(19) World Intellectual Property Organization International Bureau



(43) International Publication Date 3 October 2002 (03.10.2002)

PCT

(10) International Publication Number WO 02/077272 A2

(51) International Patent Classification7:

C12Q 1/68

(21) International Application Number: PCT/EP02/03401

(22) International Filing Date: 26 March 2002 (26.03.2002)

(25) Filing Language:

λ

English

(26) Publication Language:

English

(30) Priority Data: 60/278,333

26 March 2001 (26.03.2001) US

(71) Applicant (for all designated States except US): EPIGE-NOMICS AG [DE/DE]; Kastanienallee 24, 10435 Berlin (DE).

(72) Inventors; and

(75) Inventors/Applicants (for US only): BERLIN, Kurt [DE/DE]; Marienkäferweg 4, 14532 Stahnsdorf (DE). BRAUN, Aron [CH/DE]; Christinenstrasse 24, 10119 Berlin (DE). DISTLER, Jürgen [DE/DE]; Hewaldstrasse 2, 10825 Berlin (DE). GÜTIG, David [DE/DE]; Kastanienallee 74, 10435 Berlin (DE). HOWE, André [DE/DE]; Trelleborger Strasse 28, 13189 Berlin (DE). MÜLLER, Jürgen [DE/DE]; Schliemannstrasse 46, 10437 Berlin (DE). OLEK, Alexander [DE/DE]; Schröderstrasse 13, 10115 Berlin (DE). PIEPENBROCK, Christian [DE/DE]; Schwartzkopffstrasse 7 B, 10115 Berlin (DE). ADORJAN, Peter [DE/DE]; Dunckerstrasse 4, 10437 Berlin (DE). GRABS, Gabi [DE/DE]; Sennockstrasse 29, 12103 Berlin (DE). LESCHE, Ralf [DE/DE]; Dänenstrasse 15, 10439 Berlin (DE). LEU, Erik [DE/DE]; Mühsamstrasse 24, 10249 Berlin (DE). LEWIN, Andre [DE/DE]; Wichertstrasse 52, 10439 Berlin (DE). LIP-SCHER, Evelyne [DE/DE]; Brehmestrasse 13, 13187 Berlin (DE). MAIER, Sabine [DE/DE]; Markelstrasse 60,

12163 Berlin (DE). MODEL, Fabian [DE/DE]; Debenzerstrasse 73, 12683 Berlin (DE). MÜLLER, Volker [DE/DE]; Dossestrasse 13, 10247 Berlin (DE). OTTO, Thomas [DE/DE]; Walliserstrasse 11, 13407 Berlin (DE). PELET, Cecile [DE/DE]; Grosse Hamburger Strasse 1, 10115 Berlin (DE). SCHWOPE, Ina [DE/DE]; Isländischestrasse 16, 10439 Berlin (DE). ZIEBARTH, Heike [DE/DE]; Paul-Robeson-Strasse 20, 10439 Berlin (DE).

- (74) Agent: SCHUBERT, Klemens; Neue Promenade 5, 10178 Berlin (DE).
- (81) Designated States (national): AE, AG, AL, AM, AT, AU, AZ, BA, BB, BG, BR, BY, BZ, CA, CH, CN, CO, CR, CU, CZ, DE, DK, DM, DZ, EC, EE, ES, FI, GB, GD, GE, GH, GM, HR, HU, ID, IL, IN, IS, JP, KE, KG, KP, KR, KZ, LC, LK, LR, LS, LT, LU, LV, MA, MD, MG, MK, MN, MW, MX, MZ, NO, NZ, OM, PH, PL, PT, RO, RU, SD, SE, SG, SI, SK, SL, TJ, TM, TN, TR, TT, TZ, UA, UG, US, UZ, VN, YU, ZA, ZM, ZW.
- (84) Designated States (regional): ARIPO patent (GH, GM, KE, LS, MW, MZ, SD, SL, SZ, TZ, UG, ZM, ZW), Eurasian patent (AM, AZ, BY, KG, KZ, MD, RU, TJ, TM), European patent (AT, BE, CH, CY, DE, DK, ES, FI, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE, TR), OAPI patent (BF, BJ, CF, CG, CI, CM, GA, GN, GQ, GW, ML, MR, NE, SN, TD, TG).

Published:

- without international search report and to be republished upon receipt of that report
- with sequence listing part of description published separately in electronic form and available upon request from the International Bureau

For two-letter codes and other abbreviations, refer to the "Guidance Notes on Codes and Abbreviations" appearing at the beginning of each regular issue of the PCT Gazette.

(54) Title: METHODS AND NUCLEIC ACIDS FOR THE ANALYSIS OF HEMATOPOIETIC CELL PROLIFERATIVE DISORDERS

(57) Abstract: The present invention relates to modified and genomic sequences, to oligonucleotides and/or PNA-oligomers for detecting the cytosine methylation state of genomic DNA, as well as to a method for ascertaining genetic and/or epigenetic parameters of genes for use in the differentiation, diagnosis, treatment and/or monitoring of hematopoietic cell proliferative disorders, or the predisposition to hematopoietic cell proliferative disorders.



Methods and nucleic acids for the analysis of hematopoietic cell proliferative disorders

Field of the invention

Class prediction is of crucial importance for most therapeutic decisions in cancer, therefore prediction of tumour class is one of the most important problems in diagnostic oncology. Different classes and sub-classes of cancers respond differently to specific types of treatment. Vital therapeutic decisions such as which therapeutic regimen is used are therefore crucially dependent on the precise recognition of tumour class. Currently tumour classification generally depends on morphological, histopathological or immunological parameters, and single molecular markers (R. W. McKenna, *Clin Chem.* 46, 1252 (2000); F. R. Appelbaum, *Semin Hematol.* 36, 401 (1999)). However, none of the diagnostic procedures currently in use is sufficient to classify tumours correctly. Different methods have to be combined, and correct classification of a tumour depends on the experience of individual pathologists.

Recently, several groups have shown that precise determination of tumour class can be achieved by microarray-based expression analysis. Golub and coworkers screened the expression levels of almost 7000 genes, between 10 and 100 of which were then shown to be sufficient to distinguish between acute lymphoblastic leukaemia (ALL) and acute myeloid leukaemia (AML) (T. R. Golub et al., Science 286, 531 (1999)). In a similar approach, Alizadeh and coworkers discovered yet unknown subclasses of diffuse large B-cell lymphoma with significant differences regarding response to therapy and disease outcome (A. A. Alizadeh et al., Nature 403, 503 (2000)). Also, it was shown that global transcript analysis can predict phenotypic characteristics of malignant melanoma (M. Bittner et al., Nature 406, 535 (2000)).

Large scale analysis using mRNA based microarrays are primarily impeded by the instability of mRNA (T. Emmert-Buck et al., Am J Pathol. 156, 1109 (2000)). Also expression changes of only a minimum of a factor 2 can be routinely and reliably detected (R. J. Lipshutz, S. P. A Fodor, T. R. Gingeras, D. J. Lockhart, Nature Genetics 21, 20 (1999); D. W. Selinger, K. J. Cheung, R. Mei, E. M. Johansson, C. S. Richmond, Nature Biotechnology 18, 1262 (2000)). Furthermore, sample preparation is complicated by the fact that expression changes occur within minutes following certain triggers. The inability to resolve the individual contributions of such influences on an expression profile, and difficulties with quantifying the gradual nature of the occurring changes complicates data analysis.

Aberrant DNA methylation within CpG islands is common in human malignancies leading to

abrogation or overexpression of a broad spectrum of genes (P. A. Jones, Cancer Res 65, 2463 (1996)). Abnormal methylation has also been shown to occur in CpG rich regulatory elements in intronic and coding parts of genes for certain tumours (M. F. Chan, G. Liang, P. A. Jones, Curr Top Microbiol Immunol 249, 75(2000)). Using restriction landmark genomic scanning, Costello and coworkers were able to show that methylation patterns are tumour-type specific (J. F. Costello et al., Nat Genet 24, 132 (2000)). Highly characteristic DNA methylation patterns could also be shown for breast cancer cell lines (T. H.-M. Huang, M. R. Perry, D. E. Laux, Hum Mol Genet 8, 459 (1999)). As does large scale mRNA expression monitoring, genome wide assessment of methylation status represents a molecular fingerprint of cancer tissues and therefore should allow tumour class prediction and discovery.

Leukemia is a malignant disease that originates in a cell in the marrow, and is characterized by the uncontrolled proliferation of developing marrow cells. The majority of leukemias are classified according to the cell type from which they develop, as either myelogenous or lymphocytic. Both classes may be acute or chronic. Acute leukemia is a rapidly progressing disease that results in the accumulation of immature, functionless cells in the marrow and blood. In many cases the marrow can no longer produce sufficient quantities of normal red and white blood cells and platelets. Anemia, a deficiency of red cells, develops in virtually all leukemia patients. The lack of normal white cells impairs the body's ability to fight infections. A shortage of platelets results in bruising and easy bleeding. Chronic leukemia progresses more slowly and permits greater numbers of more mature, functional cells to be made.

Among an estimated 31,500 new cases of leukemia in the United States in 2001, about equal proportions are acute leukemia and chronic types. Most cases occur in older adults; more than half of all cases occur after age 60. Leukemia usually strikes ten times as many adults as children. Leukemia is the most common childhood cancer and acute lymphocytic leukemia (hereinafter also referred to by the abbreviation ALL) accounts for 80 percent of the childhood leukemia cases. The most common types of leukemia in adults are acute myelogenous leukemia (hereinafter also referred to by the abbreviation AML), with an estimated 10,000 new cases annually. Acute lymphocytic leukemia will account for about 3,500 cases this year.

In addition to environmental and viral causes a variety of genes involved in cellular pathways such as cell signalling, apoptosis, cell proliferation, and senescence are thought to be involved in the development of leukemia. Examples of such genes in Bcl2, ALL1,

5-methylcytosine is the most frequent covalent base modification in the DNA of eukaryotic cells. It plays a role, for example, in the regulation of the transcription, in genetic imprinting, and in tumorigenesis. Therefore, the identification of 5-methylcytosine as a component of genetic

information is of considerable interest. However, 5-methylcytosine positions cannot be identified by sequencing since 5-methylcytosine has the same base pairing behavior as cytosine. Moreover, the epigenetic information carried by 5-methylcytosine is completely lost during PCR amplification.

A relatively new and currently the most frequently used method for analyzing DNA for 5methylcytosine is based upon the specific reaction of bisulfite with cytosine which, upon subsequent alkaline hydrolysis, is converted to uracil which corresponds to thymidine in its base pairing behavior. However, 5-methylcytosine remains unmodified under these conditions. Consequently, the original DNA is converted in such a manner that methylcytosine, which originally could not be distinguished from cytosine by its hybridization behavior, can now be detected as the only remaining cytosine using "normal" molecular biological techniques, for example, by amplification and hybridization or sequencing. All of these techniques are based on base pairing which can now be fully exploited. In terms of sensitivity, the prior art is defined by a method which encloses the DNA to be analyzed in an agarose matrix, thus preventing the diffusion and renaturation of the DNA (bisulfite only reacts with single-stranded DNA), and which replaces all precipitation and purification steps with fast dialysis (Olek A, Oswald J, Walter J. A modified and improved method for bisulphite based cytosine methylation analysis. Nucleic Acids Res. 1996 Dec 15;24(24):5064-6). Using this method, it is possible to analyze individual cells, which illustrates the potential of the method. However, currently only individual regions of a length of up to approximately 3000 base pairs are analyzed, a global analysis of cells for thousands of possible methylation events is not possible. However, this method cannot reliably analyze very small fragments from small sample quantities either. These are lost through the matrix in spite of the diffusion protection.

An overview of the further known methods of detecting 5-methylcytosine may be gathered from the following review article: Rein, T., DePamphilis, M. L., Zorbas, H., Nucleic Acids Res. 1998, 26, 2255.

To date, barring few exceptions (e.g., Zeschnigk M, Lich C, Buiting K, Doerfler W, Horsthemke B. A single-tube PCR test for the diagnosis of Angelman and Prader-Willi syndrome based on allelic methylation differences at the SNRPN locus. Eur J Hum Genet. 1997 Mar-Apr;5(2):94-8) the bisulfite technique is only used in research. Always, however, short, specific fragments of a known gene are amplified subsequent to a bisulfite treatment and either completely sequenced (Olek A, Walter J. The pre-implantation ontogeny of the H19 methylation imprint. Nat Genet. 1997

Nov;17(3):275-6) or individual cytosine positions are detected by a primer extension reaction (Gonzalgo ML, Jones PA. Rapid quantitation of methylation differences at specific sites using methylation-sensitive single nucleotide primer extension (Ms-SNuPE). Nucleic Acids Res. 1997 Jun 15;25(12):2529-31, WO Patent 9500669) or by enzymatic digestion (Xiong Z, Laird PW. COBRA: a sensitive and quantitative DNA methylation assay. Nucleic Acids Res. 1997 Jun 15;25(12):2532-4). In addition, detection by hybridization has also been described (Olek et al., WO 99 28498).

Further publications dealing with the use of the bisulfite technique for methylation detection in individual genes are: Grigg G, Clark S. Sequencing 5-methylcytosine residues in genomic DNA. Bioessays. 1994 Jun;16(6):431-6, 431; Zeschnigk M, Schmitz B, Dittrich B, Buiting K, Horsthemke B, Doerfler W. Imprinted segments in the human genome: different DNA methylation patterns in the Prader-Willi/Angelman syndrome region as determined by the genomic sequencing method. Hum Mol Genet. 1997 Mar;6(3):387-95; Feil R, Charlton J, Bird AP, Walter J, Reik W. Methylation analysis on individual chromosomes: improved protocol for bisulphite genomic sequencing. Nucleic Acids Res. 1994 Feb 25;22(4):695-6; Martin V, Ribieras S, Song-Wang X, Rio MC, Dante R. Genomic sequencing indicates a correlation between DNA hypomethylation in the 5' region of the pS2 gene and its expression in human breast cancer cell lines. Gene. 1995 May 19;157(1-2):261-4; WO 97 46705, WO 95 15373 and WO 45560.

An overview of the Prior Art in oligomer array manufacturing can be gathered from a special edition of Nature Genetics (Nature Genetics Supplement, Volume 21, January 1999), published in January 1999, and from the literature cited therein.

Fluorescently labeled probes are often used for the scanning of immobilized DNA arrays. The simple attachment of Cy3 and Cy5 dyes to the 5'-OH of the specific probe are particularly suitable for fluorescence labels. The detection of the fluorescence of the hybridized probes may be carried out, for example via a confocal microscope. Cy3 and Cy5 dyes, besides many others, are commercially available.

Matrix Assisted Laser Desorption Ionization Mass Spectrometry (MALDI-TOF) is a very efficient development for the analysis of biomolecules (Karas M, Hillenkamp F. Laser desorption ionization

of proteins with molecular masses exceeding 10,000 daltons. Anal Chem. 1988 Oct 15;60(20):2299-301). An analyte is embedded in a light-absorbing matrix. The matrix is evaporated by a short laser pulse thus transporting the analyte molecule into the vapor phase in an unfragmented manner. The analyte is ionized by collisions with matrix molecules. An applied voltage accelerates the ions into a field-free flight tube. Due to their different masses, the ions are accelerated at different rates. Smaller ions reach the detector sooner than bigger ones.

MALDI-TOF spectrometry is excellently suited to the analysis of peptides and proteins. The analysis of nucleic acids is somewhat more difficult (Gut I G, Beck S. DNA and Matrix Assisted Laser Desorption Ionization Mass Spectrometry. Current Innovations and Future Trends. 1995, 1; 147-57). The sensitivity to nucleic acids is approximately 100 times worse than to peptides and decreases disproportionally with increasing fragment size. For nucleic acids having a multiply negatively charged backbone, the ionization process via the matrix is considerably less efficient. In MALDI-TOF spectrometry, the selection of the matrix plays an eminently important role. For the desorption of peptides, several very efficient matrixes have been found which produce a very fine crystallization. There are now several responsive matrixes for DNA, however, the difference in sensitivity has not been reduced. The difference in sensitivity can be reduced by chemically modifying the DNA in such a manner that it becomes more similar to a peptide. Phosphorothioate nucleic acids in which the usual phosphates of the backbone are substituted with thiophosphates can be converted into a charge-neutral DNA using simple alkylation chemistry (Gut IG, Beck S. A procedure for selective DNA alkylation and detection by mass spectrometry. Nucleic Acids Res. 1995 Apr 25;23(8):1367-73). The coupling of a charge tag to this modified DNA results in an increase in sensitivity to the same level as that found for peptides. A further advantage of charge tagging is the increased stability of the analysis against impurities which make the detection of unmodified substrates considerably more difficult.

Genomic DNA is obtained from DNA of cell, tissue or other test samples using standard methods. This standard methodology is found in references such as Fritsch and Maniatis eds., Molecular Cloning: A Laboratory Manual, 1989.

The method according to the invention presents a novel microarray-based assay which is suited for methylation analysis of very large numbers of genes and CpG dinucleotides in parallel. Samples

obtained from patients with acute lymphoblastic leukemia (ALL) or acute myeloid leukemia (AML) can be classified based solely on DNA methylation patterns.

7

Description

The invention provide a method for the analysis of biological samples for features associated with the development of hematopoietic cell proliferative disorders, characterised in that the nucleic acid of at least one member of the group comprising the genes: ABL1, ABL1, APAF1, APC, AR, ARHI, BAK1, BAX, BCL2, CASP10, CASP8, CASP9, CCND2, CDC2, CDC25A, CDH1, CDH3, CDK 4, CDKN1A, CDKN1B (p27 Kip1), CDKN1C, CDKN2a, CDKN2B, CSNK2B, DAPK1, EGR4, ELK1, ESR1, FOS, GPIb beta, GPR37, GSK3ß, GSTP1, HIC-1, HOXA5, IGF2, MDR1, MGMT, MLH1, MOS, Humos, MPL, MYC, MYCL1, MYOD1, N33, PITX2, PML, PMS2, PRAME, PTEN, RB1, RBL2, SDC4, SFN, TCL1A, TGFBR2, TP73, WT1, N-MYC, L-MYC, C-ABL, ELK1, Tubulin, CSF1, CD1R3, CSNK2B, Me491/TD63, AR, CDK 4, Humos, CDC25A, CMYCex3 is/are contacted with a reagent or series of reagents capable of distinguishing between methylated and non methylated CpG dinucleotides within the genomic sequence of interest.

The present invention makes available a method for ascertaining genetic and/or epigenetic parameters of genomic DNA. The method is for use in the improved diagnosis, treatment and monitoring of hematopoietic cell proliferative disorders, more specifically by enabling the improved identification of and differentiation between subclasses of said disorder and the genetic predisposition to said disorders. The invention presents improvements over the state of the art in that it enables a highly specific classification of hematopoietic cell proliferative disorders, thereby allowing for improved and informed treatment of patients.

In a particularly preferred embodiment the present invention makes available methods and nucleic acids that allow the differentiation between acute lymphocytic leukemia and acute myelogenous leukemia.

Furthermore, the method enables the analysis of cytosine methylations and single nucleotide polymorphisms.

In a preferred embodiment, the method comprises the following steps:

In the first step of the method the genomic DNA sample must be isolated from sources such as cell lines or blood samples. Extraction may be by means that are standard to one skilled in the art, these include the use of detergent lysates, sonification and vortexing with glass beads. Once the nucleic acids have been extracted the genomic double stranded DNA is used in the analysis.

In a preferred embodiment the DNA may be cleaved prior to the next step of the method, this may be by any means standard in the state of the art, in particular, but not limited to, with restriction endonucleases.

In the second step of the method, the genomic DNA sample is treated in such a manner that cytosine bases which are unmethylated at the 5'-position are converted to uracil, thymine, or another base which is dissimilar to cytosine in terms of hybridization behaviour. This will be understood as 'pretreatment' hereinafter.

The above described treatment of genomic DNA is preferably carried out with bisulfite (sulfite, disulfite) and subsequent alkaline hydrolysis which results in a conversion of non-methylated cytosine nucleobases to uracil or to another base which is dissimilar to cytosine in terms of base pairing behaviour. If bisulfite solution is used for the reaction, then an addition takes place at the non-methylated cytosine bases. Moreover, a denaturating reagent or solvent as well as a radical interceptor must be present. A subsequent alkaline hydrolysis then gives rise to the conversion of non-methylated cytosine nucleobases to uracil. The chemically converted DNA is then used for the detection of methylated cytosines.

Fragments of the pretreated DNA are amplified, using sets of primer oligonucleotides according to SEQ ID NO:387 to SEQ ID NO: 534, and a, preferably heat-stable, polymerase. Because of statistical and practical considerations, preferably more than ten different fragments having a length of 100 - 2000 base pairs are amplified. The amplification of several DNA segments can be carried out simultaneously in one and the same reaction vessel. Usually, the amplification is carried out by means of a polymerase chain reaction (PCR).

The method may also be enabled by the use of alternative primers, the design of such primers is obvious to one skilled in the art. These should include at least two oligonucleotides whose sequences are each reverse complementary or identical to an at least 18 base-pair long segment of the base sequences specified in the appendix (SEQ ID NO:95 through SEQ ID NO: 386). Said primer oligonucleotides are preferably characterized in that they do not contain any CpG dinucleotides. In a particularly preferred embodiment of the method, the sequence of said primer oligonucleotides are designed so as to selectively anneal to and amplify, only the hematopoietic cell specific DNA of interest, thereby minimizing the amplification of background or non relevant DNA. In the context of the present invention, background DNA is taken to mean genomic DNA which does not have a relevant tissue specific methylation pattern, in this case, the relevant tissue being

hematopoietic cells, both healthy and diseased.

According to the present invention, it is preferred that at least one primer oligonucleotide is bound to a solid phase during amplification. The different oligonucleotide and/or PNA-oligomer sequences can be arranged on a plane solid phase in the form of a rectangular or hexagonal lattice, the solid phase surface preferably being composed of silicon, glass, polystyrene, aluminum, steel, iron, copper, nickel, silver, or gold, it being possible for other materials such as nitrocellulose or plastics to be used as well.

The fragments obtained by means of the amplification can carry a directly or indirectly detectable label. Preferred are labels in the form of fluorescence labels, radionuclides, or detachable molecule fragments having a typical mass which can be detected in a mass spectrometer, it being preferred that the fragments that are produced have a single positive or negative net charge for better detectability in the mass spectrometer. The detection may be carried out and visualized by means of matrix assisted laser desorption/ionization mass spectrometry (MALDI) or using electron spray mass spectrometry (ESI).

The amplificates obtained in the second step of the method are subsequently hybridized to an array or a set of oligonucleotides and/or PNA probes. In this context, the hybridization takes place in the manner described as follows. The set of probes used during the hybridization is preferably composed of at least 10 oligonucleotides or PNA-oligomers. In the process, the amplificates serve as probes which hybridize to oligonucleotides previously bonded to a solid phase. In a particularly preferred embodiment, the oligonucleotides are taken from the group comprising SEQ ID NO: 535 to SEQ ID NO: 1258. In a further preferred embodiment as described in the examples, the oligonucleotides are taken from the group comprising SEQ ID NO: 1211 to SEQ ID NO: 1258. The non-hybridized fragments are subsequently removed. Said oligonucleotides contain at least one base sequence having a length of 10 nucleotides which is reverse complementary or identical to a segment of the base sequences specified in the appendix, the segment containing at least one CpG or TpG dinucleotide. In a further preferred embodiment the cytosine of the CpG dinucleotide, or in the case of TpG, the thiamine, is the 5th to 9th nucleotide from the 5th-end of the 10-mer. One oligonucleotide exists for each CpG or TpG dinucleotide.

In the fifth step of the method, the non-hybridized amplificates are removed.

In the final step of the method, the hybridized amplificates are detected. In this context, it is preferred that labels attached to the amplificates are identifiable at each position of the solid phase at which an oligonucleotide sequence is located.

According to the present invention, it is preferred that the labels of the amplificates are fluorescence labels, radionuclides, or detachable molecule fragments having a typical mass which can be detected in a mass spectrometer. The mass spectrometer is preferred for the detection of the amplificates, fragments of the amplificates or of probes which are complementary to the amplificates, it being possible for the detection to be carried out and visualized by means of matrix assisted laser desorption/ionization mass spectrometry (MALDI) or using electron spray mass spectrometry (ESI). The produced fragments may have a single positive or negative net charge for better detectability in the mass spectrometer.

The aforementioned method is preferably used for ascertaining genetic and/or epigenetic parameters of genomic DNA.

In order to enable this method, the invention further provides the modified DNA of genes ABL1, ABL1, APAF1, APC, AR, ARHI, BAK1, BAX, BCL2, CASP10, CASP8, CASP9, CCND2, CDC2, CDC25A, CDH1, CDH3, CDK 4, CDKN1A, CDKN1B (p27 Kip1), CDKN1C, CDKN2a, CDKN2B, CSNK2B, DAPK1, EGR4, ELK1, ESR1, FOS, GPIb beta, GPR37, GSK3ß, GSTP1, HIC-1, HOXA5, IGF2, MDR1, MGMT, MLH1, MOS, Humos, MPL, MYC, MYCL1, MYOD1, N33, PITX2, PML, PMS2, PRAME, PTEN, RB1, RBL2, SDC4, SFN, TCL1A, TGFBR2, TP73, WT1, N-MYC, L-MYC, C-ABL, ELK1, Tubulin, CSF1, CD1R3, CSNK2B, Me491/TD63, AR, CDK 4, Humos, CDC25A, CMYCex3 as well as oligonucleotides and/or PNA-oligomers for detecting cytosine methylations within said genes. The present invention is based on the discovery that genetic and epigenetic parameters and, in particular, the cytosine methylation patterns of genomic DNA are particularly suitable for improved diagnosis, treatment and monitoring of hematopoietic cell proliferative disorders. Furthermore, the invention enables the differentiation between different subclasses of hematopoietic cell proliferative disorders or detection of a predisposition to hematopoietic cell proliferative disorders.

The nucleic acids according to the present invention can be used for the analysis of genetic and/or epigenetic parameters of genomic DNA.

11

This objective is achieved according to the present invention using a nucleic acid containing a sequence of at least 18 bases in length of the pretreated genomic DNA according to one of SEQ ID NO: 95 through SEQ ID NO: 386 and sequences complementary thereto.

The modified nucleic acid could heretofore not be connected with the ascertainment of disease relevant genetic and epigenetic parameters.

The object of the present invention is further achieved by an oligonucleotide or oligomer for the analysis of pretreated DNA, for detecting the genomic cytosine methylation state, said oligonucleotide containing at least one base sequence having a length of at least 10 nucleotides which hybridizes to a pretreated genomic DNA according to SEQ ID NO: 95 through SEQ ID NO: 386. The oligomer probes according to the present invention constitute important and effective tools which, for the first time, make it possible to ascertain specific genetic and epigenetic parameters during the analysis of biological samples for features associated with the development of hematopoietic cell proliferative disorders. Said oligonucleotides allow the improved diagnosis, treatment and monitoring of hematopoietic cell proliferative disorders and detection of the predisposition to said disorders. Furthermore, they allow the differentiation of different subclasses of hematopoietic carcinomas. The base sequence of the oligomers preferably contains at least one CpG or TpG dinucleotide. The probes may also exist in the form of a PNA (peptide nucleic acid) which has particularly preferred pairing properties. Particularly preferred are oligonucleotides according to the present invention in which the cytosine of the CpG dinucleotide is the 5th - 9th nucleotide from the 5'-end of the 13-mer; in the case of PNA-oligomers, it is preferred for the cytosine of the CpG dinucleotide to be the 4th - 6th nucleotide from the 5'-end of the 9-mer.

The oligomers according to the present invention are normally used in so called "sets" which contain at least one oligomer for each of the CpG dinucleotides within SEQ ID NO: 95 through SEQ ID NO: 386. Preferred is a set which contains at least one oligomer for each of the CpG dinucleotides, from SEQ ID NO: 535 through SEQ ID NO: 1258. Further preferred is a set comprising SEQ ID NO: 1211 to SEQ ID NO: 1258.

In the case of the sets of oligonucleotides according to the present invention, it is preferred that at least one oligonucleotide is bound to a solid phase. It is further preferred that all the oligonucleotides of one set are bound to a solid phase.

The present invention moreover relates to a set of at least 10 n (oligonucleotides and/or PNA-oligomers) used for detecting the cytosine methylation state of genomic DNA using treated versions of said genomic DNA (according to SEQ ID NO: 95 through SEQ ID NO: 386 and sequences complementary thereto). These probes enable improved diagnosis, treatment and monitoring of hematopoietic cell proliferative disorders. In particular they enable the differentiation between different sub classes of hematopoietic cell proliferative disorders and the detection of a predisposition to said disorders. In a particularly preferred embodiment the set comprises SEQ ID NO: 74 to SEQ ID NO: 1258.

The set of oligomers may also be used for detecting single nucleotide polymorphisms (SNPs) using pretreated genomic DNA according to one of SEQ ID NO: 95 through SEQ ID NO: 386.

According to the present invention, it is preferred that an arrangement of different oligonucleotides and/or PNA-oligomers (a so-called "array") made available by the present invention is present in a manner that it is likewise bound to a solid phase. This array of different oligonucleotide- and/or PNA-oligomer sequences can be characterized in that it is arranged on the solid phase in the form of a rectangular or hexagonal lattice. The solid phase surface is preferably composed of silicon, glass, polystyrene, aluminum, steel, iron, copper, nickel, silver, or gold. However, nitrocellulose as well as plastics such as nylon which can exist in the form of pellets or also as resin matrices are suitable alternatives.

Therefore, a further subject matter of the present invention is a method for manufacturing an array fixed to a carrier material for the improved diagnosis, treatment and monitoring of hematopoietic cell proliferative disorders, the differentiation between different subclasses of hematopoietic carcinomas and/or detection of the predisposition to hematopoietic cell proliferative disorders. In said method at least one oligomer according to the present invention is coupled to a solid phase. Methods for manufacturing such arrays are known, for example, from US Patent 5,744,305 by

means of solid-phase chemistry and photolabile protecting groups.

A further subject matter of the present invention relates to a DNA chip for the improved diagnosis, treatment and monitoring of hematopoietic cell proliferative disorders. Furthermore the DNA chip enables detection of the predisposition to hematopoietic cell proliferative disorders and the differentiation between different subclasses of hematopoietic carcinomas. The DNA chip contains at least one nucleic acid according to the present invention. DNA chips are known, for example, in US Patent 5,837,832.

Moreover, a subject matter of the present invention is a kit which may be composed, for example, of a bisulfite-containing reagent, a set of primer oligonucleotides containing at least two oligonucleotides whose sequences in each case correspond to or are complementary to a 18 base long segment of the base sequences specified in the appendix (SEQ ID NO: 95 through SEQ ID NO: 386), oligonucleotides and/or PNA-oligomers as well as instructions for carrying out and evaluating the described method. However, a kit along the lines of the present invention can also contain only part of the aforementioned components.

The oligomers according to the present invention or arrays thereof as well as a kit according to the present invention are intended to be used for the improved diagnosis, treatment and monitoring of hematopoietic cell proliferative disorders. Furthermore the use of said inventions extends to the differentiation between different subclasses of hematopoietic carcinomas and detection of the predisposition to hematopoietic cell proliferative disorders. According to the present invention, the method is preferably used for the analysis of important genetic and/or epigenetic parameters within genomic DNA, in particular for use in improved diagnosis, treatment and monitoring of hematopoietic cell proliferative disorders, detection of the predisposition to said disorders and the differentiation between subclasses of said disorders.

The methods according to the present invention are used, for example, for improved diagnosis, treatment and monitoring of hematopoietic cell proliferative disorders progression, detection of the predisposition to said disorders and the differentiation between subclasses of said disorders.

A further embodiment of the invention is a method for the analysis of the methylation status of genomic DNA without the need for pretreatment. In the first step of the method the genomic DNA

sample must be isolated from tissue or cellular sources. Such sources may include cell lines, histological slides, body fluids, or tissue embedded in paraffin. Extraction may be by means that are standard to one skilled in the art, these include the use of detergent lysates, sonification and vortexing with glass beads. Once the nucleic acids have been extracted the genomic double stranded DNA is used in the analysis.

In a preferred embodiment the DNA may be cleaved prior to the treatment, this may be any means standard in the state of the art, in particular with restriction endonucleases. In the second step, the DNA is then digested with one or more methylation sensitive restriction enzymes. The digestion is carried out such that hydrolysis of the DNA at the restriction site is informative of the methylation status of a specific CpG dinucleotide.

In the third step the restriction fragments are amplified. In a preferred embodiment this is carried out using a polymerase chain reaction.

In the final step the amplificates are detected. The detection may be by any means standard in the art, for example, but not limited to, gel electrophoresis analysis, hybridisation analysis, incorporation of detectable tags within the PCR products, DNA array analysis, MALDI or ESI analysis.

The analysis of the methylation status of CpG positions using either of the molecular biological techniques described above in multiple patient samples generates a large set of data points which are then analysed using the following techniques in order to define CpG marker positions fpr a particular phenotypic class, in particular by using supervised machine learning techniques.

The problem of class prediction was addressed by using a supervised machine learning method called support vector machine (SVM) (V. Vapnik, Statistical Learning Theory (Wiley, New York 1998); N. Christianini, J. Shawe-Taylor, An Introduction to Support Vector Machines (Cambridge University Press, Cambridge 2000)), which has already been successfully applied to the analysis of microarray gene expression data (M. P. Brown et al., Proc Natl Acad Sci U S A 97, 262(2000); T. Gaasterland, S. Bekiranov, Nature Genetics 24, 204(2000)). The SVM constructs an optimal discriminant between two classes of given training samples. In this case each sample is described by the methylation patterns (CG/TG ratios) at the investigated CpG sites. A supervised learning technique, such as SVM, has the advantage that it exploits the prior knowledge represented by data labels. Furthermore, SVMs are capable to account for highly non-linear interdependencies between individual dimensions (in this case CpG sites) and still avoid overfitting the data (V. Vapnik, Statistical Learning Theory (Wiley, New York 1998)).

In order to quantify the contribution of single CpG sites to the prediction of classes, we initially ranked all CpG sites according to their discriminatory power for the investigated sample classes using a two sample t-test. Then the SVM was trained using an increasing number of CpG sites in the order of their ranking. The number of CpG sites where the prediction error is minimized depends on which particular classes are to be separated from each other. This is explained by the fact that, on the one hand, the exclusion of certain CpG positions from the analysis potentially leads to a loss of information. On the other hand, as suggested by statistical learning theory (V. Vapnik, Statistical Learning Theory (Wiley, New York 1998)), for a limited number of samples the reliability of the class-predictor generally increases with decreasing number of free parameters, in our case number of CpG sites. The inherent, task-specific complexity of the classification problem determines where these two counteracting effects are balanced and class prediction achieves its optimal performance. This suggests, that a proper feature selection in genome-wide analyses not only improves class prediction J. Weston et al., in Advances in neural information processing systems (MIT Press, Cambridge, MA 2001), but also provides information about the complexity of the problem.

The present invention moreover relates to the diagnosis and/or prognosis of events which are disadvantageous or relevant to patients or individuals in which important genetic and/or epigenetic parameters within genomic DNA, said parameters obtained by means of the present invention may be compared to another set of genetic and/or epigenetic parameters, the differences serving as the basis for the diagnosis and/or prognosis of events which are disadvantageous or relevant to patients or individuals.

In the context of the present invention the term "hybridization" is to be understood as a bond of an oligonucleotide to a completely complementary sequence along the lines of the Watson-Crick base pairings in the sample DNA, forming a duplex structure.

In the context of the present invention, "genetic parameters" are mutations and polymorphisms of genomic DNA and sequences further required for their regulation. To be designated as mutations

are, in particular, insertions, deletions, point mutations, inversions and polymorphisms and, particularly preferred, SNPs (single nucleotide polymorphisms).

In the context of the present invention, "epigenetic parameters" are, in particular, cytosine methylations and further modifications of DNA bases of genomic DNA and sequences further required for their regulation. Further epigenetic parameters include, for example, the acetylation of histones which, cannot be directly analyzed using the described method but which, in turn, correlates with the DNA methylation.

In the following, the present invention will be explained in greater detail on the basis of the sequences and examples without being limited thereto.

SEQ ID NO: 1 to SEQ ID NO: 73 represent 5' and/or regulatory regions of the genomic DNA of genes: ABL1, ABL1, APAF1, APC, AR, ARHI, BAK1, BAX, BCL2, CASP10, CASP8, CASP9, CCND2, CDC2, CDC25A, CDH1, CDH3, CDK 4, CDKN1A, CDKN1B (p27 Kip1), CDKN1C, CDKN2a, CDKN2B, CSNK2B, DAPK1, EGR4, ELK1, ESR1, FOS, GPIb beta, GPR37, GSK3ß, GSTP1, HIC-1, HOXA5, IGF2, MDR1, MGMT, MLH1, MOS, Humos, MPL, MYC, MYCL1, MYOD1, N33, PITX2, PML, PMS2, PRAME, PTEN, RB1, RBL2, SDC4, SFN, TCL1A, TGFBR2, TP73, WT1, N-MYC, L-MYC, C-ABL, ELK1, Tubulin, CSF1, CD1R3, CSNK2B, Me491/TD63, AR, CDK 4, Humos, CDC25A, CMYCex3. These sequences are derived from Genbank and will be taken to include all minor variations of the sequence material which are currently unforeseen, for example, but not limited to, minor deletions and SNPs.

SEQ ID NO: 95 to SEQ ID NO: 386 exhibit the pretreated sequence of DNA derived from genes: ABL1, ABL1, APAF1, APC, AR, ARHI, BAK1, BAX, BCL2, CASP10, CASP8, CASP9, CCND2, CDC2, CDC25A, CDH1, CDH3, CDK 4, CDKN1A, CDKN1B (p27 Kip1), CDKN1C, CDKN2a, CDKN2B, CSNK2B, DAPK1, EGR4, ELK1, ESR1, FOS, GPIb beta, GPR37, GSK3ß, GSTP1, HIC-1, HOXA5, IGF2, MDR1, MGMT, MLH1, MOS, Humos, MPL, MYC, MYCL1, MYOD1, N33, PITX2, PML, PMS2, PRAME, PTEN, RB1, RBL2, SDC4, SFN, TCL1A, TGFBR2, TP73, WT1, N-MYC, L-MYC, C-ABL, ELK1, Tubulin, CSF1, CD1R3, CSNK2B, Me491/TD63, AR, CDK 4, Humos, CDC25A, CMYCex3. These sequences will be taken to include all minor variations of the sequence material which are currently unforeseen, for example, but not limited to, minor deletions and SNPs.

SEQ ID NO: 387 to SEQ ID NO: 534 exhibit the sequence of primer oligonucleotides for the amplification of pretreated DNA according to SEQ ID NO: 95 to SEQ ID NO: 386.

SEQ ID NO: [&IDOLIGOFIRST] to SEQ ID NO: 1258 exhibit the sequence of oligomers which are useful for the analysis of CpG positions within genomic DNA according to SEQ ID NO: 1 to SEQ ID NO: 73.

SEQ ID NO: 1211 to SEQ ID NO: 1258 exhibit the sequence of oligomers which are useful for the analysis of CpG positions within genomic DNA according to SEQ ID NO: 1 to SEQ ID NO: 73.

Examples 1 and 2: Digital Phenotype

In the following examples, multiplex PCR was carried out on samples from patients with acute lymphocytic leukemia and acute myelogenous leukemia. Each sample was treated in the manner described below in Example 1 in order to deduce the methylation status of CpG positions, the CpG methylation information for each sample was collated and then used in an analysis, as detailed in Example 2. An alternative method for the analysis of CpG methylation status is further described in Example 3.

Example 1

A total of 17 ALL samples and 8 AML samples were included. T- and B-cells from 8 healthy donors, MACS-sorted with anti-CD4 and anti-CD19 antibodies, respectively, served as controls. 11 genes were randomly selected from a panel of genes representing different pathways associated with tumourigenesis. Two of the 11 selected genes are located on the X-chromosome. A total of 81 CpG dinucleotides located in CpG rich regions of the promoters, intronic and coding sequences of these genes were evaluated for methylation status.

In order to allow sequence specific distinction of the methylated from the unmethylated state of CpG dinucleotides by hybridisation analysis, total DNA from all samples was bisulphite treated converting all unmethylated cytosines to uracil whereas methylated cytosines were conserved (M. Frommer et al., Proc Natl Acad Sci USA 89, 1827 (1992)). Regions of interest were then amplified by PCR using fluorescently labelled primers converting originally unmethylated CpG dinucleotides to TG and conserving originally methylated CpG sites. Primers were designed complementary to DNA segments containing no CpG dinucleotides. This allowed unbiased amplification of both methylated and unmethylated alleles in one reaction. All PCR products performed on an individual sample were mixed and hybridised to glass slides carrying for each CpG position a pair of immobilised oligonucleotides. Each of these detection oligonucleotides was designed to hybridise to the bisulphite converted sequence around one CpG site which was either originally unmethylated (TG) or methylated (CG). Hybridisation conditions were selected to allow the detection of the single nucleotide differences between the TG and CG variants. Ratios for the two signals were calculated based on comparison of intensity of the fluorescent signals. Sensitivity for detection of methylation changes was determined using artificially up- and downmethylated DNA fragments mixed at different ratios. For each of those mixtures, a series of experiments was conducted to define the range of CG/TG ratios that corresponds to varying degrees of methylation at each of the

CpG sites tested.

In the first step the genomic DNA was isolated from the cell samples using the Wizzard kit from (Promega).

The isolated genomic DNA from the samples are treated using a bisulfite solution (hydrogen sulfite, disulfite). The treatment is such that all non methylated cytosines within the sample are converted to thiamidine, conversely 5-methylated cytosines within the sample remain unmodified. Bisulphite treatment of genomic DNA was done with minor modifications as described by A. Olek, J. Oswald, J. Walter, *Nucleic Acid Res.* 24, 5064 (1996). Genomic DNA was digested with *MssI* (MBI Fermentas, St. Leon-Rot, Germany) prior to the modification by bisulphite.

The treated nucleic acids were then amplified using multiplex PCRs, amplifying 8 fragments per reaction with Cy5 fluorescently labelled primers. PCR primers used are described in Table 1. PCR conditions were as follows.

For the PCR amplification of the bisulphite treated sense strand of the genes used for class prediction and discovery the primers were designed according to the guidelines of Clark and Frommer (S. J. Clark, M. Frommer, in *Laboratory Methods for the Detection of Mutations and Polymorphisms in DNA*, G. R. Taylor ed. (CRC Press, Boca Raton 1997)). CpG sites from the following genes were analysed: ABL1, ABL1, APAF1, APC, AR, ARHI, BAK1, BAX, BCL2, CASP10, CASP8, CASP9, CCND2, CDC2, CDC25A, CDH1, CDH3, CDK 4, CDKN1A, CDKN1B (p27 Kip1), CDKN1C, CDKN2a, CDKN2B, CSNK2B, DAPK1, EGR4, ELK1, ESR1, FOS, GPIb beta, GPR37, GSK3ß, GSTP1, HIC-1, HOXA5, IGF2, MDR1, MGMT, MLH1, MOS, Humos, MPL, MYC, MYCL1, MYOD1, N33, PITX2, PML, PMS2, PRAME, PTEN, RB1, RBL2, SDC4, SFN, TCL1A, TGFBR2, TP73, WT1, N-MYC, L-MYC, C-ABL, ELK1, Tubulin, CSF1, CD1R3, CSNK2B, Me491/TD63, AR, CDK 4, Humos, CDC25A, CMYCex3.

10 ng DNA was used as template DNA for the PCR reactions. The template DNA, 12.5 pmol or 40 pmol (CY5-labelled) of each primer, 0.5-2 U Taq polymerase (HotStarTaq, Qiagen, Hilden, Germany) and 1 mM dNTPs were incubated with the reaction buffer supplied with the enzyme in a total volume of 20 µl. After activation of the enzyme (15 min, 96 °C) the incubation times and temperatures were 95°C for 1 min followed by 34 cycles (95°C for 1 min, annealing temperature (see Supplementary information) for 45 sec, 72°C for 75 sec) and 72°C for 10 min.

See Table 2 for further details of PCR products.

All PCR products from each individual sample were then hybridised to glass slides carrying a pair

of immobilised oligonucleotides for each CpG position under analysis. Each of these detection oligonucleotides was designed to hybridise to the bisulphite converted sequence around one CpG site which was either originally unmethylated (TG) or methylated (CG). See Table 3 for further details of all hybridisation oligonucleotides used (both informative and non-informative) Hybridisation conditions were selected to allow the detection of the single nucleotide differences between the TG and CG variants.

Oligonucleotides with a C6-amino modification at the 5'end were spotted with 4-fold redundancy on activated glass slides (T. R. Golub *et al.*, *Science* 286, 531 (1999)). For each analysed CpG position two oligonucleotides N₍₂₋₁₆₎-CG-N₍₂₋₁₆₎ and N₍₂₋₁₆₎-TG-N₍₂₋₁₆₎, reflecting the methylated and non methylated status of the CpG dinucleotides, were spotted and immobilised on the glass array. The oligonucleotide microarrays representing 81 CpG sites were hybridised with a combination of up to 11 Cy5-labelled PCR fragments as described earlier (D. Chen, Z. Yan, D. L. Cole, G. S. Srivatsa, *Nucleic Acid Res* 27, 389 (1999)).. Subsequently, the fluorescent images of the hybridised slides were obtained using a GenePix 4000 microarray scanner (Axon Instruments). Hybridisation experiments were repeated at least three times for each sample.

The sensitivity of the method for detection of methylation changes was determined using artificially up- and downmethylated DNA fragments mixed at different ratios. For each of those mixtures, a series of experiments was conducted to define the range of CG/TG ratios that corresponds to varying degrees of methylation at each of the CpG sites tested. In Fig. 1B results for two CpG positions located in exon 14 of the human factor VIII gene are shown as examples. For the mixtures of 3:0, 2:1, 1:2 and 0:3 the degree of methylation of the individual CpG sites could safely be distinguished.

To verify the detection of methylation changes in the real data set two X-chromosomal genes were included in the gene set. Because one of the two X-chromosomes in females becomes inactivated by methylation we can expect a higher degree of methylation of X-chromosomal genes in females compared to males. In Fig. 2A CpGs are ranked according to the significance of the difference between male and female methylation levels. As expected, the X-chromosomal genes (ELK1, AR) show a significantly higher methylation for females. This clearly demonstrates that the method really detects changes in methylation.

Support Vector Machines

In our case, the task of cancer classification consists of constructing a machine that can predict the leukemia subtype (ALL or AML) from a patients methylation pattern. For every patient sample this pattern is given as a vector of average (Every hybridisation experiment was at least 3 times repeated and the results averaged.) log CG/TG ratios at 81 CpG positions. Based on a given set of training examples $X = \{x^i : x^i : i : R^n\}$ with known diagnosis $Y = \{y^i : y^i : \{ALL,AML\}\}$ a discriminant function $f : R^n \otimes \{ALL,AML\}$, where n is the number of CpGs, has to be learned. The number of misclassifications of f on the training set $\{X,Y\}$ is called training error and is usually minimised by the learning machine during the training phase. However, what is of practical interest is the capability to predict the class of previously unseen samples, the so called generalisation performance of the learning machine. This performance is usually estimated by the test error, which is the number of misclassifications on an independent test set $\{X',Y'\}$.

The major problem of training a learning machine with good generalisation performance is to find a discriminant function f which on the one hand is complex enough to capture the essential properties of the data distribution, but which on the other hand avoids over-fitting the data. The Support Vector Machine (SVM) tries to solve this problem by constructing a linear discriminant that separates the training data and maximises the distance to the nearest points of the training set. This maximum margin separating hyperplane minimises the ratio between the radius of the minimum enclosing sphere of the training set and the margin between hyperplane and training points. This corresponds to minimising the so called radius margin bound on the expected probability of a test error and promises good generalisation performance (Vapnik V. "Statistical Learning Theory." Wiley, New York (1998)).

Of course there are more complex classification problems, where the dependence between class labels y^i and features x^i is not linear and the training set can not be separated by a hyperplane. In order to allow for non-linear discriminant functions the input space can be non-linearly mapped into a potentially higher dimensional feature space by a mapping function $F: x^i \otimes F(x^i)$. Because the SVM algorithm in its dual formulation uses only the inner product between elements of the input

WO 02/077272

space, the knowledge of the kernel function $K(x^i, x^j) = \langle F(x^i) \cdot F(x^j) \rangle$ is sufficient to train the SVM. It is not necessary to explicitly know the mapping F and a non-linear SVM can be trained efficiently by computing only the kernel function. Here we will only use the linear kernel $K(x^i, x^j) = \langle x^i \cdot x^j \rangle$

and the quadratic kernel $K(x^i, x^j) = (\langle x^{i x_j} \rangle + 1)^2$. In the next section we will compare SVMs trained on different feature sets. In order to evaluate the prediction performance of these SVMs we used a cross-validation method (Bishop, C.M. "Neural networks for pattern recognition." Oxford University Press, New York (1995)). For each classification task, the samples were partitioned into 8 groups of approximately equal size. Then the SVM predicted the class for the test samples in one group after it had been trained using the 7 other groups. The number of misclassifications was counted over 8 runs of the SVM algorithm for all possible choices of the test group. To obtain a reliable estimate for the test error the number of misclassifications were averaged over 50 different partitionings of the samples into 8 groups.

Feature Selection

The simplest way for applying a SVM to our methylation data is to use every CpG position as a separate dimension, not making any assumption about the interdependence of CpG sites from the same gene. On the leukemia subclassification task the SVM with linear kernel trained on this 81 dimensional input space had an average test error of 16%. Using a quadratic kernel did not significantly improve the results.

An obvious explanation for this relatively poor performance is that we have only 25 data points (even less in the training set) in a 81 dimensional space. Finding a separating hyperplane under these conditions is a heavily under-determined problem. And as it turns out, the SVM technique of maximising the margin is not sufficient to find the solution with optimal generalisation properties. It is necessary to reduce the dimensionality of the input space while retaining the relevant information for classification. This should be possible because it can be expected that only a minority of CpG positions has any connection with the two subtypes of leukemia.

The probably most popular method for dimension reduction is principle component analysis (PCA) (Bishop, C.M. "Neural networks for pattern recognition." Oxford University Press, New York (1995)). For a given training set X, PCA constructs a set of orthogonal vectors (principle components) which correspond to the directions of maximum variance. The projection of X onto the first k principle components gives the 2-norm optimal representation of X in a k-dimensional orthogonal subspace. Because this projection does not explicitly use the class information Y, PCA is an unsupervised learning technique.

In order to reduce the dimension of the input space for the SVM we performed a PCA on the combined training and test set {X,X'} and projected both sets on the first k principle components. This gives considerably better results than performing PCA only on the training set X and is justified by the fact that no label information is used. However, the generalisation results for k=2 and k=5, were even worse than for the SVM without feature selection. The reason for this is that PCA does not necessarily extract features that are important for the discrimination between ALL and AML. It first picks the features with the highest variance, which are in this case discriminating between cell lines and primary patient tissue (see Fig. 4, right, second line), i.e. subgroups that are not relevant to the classification task. As is shown in Fig. 6, features carrying information about the leukemia subclasses appear only from the 9th principle component on. The generalisation performance including the 9th component is significantly better than for a SVM without feature selection. However, it seems clear that a supervised feature selection method, which takes the class labels of the training set into account, should be more reliable and give better generalisation.

Fig. 4A shows the methylation profiles of the best 20 CpGs according to the Fisher criterion.

Fisher Criterion and t-Test

A classical measure to asses the degree of separation between two classes is given by the Fisher criterion cite (Bishop, C.M. "Neural networks for pattern recognition." Oxford University Press, New York (1995)). In our case it gives the discriminative power of the kth CpG as

where $(m_k)^{ALL/AML}$ is the mean and $(s_k)^{ALL/AML}$ is the standard deviation of all $(k_k)^i$ with $y_i = ALL/AML$. The Fisher criterion gives a high ranking for CpGs where the two classes are far apart compared to the within class variances. A very similar criterion was used by Golub and coworkers for their ALL/AML classification based on mRNA expression data (T. R. Golub *et al.*, *Science* 286, 531 (1999)). Fig. 4B hows the methylation profiles of the best 20 CpGs according to the Fisher criterion.

Another approach to rank CpGs by their discriminative power is to use a test statistic for computing the significance of class differences. Here we assumed normal distribution of the two classes and used a two sample t-test to rank the CpGs according to the significance of the difference between the class means (Mendenhall, W., Sincich, T., "Statistics for engineering and the sciences." Fig. 4 (right, first line) shows the ranking, which is very similar to the Fisher criterion because a large mean difference and a small within class variance are again the important factors.

In order to improve classification performance we trained SVMs on the k highest ranking CpGs according to the Fisher criterion or t-test. Fig. 5 shows a trained SVM on the best two CpGs from the Fisher criterion. The test errors for k=2 and k=5 are given in the table above. The results show a dramatic improvement of generalisation performance. Using Fisher criterion for feature selection and k=5 CpGs the test error was decreased to 3% compared to 16% for the SVM without feature selection. Fig. 6 shows the dependence of generalisation performance from the selected dimension k and indicates that especially the Fisher criterion gives dimension independent good generalisation for reasonable small k.

Although the two described CpG ranking methods give very good generalisation, they have some potential drawbacks. One problem is that they can only detect linear dependencies between features and class labels. A simple XOR or even OR combination of two CpGs would be completely missed. Another drawback is that redundant features are not removed. In our case there are usually several CpGs from the same gene which have a high likelihood of comethylation. This can result in a large set of high ranking features which carry essentially the same information. Although the good results seem to indicate that the described problems do not appear in our data set, they should be considered.

Backward Elimination

PCA, Fisher criterion and t-test construct or rank features independent of the learning machine that does the actual classification and are therefore called filter methods (Blum, A., Langley, P., "Selection of relevant features and examples in machine learning." Artificial Intelligence 97, 245-271 (1997)). Another approach is to use the learning machine itself for feature selection. These techniques are called wrapper methods and try to identify the features that are important for the generalisation capability of the machine. Here we propose to use the features that are important for achieving a low training error as a simple approximation. In the case of a SVM with linear kernel these features are easily identified by looking at the normal vector w of the separating hyperplane. The smaller the angle between a feature basis vector and the normal vector the more important is the feature for the separation. Features orthogonal to the normal vector have obviously no influence on the discrimination at all. This means the feature ranking is simply given by the components of the normal vector as $(w_k)^2$. Of course this ranking is not very realistic because the SVM solution on the full feature set is far from optimal as we demonstrated in the last subsections. A simple heuristic is to assume that the feature with the smallest $(w_k)^2$ is really unimportant for the solution and can be safely removed from the feature set.

Then the SVM can be retrained on the reduced feature set and the procedure is repeated until the feature set is empty. Such a successive feature removal is called backward elimination (Blum, A., Langley, P. "Selection of relevant features and examples in machine learning." Artificial Intelligence 97, 245-271 (1997)). The resulting CpG ranking on our data set is shown in Fig. 4D and differs considerably from the Fisher and t-test rankings. It seems backward elimination is able to remove redundant features. However, as shown in the table above and Fig. 6 the generalisation results are not better than for the Fisher criterion. Furthermore, backward elimination seems to be more dimension dependent and it is computationally more expensive. It follows that at least for this data set the simple Fisher criterion is the preferable feature selection technique.

A canonical way to construct a wrapper method for feature selection is to evaluate the generalisation performance of the learning machine on every possible feature subset. Cross-validation on the training set can be used to estimate the generalisation of the machine on a given feature set. What makes this exhaustive search of the feature space practically useless is the enormous number of 2ⁿ different feature combinations and there are numerous heuristics to search the feature space more efficiently (e.g. backward elimination) (Blum, A., Langley, P. "Selection of relevant features and examples in machine learning." Artificial Intelligence 97, 245-271 (1997)).

Here we only want to demonstrate that there are no higher order correlations between features and class labels in our data set. In order to do this we exhaustively searched the space of all two feature combinations. For every of the 3240 two CpG combinations of k = 2 CpGs at the same time without repetitions, out of all together n = 81 CpGs we computed the leave-one-out cross-validation error of a SVM with quadratic kernel on the training set. From all CpG pairs with minimum leave-one-out error we selected the one with the smallest radius margin ratio. This pair was considered to be the optimal feature combination and was used to evaluate the generalisation performance of the SVM on the test set.

The average test error of the exhaustive search method was with 6% the same as the one of the Fisher criterion in the case of two features and a quadratic kernel. For five features the exhaustive computation is already infeasible. In the absolute majority of cross-validation runs the CpGs selected by exhaustive search and Fisher criterion were identical. In some cases suboptimal CpGs were chosen by the exhaustive search method. These results clearly demonstrate that there are no higher order feature combinations in our data set that are important for an ALL/AML discrimination.

Example 2

The data obtained according to Example 1 is sorted into a ranked matrix (Figures 3 to 8) according to CpG methylation differences between the two classes of diseased and healthy hematopoietic cells using an algorithim. The most significant CpG positions are at the bottom of the matrix with significance decreasing towards the top. On the right side of the matrix p values for the individual CpG positions are shown. The p values are the probabilities that the observed distribution occurred by chance in the data set.

As discussed above, for selected distinctions, we trained a learning algorithm (support vector machine, SVM). The SVM (as discussed by F. Model, P. Adorjan, A. Olek, C. Piepenbrock, Feature selection for DNA methylation based cancer classification. Bioinformatics. 2001 Jun; 17 Suppl 1:S157-64) constructs an optimal discriminant between two classes of given training samples. In this case each sample is described by the methylation patterns (CG/TG ratios) at the investigated CpG sites. The SVM was trained on a subset of samples of each class, which were presented with the diagnosis attached. Independent test samples, which were not shown to the SVM before were then presented to evaluate, if the diagnosis can be predicted correctly based on the predictor created in the training round.

A cross validation method was used to evaluate the prediction performance of the SVM (C. M. Bishop, Neural networks for pattern recognition (Oxford University Press, New York 1995)). For each classification task, the samples were partitioned into 8 groups of approximately equal size. Then the SVM predicted the class for the test samples in one group after it had been trained using the 7 other groups. The number of misclassifications was counted over 8 runs of the SVM algorithm for all possible choices of the test group. To obtain a reliable estimate for the test error, i.e. the probability of misclassification for a previously unknown sample, the number of misclassifications were averaged over 10 different partitionings of the samples into 8 groups.

The SVM was trained to recognise the difference between healthy CD19+ B cells and CD4+ T cells, and T and B cell leukaemias (from both patient samples and cell lines). Samples could be classified with 15% test error using the two most informative CpG positions. Remarkably, the test error could be reduced to 4 % by including a total of 54 CpG positions into the analysis. Individual CpG sites were ranked according to their contribution to the decision of the support vector machine, showing that the decision between healthy T and B cells and ALL was primarily based on CpG sites located in intron 1 of the CDK4 gene and the coding sequence of the c-MOS oncogene, but CpG sites located in regulatory regions of other genes contributed significantly (Table 1 and Figures 2B and 2C). Compared to the healthy group, the tumour cells consistently showed relative hypermethylation

of these particular CpG sites. Also, classification could be achieved between ALL patient samples and healthy donor B and T cells with a test error of only 13% using two CpG positions. This low-dimensional classification was based on methylation status of CpG sites from CSNK2B and CDC25A, both of which were hypermethylated in ALL patients. The test error could be further improved to 5% by increasing the number of CpG positions used for classification to a total of 31.

In a next step, we tried to distinguish the two classes of leukaemia. Again the support vector machine was initially presented with a training set of AML and ALL samples (from both patient samples and cell lines) with the class information attached. Previously not included samples could then be classified into one of the classes with a test error of 5% when using two CpG positions. The optimal number of CpG sites in this case was calculated to be 6, decreasing the test error to only 1% (Table 1 and Fig 3B). The most informative CpG sites were located in intron 1 of the CDK4 gene and the promoter region of CSNK2B. The overall degree of methylation at these sites was reproducibly higher in ALL than in AML cells. However, the CpG sites in the CDK4 gene contributing information to this particular decision were different from those in the same gene distinguishing between ALL and healthy lymphocytes. This shows clearly that in some cases, different CpG sites within one cluster can contribute independent information, with different CpG sites potentially answering questions on different aspects of a phenotype.

There is evidence that CpG islands of several tissue-specific genes become methylated in cell lines (F. Antequera, J. Boyes, A. Bird, *Cell* 62, 503 (1990)). Hence, it is important to rule out that the above classifications are based on cell line specific artefacts rather than on disease specific changes. We therefore grouped leukaemia cell lines in one and primary leukaemia samples into another class. Classification was possible with at test error of 28% based on 2 CpG positions and 12% using 30 CpG sites (Table 1). Importantly, the CpG sites that contributed most to this distinction were not informative in any of the earlier classifications, excluding that any of them is based merely on cell line specific markers. This knowledge then enabled us to identify disease specific methylation differences between healthy T cells and T-ALL cell lines. Classification was performed with a test error of 0% based on two CpG positions, which could not distinguish cell lines from non cell lines. Decisive CpG sites were located in the coding sequence of the c-MOS gene (Table 1). In conclusion, although our results confirm the occurrence of cell line specific methylation changes, these do not interfere with classification of tumour types.

Figure 7 shows the use of a selection of genes from the panel for use in differentiating between acute lymphocytic leukemia and acute myelogenous leukemia. Using the 2 most significant CpG positions, test error was calculated to be 5%. The optimal number of CpG sites in this case was calculated to be 6, decreasing the test error to 1%. The most informative CpG sites were located in intron 1 of the CDK4 gene and the promoter region of CSNK2B. The overall degree of methylation at these sites was reproducibly higher in ALL than in AML cells. However, the CpG sites in the CDK4 gene contributing information to this particular decision were different from those in the same gene distinguishing between ALL and healthy lymphocytes. This shows clearly that in some cases, different CpG sites within one cluster can contribute independent information, with different CpG sites potentially answering questions on different aspects of a phenotype.

Figure 8 shows the use of an alternative selection of genes from the panel for use in differentiating between acute lymphocytic leukemia and acute myelogenous leukemia. The test error in this case was 16.3%.

In all classes classification results were comparable to mRNA based assays regarding test error, safety and reproducibility. However, classification was achieved with significantly fewer, randomly chosen genes, suggesting that methylation data may actually contain more information per dimension than expression profiles.

More importantly, in expression profiling signal intensities strongly depend on both the absolute and relative amounts of the different mRNA-species, comparison between independent experiments is challenging at best. Hence, although signals are generally calibrated against those derived from housekeeping genes, there remains an inability of mRNA profiling to detect subtle changes in expression levels (R. J. Lipshutz, S. P. A Fodor, T. R. Gingeras, D. J. Lockhart, *Nature Genetics* 21, 20 (1999)). In our approach, a CG/TG ratio is used as an internal calibration, and thus the amount of probe hybridised to the slide does not influence the results. This greatly improves the comparability of the results and therefore enables the screening of larger populations, as for example needed in multicenter trials and prospective studies.

Example 3: Identification of the methylation status of a CpG site within the gene N33 (NM 006765).

A fragment of the gene N33 (Seq ID NO: &[GENEID_2188]) was PCR amplified using primers AAAGCCGCTGCCATCC and TTTCGGCGACGGTAGG. The resultant fragment (548 bp in length) contained an informative CpG at position 446. The amplificate DNA was digested with the restriction endonuclease *FauI*, recognition site CCCGC. Hydrolysis by said endonuclease is blocked by methylation of the CpG at position 446 of the amplificate. The digest was used as a control.

Genomic DNA was isolated from sample using the DNA wizzard DNA isolation kit (Promega). Each sample was digested using *AvaI* according to manufacturer's recommendations (New England Biolabs).

10 ng of each genomic digest was then amplified using PCR primers AAAGCCGCTGCCATCC and TTTCGGCGACGGTAGG. The PCR reactions were performed using a thermocycler (Eppendorf GmbH) using 10 ng of DNA, 6 pmole of each primer, 200 μM of each dNTP, 1.5 mM MgCl2 and 1 U of HotstartTaq (Qiagen AG). The other conditions were as recommended by the Taq polymerase manufacturer. Using the above mentioned primers, gene fragments were amplified by PCR performing a first denaturation step for 14 min at 96 °C, followed by 30 - 45 cycles (step 2: 60 sec at 96°C, step 3: 45 sec at 52 °C, step 4: 75 sec at 72 °C) and a subsequent final elongation of 10 min at 72 °C. The presence of PCR products was analysed by agrarose gel electrophoresis.

PCR products were detectable with Faul hydrolyzed DNA isolated wherein the CpG position in question was upmethylated, when step 2 to step 4 of the cycle program were repeated 34, 37, 39, 42 and 45 fold. In contrast PCR products were only detectable with Faul hydrolyzed DNA isolated from downmethylated DNA (and control DNA) when step 2 to step 4 of the cycle program were repeated 39, 42 and 45 fold. These results were incorporated into a CpG methylation matrix analysis as described in Example 2.

<u>Tables</u>

Table 1

Class 1		Class 2		Test Error 2 CpGs	Significance 2 CpGs	Optimal CpG Number	Test Error optimal CpGs	Significance optimal CpGs
Healthy-B Healthy-T	13	B-ALL-cl T-ALL-cl ALL	17	15%	1.3E-4	54	4%	3.2E-7
B-ALL-ci T-ALL-ci ALL	17	AML AML-cl	8	5%	3.7E-4	6	1%	5.4E-10
Female-nocl	14	Male-nocl	7	5%	7.4E-4	2	5%	7.4E-4
Healthy-T	6	T-ALL-cl	6	0%	2.5E-3	2	0%	2.5E-3
Healthy-B Healthy-T	13	AML AML-ci	8	16%	6.9E-3	22	0%	2.3E-6
Primary leukaemia samples	8	CellLines	17	28%	2.2E-2	30	12%	7.9E-5
Healthy-B Healthy-T	13	ALL	5	13%	8.2E-3	31	5%	9.1E-5
Healthy-B	7	B-ALL-cl	6	27%	5.5E-2	11	8%	1.8E-3
Female-cl	9	Male-cl	8	52%	5.3E-1	42	40%	1.9E-1
B-ALL-cl	6	T-ALL-cl	6	58%	5.7E-1	70	51%	3.9E-1
Healthy-B	7	Healthy-T	6	63%	6.8E-1	11	49%	3.8E-1

Table 2: PCR primers and products

No:	Gene:	Primer:	Amplificate
		·	Length:
1	MDR1	TAAGTATGTTGAAGAAAGATTATTGTAG	633
	(SEQ ID NO: 37)	(SEQ ID NO: 388)	
		TAAAAACTATCCCATAATAACTCCCAAC	
<u>.</u>		(SEQ ID NO: 387)	
2	CSNK2B	GGGGAAATGGAGAAGTGTAA	524
	(SEQ ID NO: 67)	(SEQ ID NO: 512)	
		CTACCAATCCCAAAATAACC	
		(SEQ ID NO: 511)	
3	EGR4	AGGGGGATTGAGTGTTAAGT	293
	(SEQ ID NO: 26)	(SEQ ID NO: 392)	
		CCCAAACATAAACACAAAAT	
		(SEQ ID NO: 391)	
4	AR	GTAGTAGTAGTAAGAGA	460
	(SEQ ID NO: 69)	(SEQ ID NO: 532)	
		ACCCCCTAAATAATTATCCT	
		(SEQ ID NO: 531)	

ì

No:	Gene:	Primer:	Amplificate
			Length:
5	CDK 4	GGTAGTTGGTTATATGGTGAGG	748
	(SEQ ID NO: 70)	(SEQ ID NO: 395)	
		TCACACTCTTAAAAACCACAAAA	
		(SEQ ID NO: 396)	
6	Humos	TGATTGGGAGTAGGTGTTT	523
	(SEQ ID NO: 71)	(SEQ ID NO: 397)	
		CAAATCTTCCAACTTCTCAAA	
		(SEQ ID NO: 398)	
7	RB1	TTTAAGTTTGTTTTTGTTTTGGT	719
	(SEQ ID NO: 52)	(SEQ ID NO: 399)	
		TCCTACTCTAAATCCTCCTCAA	
		(SEQ ID NO: 400)	
8	CDC25A	TTGGGAGTTTTTATTGATTTTT	445
	(SEQ ID NO: 72)	(SEQ ID NO: 401)	
		ACAACCTAAAAATTAAATCCAAA	
		(SEQ ID NO: 402)	·
9	GPIb beta	GGTGATAGGAGAATAATGTTGG	379
	(SEQ ID NO: 30)	(SEQ ID NO: 403)	
		TCTCCCAACTACAACCAAAC	
		(SEQ ID NO: 404)	
10	MYOD1	ATTAGGGGTATAGAGGAGTATTGA	883
	(SEQ ID NO: 45)	(SEQ ID NO: 405)	
		CTTACAAACCCACAATAAACAA	
		(SEQ ID NO: 406)	
11	CDH3	GTTTAGAAGTTTAAGATTAG	611
	(SEQ ID NO: 17)	(SEQ ID NO: 407)	
		CAAAAACTCAACCTCTATCT	
		(SEQ ID NO: 408)	
12	WT1	AAAGGAAATTAAGTGTTGT	747
	(SEQ ID NO: 59)	(SEQ ID NO: 410)	
	}	TAACTACCCTCAACTTCCC	
		(SEQ ID NO: 409)	

No:	Gene:	Primer:	Amplificate
			Length:
13	MYCL1	AGGTTTGGGTTATTGAGTTT	491
	(SEQ ID NO: 44)	(SEQ ID NO: 411)	
		CATTATTTCCTAACTACCTTATATCTC	
		(SEQ ID NO: 412)	
14	ABL1	GGTTGGGAGATTTAATTTATT	967
	(SEQ ID NO: 2)	(SEQ ID NO: 414)	
		ACCAATCCAAACTTTTCCTT	
		(SEQ ID NO: 413)	•
15	ELK1	AAGTGTTTTAGTTTTTAATGGGTA	966
	(SEQ ID NO: 27)	(SEQ ID NO: 415)	
		CAAACCCAAAACTCACCTAT	
		(SEQ ID NO: 416)	
16	ABL1	GTTAGGAGGGGTTAAGG	291
	(SEQ ID NO: 1)	(SEQ ID NO: 417)	
	.	CCAACTTCAAACAAATCTCC	
		(SEQ ID NO: 418)	
. 17	APC	TCAACTACCATCAACTTCCTTA	
	(SEQ ID NO: 4)	(SEQ ID NO: 419)	
	,	AATTTATTTTAGTGTTGTAGTGGG	
		(SEQ ID NO: 420)	
18	ARHI	GTGAGTTTTTGGGGTGTTTA	444
	(SEQ ID NO: 6)	(SEQ ID NO: 421)	
		TCAATCTTACTTTCACACTACATAA	
		(SEQ ID NO: 422)	
19	BCL2	GTATTTATGTTAAGGGGGAAA	640
	(SEQ ID NO: 9)	(SEQ ID NO: 423)	
		AAAAACCACAATCCTCCC	
		(SEQ ID NO: 424)	
20	CCND2	TTTTGGTATGTAGGTTGGATG	426
	(SEQ ID NO: 13)	(SEQ ID NO: 425)	
		CCTAACCTCCTTCCTTTAACT	
		(SEQ ID NO: 426)	

PCT/EP02/03401

No:	Gene:	Primer:	Amplificate
21	CDH1	CAAATAAACCCTCAACCAATC	Length:
	(SEQ ID NO: 16)	(SEQ ID NO: 427)	,,,
		TGGAGGGGTAGGAAAGT	
		(SEQ ID NO: 428)	
22	CDKN1A	GGATTAGTGGGAATAGAGGTG	408
	(SEQ ID NO: 19)	(SEQ ID NO: 429)	,
		AAACCCAAACTCCTAACTACC	
		(SEQ ID NO: 430)	
23	CDKN1B (p27	GTGGGGAGGTAGTTGAAGA	478
,	Kip1)	(SEQ ID NO: 431)	
	(SEQ ID NO: 20)	ATACACCCTAACCCAAAAT	
		(SEQ ID NO: 432)	
24	CDKN2a	TTGAAAATTAAGGGTTGAGG	598
	(SEQ ID NO: 22)	(SEQ ID NO: 433)	
		CACCCTCTAATAACCAACCA	
		(SEQ ID NO: 434)	
25	CDKN2a	GGGGTTGGTTATTAGA	256
	(SEQ ID NO: 22)	(SEQ ID NO: 435)	
		AACCCTCTACCCACCTAAAT	
		(SEQ ID NO: 436)	
26	CDKN2B	GGTTGGTTGAAGGAATAGAAAT	708
	(SEQ ID NO: 23)	(SEQ ID NO: 437)	
		CCCACTAAACATACCCTTATTC	
		(SEQ ID NO: 438)	
27	DAPK1	AACCCTTTCTTCAAATTACAAA	348
	(SEQ ID NO: 25)	(SEQ ID NO: 439)	
		TGATTGGGTTTTAGGGAAATA	
		(SEQ ID NO: 440)	
28	FOS	TTTTTGGGGTTTAGTTTAGAAT	308
	(SEQ ID NO: 29)	(SEQ ID NO: 441)	
		AACCTTCATCCCCTAACCT	
		(SEQ ID NO: 442)	. ,

No:	Gene:	Primer:	Amplificate
			Length:
29	GSTP1	ATTTGGGAAAGAGGGAAAG	300
	(SEQ ID NO: 33)	(SEQ ID NO: 443)	
		TAAAAACTCTAAACCCCATCC	
		(SEQ ID NO: 444)	
30	HIC-1	TGGGTTGGAGAAGAGTTTA	280
	(SEQ ID NO: 34)	(SEQ ID NO: 445)	
		TCATATTTCCAAAAACACACC	
		(SEQ ID NO: 446)	
31	IGF2	CCCTTCCCCTTAACTAAACT	364
	(SEQ ID NO: 36)	(SEQ ID NO: 447)	
		AATTTGGGTTAGGTTTGGA	
		(SEQ ID NO: 448)	
32	MGMT	AAGGTTTTAGGGAAGAGTGTTT	636
	(SEQ ID NO: 38)	(SEQ ID NO: 449)	
		ACCTTTTCCTATCACAAAAATAA	
		(SEQ ID NO: 450)	
33	MLH1	TAAGGGAGAGGAGTTT	545
	(SEQ ID NO: 39)	(SEQ ID NO: 451)	
		ACCAATTCTCAATCATCTCTTT	
		(SEQ ID NO: 452)	·
34	MOS	ACCCTACAACAATCCCTCA	343
	(SEQ ID NO: 40)	(SEQ ID NO: 453)	
		TGGTTTTTAGGTTATTGGATTT	
		(SEQ ID NO: 454)	
35	MPL	TGGGGAGATTGATTTGAGTAT	594
	(SEQ ID NO: 42)	(SEQ ID NO: 455)	
		AATTCCTTCAATAACATACCCTT	
		(SEQ ID NO: 456)	
36	MYC	AGAGGAGTAAAAGAAAATGGT	712
	(SEQ ID NO: 43)	(SEQ ID NO: 457)	
;		CCAAATAAACAAAATAACCTCC	
		(SEQ ID NO: 458)	

No:	Gene:	Primer:	Amplificate
			Length:
37	N33	TTTTAGATTGAGGTTTTAGGGT	498
	(SEQ ID NO: 46)	(SEQ ID NO: 459)	
		ATCCATTCTACCTCCTTTTTCT	
		(SEQ ID NO: 460)	
38	PMS2	GATTGAGATTATTTTGGGTTTT	561
	(SEQ ID NO: 49)	(SEQ ID NO: 461)	
		ACTTAAACCTTCCCTCTCCAC	
		(SEQ ID NO: 462)	
39	PTEN	TTTTAGGTAGTTATATTGGGTATGTT	346
	(SEQ ID NO: 51)	(SEQ ID NO: 463)	
		TCAACTCTCAAACTTCCATCA	
		(SEQ ID NO: 464)	
40	RBL2	GAAAATGGGTGTGTGG	112
	(SEQ ID NO: 53)	(SEQ ID NO: 465)	
		TACAAATAAAAACAAATCCCCT	
		(SEQ ID NO: 466)	
41	SFN	GAAGAGAGAGGAGGTA	489
	(SEQ ID NO: 55)	(SEQ ID NO: 468)	
		CTATCCAACAACCCAACA	
		(SEQ ID NO: 467)	
42	TGFBR2	GTAATTTGAAGAAAGTTGAGGG	296
	(SEQ ID NO: 57)	(SEQ ID NO: 469)	
		CCAACAACTAAACAAAACCTCT	
<u>.</u>		(SEQ ID NO: 470)	
43	TP73	AGTAAATAGTGGGTGAGTTATGAA	607
	(SEQ ID NO: 58)	(SEQ ID NO: 472)	
:		GAAAAACCTCTAAAAACTACTCTCC	
		(SEQ ID NO: 471)	
44	CDKN1C	GGGGAGGTAGATATTTGGATAA	300
	(SEQ ID NO: 21)	(SEQ ID NO: 473)	
		AACTACACCATTTATATTCCCAC	
		(SEQ ID NO: 474)	

No:	Gene:	Primer:	Amplificate
			Length:
45	GSK3ß	TAAGTGATAAAGGAAGGA	243
	(SEQ ID NO: 32)	(SEQ ID NO: 475)	
		CCTTCAAACCCCAAACAA	
		(SEQ ID NO: 476)	_
46	CDC2	ATTAGAAGTGAAAGTAATGGAATTT	418
	(SEQ ID NO: 14)	(SEQ ID NO: 477)	
		TCAATTTCCAAAAACCAAC	
		(SEQ ID NO: 478)	
47	ESR1	AGGGGAATTAAATAGAAAGAG	662
	(SEQ ID NO: 28)	(SEQ ID NO: 479)	
		CAATAAAACCATCCCAAATACT	
		(SEQ ID NO: 480)	•
48	APAF1	AGATATGTTTGGAGATTTTAGGA	674
	(SEQ ID NO: 3)	(SEQ ID NO: 481)	
		AACTCCCCACCTCTAATTCTAT	
		(SEQ ID NO: 482)	
49	BAK1	AATTAGGGATGGGAAAAGTAGT	558
	(SEQ ID NO: 7)	(SEQ ID NO: 483)	
		AAACATAACAAAATCAAATCCC	
		(SEQ ID NO: 484)	
50	BAX	AAATAAATAGAAAAGTAGGTTTGGC	716
	(SEQ ID NO: 8)	(SEQ ID NO: 485)	
•		TTCTACCCCTCAATACTTAAAAA	•
		(SEQ ID NO: 486)	
51	CASP10	TCTTCCCAAACAAATAAAAACT	468
	(SEQ ID NO: 10)	(SEQ ID NO: 487)	
		GTTTTGAGGTTAAATGAGTGGT	
		(SEQ ID NO: 488)	
52	CASP8	AGTGGATTTGGAGTTTAGATGT	431
	(SEQ ID NO: 11)	(SEQ ID NO: 489)	
		AACAAAATAAAAACTTCTCCCA	
		(SEQ ID NO: 490)	

No:	Gene:	Primer:	Amplificate
	}		Length:
53	CASP9	GGAAGAGTTGTAGGTGGATTAG	424
	(SEQ ID NO: 12)	(SEQ ID NO: 491)	
		CCACATTTTCCCAATAAATACT	·
		(SEQ ID NO: 492)	
54	TCL1A	CCACAACACATCAAACCTATAA	491
	(SEQ ID NO: 56)	(SEQ ID NO: 493)	
		TTTGGTGGTAAGTGAGGGT	
		(SEQ ID NO: 494)	
55	PML	AATTCCCAAACAAACTTTTACA	459
	(SEQ ID NO: 48)	(SEQ ID NO: 495)	
		TGGAGAGGAAGTGAGATAGAAG	
		(SEQ ID NO: 496)	
56	HOXA5	AAACCCCAAACAACCTCTAT	392
	(SEQ ID NO: 35)	(SEQ ID NO: 498)	
		GAAGGGGAAAGTTATTTAGTTA	
	·	(SEQ ID NO: 497)	
57	SDC4	CCTAACTACCCTCATTCCTTT	269
	(SEQ ID NO: 54)	(SEQ ID NO: 499)	
		AGTTGGGGAAATTAAGGTTTAG	
		(SEQ ID NO: 500)	
58	PITX2	TCCTCAACTCTACAAACCTAAAA	408
	(SEQ ID NO: 47)	(SEQ ID NO: 501)	
		GTAGGGGAGGGAAGTAGATGT	
		(SEQ ID NO: 502)	
59	GPR37	ACTTATTTTCTTTTCCTCTAAAAAC	489
	(SEQ ID NO: 31)	(SEQ ID NO: 503)	
		TATGGTTTGGTGAGGGTATATT	
		(SEQ ID NO: 504)	
60	PRAME	ACCACTCAAAAACAAACCTTTA	303
	(SEQ ID NO: 50)	(SEQ ID NO: 505)	
		ATGGGTATGGTGGTTAAGAGTA	
		(SEQ ID NO: 506)	

No:	Gene:	Primer:	Amplificate
			Length:
61	CDK 4	TTTTGGTAGTTGGTTATATG	474
	(SEQ ID NO: 70)	(SEQ ID NO: 508)	
		AAAAATAACACAATAACTCA	
		(SEQ ID NO: 507)	
62	CDC25A	AGAAGTTGTTTATTGATTGG	272
	(SEQ ID NO: 72)	(SEQ ID NO: 510)	
		AAAATTAAATCCAAACAAAC	
		(SEQ ID NO: 509)	
63	CSNK2B	GGGAAATGGAGAAGTGTAA	524
	(SEQ ID NO: 67)	(SEQ ID NO: 512)	
	1	CTACCAATCCCAAAATAACC	
		(SEQ ID NO: 511)	
64	C-ABL	ATTTATTTTG	568
	(SEQ ID NO: 62)	(SEQ ID NO: 514)	
		ATCCCAAACTAAAAAAAAACAAACCACC	
		(SEQ ID NO: 513)	
65	L-MYC	AAACCCAAAAACTAAAAAA	461
	(SEQ ID NO: 61)	(SEQ ID NO: 515)	-
		ATGTTATTTAGGTTTGGGTTATTG	
		(SEQ ID NO: 516)	
66	CSF1	GAGAATGGGGTGG	
	(SEQ ID NO: 65)	(SEQ ID NO: 518)	
		TATAAAAAATCACCCTAACC	
		(SEQ ID NO: 517)	
67	Me491/TD63	ACCTAAATCA	
	(SEQ ID NO: 68)	(SEQ ID NO: 519)	
		ATGTGAGTTAGGGATTTTTT	
-25-		(SEQ ID NO: 520)	
68	CD1R3	ATTATGGTTGGAATTGTAAT	413
	(SEQ ID NO: 66)	(SEQ ID NO: 522)	
		ACAAAACAACAACACCCC	
		(SEQ ID NO: 521)	

No:	Gene:	Primer:	Amplificate
			Length:
69	N-MYC	GGGAGGAGTATATTTTG	368
	(SEQ ID NO: 60)	(SEQ ID NO: 524)	
		AAAATATCCTC	
		(SEQ ID NO: 523)	
70	Tubulin	CTTATATACCCTCCC	453
	(SEQ ID NO: 64)	(SEQ ID NO: 525)	
		GGGTTTAGAGTTAAGATTG	·
		(SEQ ID NO: 526)	
71	CMYCex3	AGGAGGAATAAGAAGATGAG	
	(SEQ ID NO: 73)	(SEQ ID NO: 528)	
		TTCAATCTCAAAACTCAACC	
		(SEQ ID NO: 527)	
72	ELK1	GAAAAAAAAACCCAATAT	523
	(SEQ ID NO: 63)	(SEQ ID NO: 529)	
		ATGGTTTTGTTTAAT	
		(SEQ ID NO: 530)	
73	AR	GTAGTAGTAGTAAGAGA	460
	(SEQ ID NO: 69)	(SEQ ID NO: 532)	
		ACCCCTAAATAATTATCCT	
		(SEQ ID NO: 531)	
74	Humos	TTTATTGATTGGGAGTAGGT	494
	(SEQ ID NO: 71)	(SEQ ID NO: 534)	
		CTAATTTTACAAACATCCTA	
		(SEQ ID NO: 533)	

Table 3: Hybridisation oligonucleotides

No:	Gene	Oligo:
1	MDR1	TTGGTGGTCGTTTTAAGG
	(SEQ ID NO: 37)	(SEQ ID NO: 535)
2	MDR1	TTGGTGGTTGTTTAAGG
	(SEQ ID NO: 37)	(SEQ ID NO: 536)

No:	Gene	Oligo:
3	MDR1	TTGAAAGACGTGTTTATA
	(SEQ ID NO: 37)	(SEQ ID NO: 537)
4	MDR1	TTGAAAGATGTGTTTATA
	(SEQ ID NO: 37)	(SEQ ID NO: 538)
5	MDR1	AGGTGTAACGGAAGTTAG
	(SEQ ID NO: 37)	(SEQ ID NO: 539)
6	MDR1	AGGTGTAATGGAAGTTAG
	(SEQ ID NO: 37)	(SEQ ID NO: 540)
7	MDR1	TAGTTTTTCGAGGAATTA
	(SEQ ID NO: 37)	(SEQ ID NO: 541)
8	MDR1	TAGTTTTTGAGGAATTA
	(SEQ ID NO: 37)	(SEQ ID NO: 542)
9	CSNK2B	AGGAGTTTCGGAGGAAAT
	(SEQ ID NO: 67)	(SEQ ID NO: 1235)
10	CSNK2B	AGGAGTTTTGGAGGAAAT
	(SEQ ID NO: 67)	
11	CSNK2B	(SEQ ID NO: 1236) GAGAGTTGCGGAAAGAGA
	(SEQ ID NO: 67)	
12	CSNK2B	(SEQ ID NO: 543) GAGAGTTGTGGAAAGAGA
13	(SEQ ID NO: 67) CSNK2B	(SEQ ID NO: 544) GGGTTTTCGTGATAGT
14	(SEQ ID NO: 67) CSNK2B	(SEQ ID NO: 545) GGGTTTTTGTGATAGT
15	(SEQ ID NO: 67) CSNK2B	(SEQ ID NO: 546) TAGGTTAGCGTATTGGGA
15		
16	(SEQ ID NO: 67) CSNK2B	(SEQ ID NO: 1117) TAGGTTAGTGTATTGGGA
10		
17	(SEQ ID NO: 67) EGR4	(SEQ ID NO: 1118) GGTGGGAAGCGTATTTAT
11		
18	(SEQ ID NO: 26) EGR4	(SEQ ID NO: 1213)
10		GGTGGGAAGTGTATTTAT
19	(SEQ ID NO: 26) EGR4	(SEQ ID NO: 1214)
19		AATAATAACGTTATAGTT
20	(SEQ ID NO: 26) EGR4	(SEQ ID NO: 549)
40		AATAATATGTTATAGTT
21	(SEQ ID NO: 26) EGR4	(SEQ ID NO: 550)
21		TTATAGTTCGAGTTTTTT
	(SEQ ID NO: 26)	(SEQ ID NO: 1225)

No:	Gene	Oligo:
22	EGR4	TTATAGTTTGAGTTTTT
	(SEQ ID NO: 26)	(SEQ ID NO: 1226)
23	EGR4	GGAGTTTTCGGTATATAT
	(SEQ ID NO: 26)	(SEQ ID NO: 551)
24	EGR4	GGAGTTTTTGGTATATAT
	(SEQ ID NO: 26)	(SEQ ID NO: 552)
25	AR	TGTTATTTCGAGAGAGGT
	(SEQ ID NO: 69)	(SEQ ID NO: 553)
26	AR	TGTTATTTTGAGAGAGGT
1	(SEQ ID NO: 69)	(SEQ ID NO: 554)
27	AR	AGAGGTTGCGTTTTAGAG
	(SEQ ID NO: 69)	(SEQ ID NO: 555)
28	AR	AGAGGTTGTGTTTTAGAG
	(SEQ ID NO: 69)	(SEQ ID NO: 556)
29	AR	GTAGTATTCGAAGGTAGT
	(SEQ ID NO: 69)	(SEQ ID NO: 1255)
30	AR	GTAGTATTTGAAGGTAGT
	(SEQ ID NO: 69)	(SEQ ID NO: 1256)
31	AR	GGAGGTTTCGGGGGTTTT
	(SEQ ID NO: 69)	(SEQ ID NO: 557)
32	AR	GGAGGTTTTGGGGGTTTT
	(SEQ ID NO: 69)	(SEQ ID NO: 558)
33	CDK 4	GTATGGGGTCGTAGGAAT
	(SEQ ID NO: 70)	(SEQ ID NO: 559)
34	CDK 4	GTATGGGGTTGTAGGAAT
	(SEQ ID NO: 70)	(SEQ ID NO: 560)
35	CDK 4	GGAAGGGTCGTTTAAGGG
	(SEQ ID NO: 70)	(SEQ ID NO: 1147)
36	CDK 4	GGAAGGGTTGTTTAAGGG
	(SEQ ID NO: 70)	(SEQ ID NO: 1148)
37	CDK 4	GGGTTGGCGTGAGGTA
	(SEQ ID NO: 70)	(SEQ ID NO: 563)
38	CDK 4	GGGTTGGTGAGGTA
	(SEQ ID NO: 70)	(SEQ ID NO: 564)
39 ·	CDK 4	AGGATTTTCGATGTAAGG
	(SEQ ID NO: 70)	(SEQ ID NO: 565)
40	CDK 4	AGGATTTTTGATGTAAGG
	(SEQ ID NO: 70)	(SEQ ID NO: 566)

77	<i>a</i>	OI!
No: 41	Gene CDK 4	Oligo: GGGTTTTACGTGGTTGGA
41		
40	(SEQ ID NO: 70) CDK 4	(SEQ ID NO: 567)
42		GGGTTTTATGTGGTTGGA
4.0	(SEQ ID NO: 70)	(SEQ ID NO: 568)
43	Humos	GAGTTTAACGTAGTAAGG
	(SEQ ID NO: 71)	(SEQ ID NO: 1221)
44	Humos	GAGTTTAATGTAGTAAGG
	(SEQ ID NO: 71)	(SEQ ID NO: 1222)
45	Humos	TATGGAGTTCGGTGGTAA
	(SEQ ID NO: 71)	(SEQ ID NO: 569)
46	Humos	TATGGAGTTTGGTGGTAA
	(SEQ ID NO: 71)	(SEQ ID NO: 570)
47	Humos	TTTATTGTCGTATTGGAG
	(SEQ ID NO: 71)	(SEQ ID NO: 571)
48	Humos	TTTATTGTTGTATTGGAG
	(SEQ ID NO: 71)	(SEQ ID NO: 572)
49	Humos	GTTGTGAACGGTTTGTTT
	(SEQ ID NO: 71)	(SEQ ID NO: 573)
50	Humos	GTTGTGAATGGTTTGTTT
	(SEQ ID NO: 71)	(SEQ ID NO: 574)
51	RB1	TTAGATTTCGGGATAGGG
	(SEQ ID NO: 52)	(SEQ ID NO: 575)
52	RB1	TTAGATTTTGGGATAGGG
	(SEQ ID NO: 52)	(SEQ ID NO: 576)
53	RB1	TATAGTTTCGTTAAGTGT
	(SEQ ID NO: 52)	(SEQ ID NO: 577)
54	RB1	TATÀGTTTTGTTAAGTGT
	(SEQ ID NO: 52)	(SEQ ID NO: 578)
55	RB1	GTGTATTTCGGTTTGGAG
	(SEQ ID NO: 52)	(SEQ ID NO: 579)
56	RB1	GTGTATTTTGGTTTGGAG
	(SEQ ID NO: 52)	(SEQ ID NO: 580)
57	RB1	TTGGAAGGCGTTTGGATT
	(SEQ ID NO: 52)	(SEQ ID NO: 581)
58	RB1	TTGGAAGGTGTTTGGATT
	(SEQ ID NO: 52)	(SEQ ID NO: 582)
59	RB1	TGGATTTACGTTAGGTTT
	(SEQ ID NO: 52)	(SEQ ID NO: 583)

No:	Gene	Oligo:
60	RB1	TGGATTTATGTTAGGTTT
	(SEQ ID NO: 52)	(SEQ ID NO: 584)
61	CDC25A	GTGTAGGTCGGTTTGGTT
	(SEQ ID NO: 72)	(SEQ ID NO: 1159)
62	CDC25A	GTGTAGGTTGGTT
	(SEQ ID NO: 72)	(SEQ ID NO: 1160)
63	CDC25A	TTGTTATTCGGAGTTGGG
	(SEQ ID NO: 72)	(SEQ ID NO: 1163)
64	CDC25A	TTGTTATTTGGAGTTGGG
	(SEQ ID NO: 72)	(SEQ ID NO: 1164)
65	CDC25A	GGAGAATAGCGAAGATAG
	(SEQ ID NO: 72)	(SEQ ID NO: 589)
66	CDC25A	GGAGAATAGTGAAGATAG
	(SEQ ID NO: 72)	
67	CDC25A	(SEQ ID NO: 590) GAAAGGTCGGTTTGGT
68	(SEQ ID NO: 72) CDC25A	(SEQ ID NO: 591) GAAAGGTTGGTT
69	(SEQ ID NO: 72) GPIb beta	(SEQ ID NO: 592) TTTGAGAGCGGGTGGGAG
70	(SEQ ID NO: 30) GPIb beta	(SEQ ID NO: 593) TTTGAGAGTGGGTGGGAG
1	·	
71	(SEQ ID NO: 30) GPIb beta	(SEQ ID NO: 594) GTGGGAGCGGAAGTTTGA
'1		
72	(SEQ ID NO: 30) GPIb beta	(SEQ ID NO: 595)
1 /2		GTGGGAGTGGAAGTTTGA
72	(SEQ ID NO: 30) GPIb beta	(SEQ ID NO: 596)
73		GGTTAGGTCGTAGTATTG
71	(SEQ ID NO: 30)	(SEQ ID NO: 597)
74	GPIb beta	GGTTAGGTTGTAGTATTG
7.	(SEQ ID NO: 30)	(SEQ ID NO: 598)
75	GPIb beta	ATGGGTTTCGGTGAGTTT
	(SEQ ID NO: 30)	(SEQ ID NO: 599)
76	GPIb beta	ATGGGTTTTGGTGAGTTT
	(SEQ ID NO: 30)	(SEQ ID NO: 600)
77	GPIb beta	GGGTTTTCGGTTGTTTT
	(SEQ ID NO: 30)	(SEQ ID NO: 601)
78	GPIb beta	GGGTTTTTGGTTGTTTTT
	(SEQ ID NO: 30)	(SEQ ID NO: 602)

No:	Gene	Oligo:
79	MYOD1	ATAGTAGTCGGGTGTTGG
ļ	(SEQ ID NO: 45)	(SEQ ID NO: 603)
80	MYOD1	ATAGTAGTTGGGTGTTGG
	(SEQ ID NO: 45)	(SEQ ID NO: 604)
81	MYOD1	GTGTTAGTCGTTTAGGGT
	(SEQ ID NO: 45)	(SEQ ID NO: 605)
82	MYOD1	GTGTTAGTTGTTTAGGGT
	(SEQ ID NO: 45)	(SEQ ID NO: 606)
83	MYOD1	TAGTTGTTCGTTTGGGTT
	(SEQ ID NO: 45)	(SEQ ID NO: 607)
84	MYOD1	TAGTTGTTTGGGTT
	(SEQ ID NO: 45)	(SEQ ID NO: 608)
85	MYOD1	GGTTATTACGGATAAATA
	(SEQ ID NO: 45)	(SEQ ID NO: 609)
86	MYOD1	GGTTATTATGGATAAATA
	(SEQ ID NO: 45)	
87	MYOD1	(SEQ ID NO: 610) AATTAGGTCGGATAGGAG
88	(SEQ ID NO: 45) MYOD1	(SEQ ID NO: 611) AATTAGGTTGGATAGGAG
V)	
89	(SEQ ID NO: 45) CDH3	(SEQ ID NO: 612) AAATTAGTCGGGTGTGGT
90	(SEQ ID NO: 17) CDH3	(SEQ ID NO: 613) AAATTAGTTGGGTGTGT
70		
91	(SEQ ID NO: 17) CDH3	(SEQ ID NO: 614) TGTGGTGGCGTAAGTTTG
71		
92	(SEQ ID NO: 17) CDH3	(SEQ ID NO: 615) TGTGGTGGTGTAAGTTTG
93	(SEQ ID NO: 17) CDH3	(SEQ ID NO: 616)
75		GGAGTTTTCGTTTTTAGT
94	(SEQ ID NO: 17) CDH3	(SEQ ID NO: 617)
7 1		GGAGTTTTTGTTTTTAGT
95	(SEQ ID NO: 17) CDH3	(SEQ ID NO: 618)
73		TAGAATTGCGAGATAGAG
96	(SEQ ID NO: 17) CDH3	(SEQ ID NO: 619)
<i>7</i> 0		TAGAATTGTGAGATAGAG
97	(SEQ ID NO: 17) WT1	(SEQ ID NO: 620)
71		TAGTGAGACGAGGTTTTT
	(SEQ ID NO: 59)	(SEQ ID NO: 621)

No:	Gene	Oligo:
98	WT1	TAGTGAGATGAGGTTTTT
ľ	(SEQ ID NO: 59)	(SEQ ID NO: 622)
99	WT1	TATATTGGCGAAGGTTAA
	(SEQ ID NO: 59)	(SEQ ID NO: 623)
100	WT1	TATATTGGTGAAGGTTAA
	(SEQ ID NO: 59)	(SEQ ID NO: 624)
101	WT1	TGTTATATCGGTTAGTTG
	(SEQ ID NO: 59)	(SEQ ID NO: 625)
102	WT1	TGTTATATTGGTTAGTTG
	(SEQ ID NO: 59)	(SEQ ID NO: 626)
103	WT1	TGTTTGGTCGGGTTTGGG
	(SEQ ID NO: 59)	(SEQ ID NO: 627)
104	WT1	TGTTTGGTTGGG
	(SEQ ID NO: 59)	(SEQ ID NO: 628)
105	WT1	TTTAGTTTCGATTTTTGG
	(SEQ ID NO: 59)	(SEQ ID NO: 629)
106	WT1	TTTAGTTTTGG
}	(SEQ ID NO: 59)	(SEQ ID NO: 630)
107	MYCL1	TTGAGGGTCGTTAGGTGG
į	(SEQ ID NO: 44)	(SEQ ID NO: 631)
108	MYCL1	TTGAGGGTTGTTAGGTGG
	(SEQ ID NO: 44)	(SEQ ID NO: 632)
109	MYCL1	TTTTAGTTCGGAGTGGGT
	(SEQ ID NO: 44)	(SEQ ID NO: 633)
110	MYCL1	TTTTAGTTTGGAGTGGGT
	(SEQ ID NO: 44)	(SEQ ID NO: 634)
111	MYCL1	AGTTTAGTCGGTTTGGTAT
	(SEQ ID NO: 44)	(SEQ ID NO: 635)
112	MYCL1	AGTTTAGTTGGTAT
	(SEQ ID NO: 44)	(SEQ ID NO: 636)
113	MYCL1	GGGGTTATCGGGGATTGA
	(SEQ ID NO: 44)	(SEQ ID NO: 637)
114	MYCL1	GGGGTTATTGGGGATTGA
	(SEQ ID NO: 44)	(SEQ ID NO: 638)
115	ABL1	GTTTGGGTCGCGAGAGTT
	(SEQ ID NO: 2)	(SEQ ID NO: 639)
116	ABL1	GTTTGGGTTGTGAGAGTT
	(SEQ ID NO: 2)	(SEQ ID NO: 640)

No:	Gene	Oligo:
117	ABL1	GGATATTGCGGGTGGTTT
1	(SEQ ID NO: 2)	(SEQ ID NO: 641)
118	ABL1	GGATATTGTGGGTGGTTT
	(SEQ ID NO: 2)	(SEQ ID NO: 642)
119	ABL1	TTTGGTTTCGTTGTGAG
	(SEQ ID NO: 2)	(SEQ ID NO: 643)
120	ABL1	TTTGGTTTTGTTGGAG
	(SEQ ID NO: 2)	(SEQ ID NO: 644)
121	ABL1	GGGTTTTGCGGTTGAGGA
	(SEQ ID NO: 2)	(SEQ ID NO: 645)
122	ABL1	GGGTTTTGTGGTTGAGGA
	(SEQ ID NO: 2)	(SEQ ID NO: 646)
123	ELK1	TTTGTTTTCGTTGAGTAG
	(SEQ ID NO: 27)	(SEQ ID NO: 647)
124	ELK1	TTTGTTTTGTTGAGTAG
	(SEQ ID NO: 27)	(SEQ ID NO: 648)
125	ELK1	TTTATTTTCGTTTTTGGG
	(SEQ ID NO: 27)	(SEQ ID NO: 649)
126	ELK1	TTTÄTTTTTTTTTGGG
	(SEQ ID NO: 27)	(SEQ ID NO: 650)
127	ELK1	GAAGGGTTCGTTTTTAA
	(SEQ ID NO: 27)	(SEQ ID NO: 651)
128	ELK1	GAAGGGTTTGTTTTTAA
	(SEQ ID NO: 27)	(SEQ ID NO: 652)
129	APC	TATTAGAGCGTTTTAAAG
	(SEQ ID NO: 4)	(SEQ ID NO: 653)
130	APC	TATTAGAGTGTTTTAAAG
	(SEQ ID NO: 4)	(SEQ ID NO: 654)
131	APC	GTTTTTTCGATTTGGGT
16-	(SEQ ID NO: 4)	(SEQ ID NO: 655)
132	APC	GTTTTTTTGATTTGGGT
	(SEQ ID NO: 4)	(SEQ ID NO: 656)
133	ARHI	TGTTGTTGCGTAGTAGAA
10.	(SEQ ID NO: 6)	(SEQ ID NO: 657)
134	ARHI	TGTTGTTGTGTAGAA
10.5	(SEQ ID NO: 6)	(SEQ ID NO: 658)
135	ARHI	GAATTATTCGTAGTTTTG
	(SEQ ID NO: 6)	(SEQ ID NO: 659)

No:	Gene	Oligo:
136	ARHI	GAATTATTTGTAGTTTTG
	(SEQ ID NO: 6)	(SEQ ID NO: 660)
137	ARHI	GATAGTTTCGTTGGGAAG
	(SEQ ID NO: 6)	(SEQ ID NO: 661)
138	ARHI	GATAGTTTTGTTGGGAAG
Ĭ	(SEQ ID NO: 6)	(SEQ ID NO: 662)
139	ARHI	TAGAAGAACGAGGTTTGA
	(SEQ ID NO: 6)	(SEQ ID NO: 663)
140	ARHI	TAGAAGAATGAGGTTTGA
	(SEQ ID NO: 6)	(SEQ ID NO: 664)
141	ARHI	TATAGTTGCGTAGGTAGT
	(SEQ ID NO: 6)	(SEQ ID NO: 665)
142	ARHI	TATAGTTGTGTAGT
	(SEQ ID NO: 6)	(SEQ ID NO: 666)
143	BCL2	AGAGGTGTCGTTGGTTTT
	(SEQ ID NO: 9)	(SEQ ID NO: 667)
144	BCL2	AGAGGTGTTGTTT
	(SEQ ID NO: 9)	(SEQ ID NO: 668)
145	BCL2	TAAGTTGTCGTAGAGGGG
	(SEQ ID NO: 9)	(SEQ ID NO: 669)
146	BCL2	TAAGTTGTTGTAGAGGGG
	(SEQ ID NO: 9)	(SEQ ID NO: 670)
147	BCL2	AGGGGTTACGAGTGGGAT
	(SEQ ID NO: 9)	(SEQ ID NO: 671)
148	BCL2	AGGGGTTATGAGTGGGAT
	(SEQ ID NO: 9)	(SEQ ID NO: 672)
149	BCL2	TTTTGTTACGGTGGTA
	(SEQ ID NO: 9)	(SEQ ID NO: 673)
150	BCL2	TTTTGTTATGGTGGTGGA
	(SEQ ID NO: 9)	(SEQ ID NO: 674)
151	CCND2	AGTTGGGTCGGGTTAGTT
	(SEQ ID NO: 13)	(SEQ ID NO: 675)
152	CCND2	AGTTGGGTTAGTT
	(SEQ ID NO: 13)	(SEQ ID NO: 676)
153	CCND2	TTTAATAACGAGAGGGGA
	(SEQ ID NO: 13)	(SEQ ID NO: 677)
154	CCND2	TTTAATAATGAGAGGGA
	(SEQ ID NO: 13)	(SEQ ID NO: 678)

No:	Gene	Oligo:
155	CCND2	TTAGTTTGCGTTATCGTT
	(SEQ ID NO: 13)	(SEQ ID NO: 679)
156	CCND2	TTAGTTTGTGTTATTGTT
1	(SEQ ID NO: 13)	(SEQ ID NO: 680)
157	CCND2	TTTTAGAGCGGAGAAGAG
	(SEQ ID NO: 13)	(SEQ ID NO: 681)
158	CCND2	TTTTAGAGTGGAGAAGAG
	(SEQ ID NO: 13)	(SEQ ID NO: 682)
159	CCND2	GGTAGTTTCGAGGTTTTG
	(SEQ ID NO: 13)	(SEQ ID NO: 683)
160	CCND2	GGTAGTTTTGAGGTTTTG
	(SEQ ID NO: 13)	(SEQ ID NO: 684)
161	CDH1	AGTTTCGACGTTATTGAG
	(SEQ ID NO: 16)	(SEQ ID NO: 685)
162	CDH1	AGTTTTGATGTTATTGAG
	(SEQ ID NO: 16)	
163	CDH1	(SEQ ID NO: 686) AGAGGTTGCGGTTTTAAG
	(SEQ ID NO: 16)	
164	CDH1	(SEQ ID NO: 687) AGAGGTTGTGGTTTTAAG
	(SEQ ID NO: 16)	
165	CDH1	(SEQ ID NO: 688) AGTAGCGTCGAGAGGTTG
	(SEQ ID NO: 16)	
166	CDH1	(SEQ ID NO: 689) AGTAGTGTTGAGAGGTTG
	(SEQ ID NO: 16)	
167	CDH1	(SEQ ID NO: 690) AGGGGATTCGGGGTATTT
	(SEQ ID NO: 16)	
168	CDH1	(SEQ ID NO: 691) AGGGGATTTGGGGTATTT
		·
169	(SEQ ID NO: 16) CDKN1A	(SEQ ID NO: 692) AGGTGTATCGTTTTATA
	(SEQ ID NO: 19)	
170	CDKN1A	(SEQ ID NO: 693) AGGTGTATTGTTTTATA
_, 5		•
171	(SEQ ID NO: 19) CDKN1A	(SEQ ID NO: 694) TGGGTTAGCGGTGAGTTA
- · -		
172	(SEQ ID NO: 19) CDKN1A	(SEQ ID NO: 695) TGGGTTAGTGGTGAGTTA
_ · - -		·
173	(SEQ ID NO: 19) CDKN1A	(SEQ ID NO: 696) GTTTATTTCGTGGGGAAA
5		
	(SEQ ID NO: 19)	(SEQ ID NO: 697)

No:	Gene	Oligo:
174	CDKN1A	GTTTATTTTGTGGGGAAA
	(SEQ ID NO: 19)	(SEQ ID NO: 698)
175	CDKN1A	TTGGAATTCGGTTAGGTT
	(SEQ ID NO: 19)	(SEQ ID NO: 699)
176	CDKN1A	TTGGAATTTGGTTAGGTT
	(SEQ ID NO: 19)	(SEQ ID NO: 700)
177	CDKN1B (p27	AAGAGAAACGTTGGAATA
	Kip1)	(SEQ ID NO: 701)
	(SEQ ID NO: 20)	
178	CDKN1B (p27	AAGAGAAATGTTGGAATA
	Kip1)	(SEQ ID NO: 702)
	(SEQ ID NO: 20)	•
179	CDKN1B (p27	TTTGATTTCGAGGGGAGT
	Kip1)	(SEQ ID NO: 703)
	(SEQ ID NO: 20)	
180	CDKN1B (p27	TTTGATTTTGAGGGGAGT
	Kip1)	(SEQ ID NO: 704)
	(SEQ ID NO: 20)	
181	CDKN1B (p27	GTATTTGGCGGTTGGATT
	Kip1)	(SEQ ID NO: 705)
	(SEQ ID NO: 20)	
182	CDKN1B (p27	GTATTTGGTGGTTT
	Kip1)	(SEQ ID NO: 706)
	(SEQ ID NO: 20)	
183	CDKN1B (p27	TATAATTTCGGGAAAGAA
	Kip1)	(SEQ ID NO: 707)
	(SEQ ID NO: 20)	
184	CDKN1B (p27	TATAATTTTGGGAAAGAA
	Kip1)	(SEQ ID NO: 708)
	(SEQ ID NO: 20)	
185	CDKN2a	AGAGTGAACGTATTTAAA
	(SEQ ID NO: 22)	(SEQ ID NO: 709)
186	CDKN2a	AGAGTGAATGTATTTAAA
	(SEQ ID NO: 22)	(SEQ ID NO: 710)
187	CDKN2a	GTTATATTCGTTAAGTGT
	(SEQ ID NO: 22)	(SEQ ID NO: 711)

No:	Gene	Oligo:
188	CDKN2a	GTTATATTTGTTAAGTGT
	(SEQ ID NO: 22)	(SEQ ID NO: 712)
189	CDKN2a	TAAGTGTTCGGAGTTAAT
	(SEQ ID NO: 22)	(SEQ ID NO: 713)
190	CDKN2a	TAAGTGTTTGGAGTTAAT
	(SEQ ID NO: 22)	(SEQ ID NO: 714)
191	CDKN2a	GTTAGTATCGGAGGAAGA
	(SEQ ID NO: 22)	(SEQ ID NO: 715)
192	CDKN2a	GTTAGTATTGGAGGAAGA
	(SEQ ID NO: 22)	(SEQ ID NO: 716)
193	CDKN2a	GGAGTTTTCGGTTGATTG
	(SEQ ID NO: 22)	(SEQ ID NO: 717)
194	CDKN2a	GGAGTTTTTGGTTGATTG
	(SEQ ID NO: 22)	(SEQ ID NO: 718)
195	CDKN2a	TTGTTTAACGTATCGAAT
	(SEQ ID NO: 22)	(SEQ ID NO: 719)
196	CDKN2a	TTGTTTAATGTATTGAAT
	(SEQ ID NO: 22)	
197	CDKN2a	(SEQ ID NO: 720) AATAGTTACGGTCGGAGG
	(SEQ ID NO: 22)	
198	CDKN2a	(SEQ ID NO: 721) AATAGTTATGGTTGGAGG
	(SEQ ID NO: 22)	
199	CDKN2B	(SEQ ID NO: 722) ATATTTAGCGAGTAGTGT
	i	
200	(SEQ ID NO: 23) CDKN2B	(SEQ ID NO: 723) ATATTTAGTGAGTAGTGT
	1	
201	(SEQ ID NO: 23) CDKN2B	(SEQ ID NO: 724) ATAGGGGGGGGGGTTTAA
	[
202	(SEQ ID NO: 23) CDKN2B	(SEQ ID NO: 725) ATAGGGGGTGGAGTTTAA
202		
203	(SEQ ID NO: 23) CDKN2B	(SEQ ID NO: 726) TTATTGTACGGGGTTTTA
200		
204	(SEQ ID NO: 23) CDKN2B	(SEQ ID NO: 727) TTATTGTATGGGGTTTTA
·		
205	(SEQ ID NO: 23) CDKN2B	(SEQ ID NO: 728) TTTTAAGTCGTAGAAGGA
206	(SEQ ID NO: 23) CDKN2B	(SEQ ID NO: 729) TTTTAAGTTGTAGAAGGA
230		
	(SEQ ID NO: 23)	(SEQ ID NO: 730)

No:	Gene	Oligo:
207	DAPK1	GTTGGAGTCGAGGTTTGA
	(SEQ ID NO: 25)	(SEQ ID NO: 731)
208	DAPK1	GTTGGAGTTTGA
j	(SEQ ID NO: 25)	(SEQ ID NO: 732)
209	DAPK1	GAAGGAGCGTATTTAT
	(SEQ ID NO: 25)	(SEQ ID NO: 733)
210	DAPK1	GAAGGAGTGTATTTAT
	(SEQ ID NO: 25)	(SEQ ID NO: 734)
211	DAPK1	TTGTTTTCGGAAATTTG
	(SEQ ID NO: 25)	(SEQ ID NO: 735)
.212	DAPK1	TTGTTTTTGGAAATTTG
	(SEQ ID NO: 25)	(SEQ ID NO: 736)
213	FOS	AATGTTTTCGTACGTAGG
	(SEQ ID NO: 29)	(SEQ ID NO: 737)
214	FOS	AATGTTTTGTATGTAGG
	(SEQ ID NO: 29)	(SEQ ID NO: 738)
215	FOS	TATATGGTCGAGAAAAT
	(SEQ ID NO: 29)	(SEQ ID NO: 739)
216	FOS	TATATGGTTGAGAAAAT
Ì	(SEQ ID NO: 29)	(SEQ ID NO: 740)
217	FOS	TTTAGTATCGTAAAGTAG
	(SEQ ID NO: 29)	(SEQ ID NO: 741)
218	FOS	TTTAGTATTGTAAAGTAG
	(SEQ ID NO: 29)	(SEQ ID NO: 742)
219	FOS	GTATTGTTCGAGTTCGAG
	(SEQ ID NO: 29)	(SEQ ID NO: 743)
220	FOS	GTATTGTTTGAG
	(SEQ ID NO: 29)	(SEQ ID NO: 744)
221	GSTP1	GGTTTTTCGGTTAGTTG
	(SEQ ID NO: 33)	(SEQ ID NO: 745)
222	GSTP1	GGTTTTTTGGTTAGTTG
	(SEQ ID NO: 33)	(SEQ ID NO: 746)
223	GSTP1	TTTTAGGGCGTTTTTTTG
	(SEQ ID NO: 33)	(SEQ ID NO: 747)
224	GSTP1	TTTTAGGGTGTTTTTTG
	(SEQ ID NO: 33)	(SEQ ID NO: 748)
225	GSTP1	GTAGTTTTCGTTATTAGT
	(SEQ ID NO: 33)	(SEQ ID NO: 749)

No:	Gene	Oligo:
226	GSTP1	GTAGTTTTTGTTATTAGT
	(SEQ ID NO: 33)	(SEQ ID NO: 750)
227	HIC-1	TTAAAACGGCGTATAGGG
	(SEQ ID NO: 34)	(SEQ ID NO: 751)
228	HIC-1	TTAAAATGGTGTATAGGG
	(SEQ ID NO: 34)	(SEQ ID NO: 752)
229	HIC-1	AGGAGATTCGAAAGTTTA
	(SEQ ID NO: 34)	(SEQ ID NO: 753)
230	HIC-1	AGGAGATTTGAAAGTTTA
	(SEQ ID NO: 34)	(SEQ ID NO: 754)
231	HIC-1	GGGTTTTACGTGGTTGTT
	(SEQ ID NO: 34)	(SEQ ID NO: 755)
232	HIC-1	GGGTTTTATGTGGTTGTT
	(SEQ ID NO: 34)	(SEQ ID NO: 756)
233	HIC-1	TTTTAGAGCGTTAGGGTT
	(SEQ ID NO: 34)	(SEQ ID NO: 757)
234	HIC-1	TTTTAGAGTGTTAGGGTT
	(SEQ ID NO: 34)	(SEQ ID NO: 758)
235	IGF2	AGTTTGAACGATGTAAGA
1	(SEQ ID NO: 36)	(SEQ ID NO: 759)
236	IGF2	AGTTTGAATGATGTAAGA
	(SEQ ID NO: 36)	(SEQ ID NO: 760)
237	IGF2	GGTTATTACGATAATTTG
	(SEQ ID NO: 36)	(SEQ ID NO: 761)
238	IGF2	GGTTATTATGATAATTTG
	(SEQ ID NO: 36)	(SEQ ID NO: 762)
239	IGF2	TTGTATGGTCGAGTTTAT
	(SEQ ID NO: 36)	(SEQ ID NO: 763)
240	IGF2	TTGTATGGTTGAGTTTAT
	(SEQ ID NO: 36)	(SEQ ID NO: 764)
241	IGF2	GATTAGGGCGGAAATAT
	(SEQ ID NO: 36)	(SEQ ID NO: 765)
242	IGF2	GATTAGGGTGGGAAATAT
	(SEQ ID NO: 36)	(SEQ ID NO: 766)
243	IGF2	TGGAGTTTACGGAGGTTT
	(SEQ ID NO: 36)	(SEQ ID NO: 767)
244	IGF2	TGGAGTTTATGGAGGTTT
	(SEQ ID NO: 36)	(SEQ ID NO: 768)

No:	Gene	Oligo:
245	MGMT	TAAGGATACGAGTTATAT
	(SEQ ID NO: 38)	(SEQ ID NO: 769)
246	MGMT	TAAGGATATGAGTTATAT
	(SEQ ID NO: 38)	(SEQ ID NO: 770)
247	MGMT	TTGGAGAGCGGTTGAGTT
	(SEQ ID NO: 38)	(SEQ ID NO: 771)
248	MGMT	TTGGAGAGTGGTTGAGTT
	(SEQ ID NO: 38)	(SEQ ID NO: 772)
249	MGMT	TAGGTTATCGGTGATTGT
	(SEQ ID NO: 38)	(SEQ ID NO: 773)
250	MGMT	TAGGTTATTGGTGATTGT
	(SEQ ID NO: 38)	(SEQ ID NO: 774)
251	MGMT	TAGGGAGCGGTTTTAGG
	(SEQ ID NO: 38)	(SEQ ID NO: 775)
252	MGMT	TAGGGAGTGGTTTTAGG
	(SEQ ID NO: 38)	(SEQ ID NO: 776)
253	MGMT	AGTAGGATCGGGATTTT
	(SEQ ID NO: 38)	(SEQ ID NO: 777)
254	MGMT	AGTAGGATTGTTT
	(SEQ ID NO: 38)	(SEQ ID NO: 778)
255	MLH1	TTGAGAAGCGTTAAGTAT
	(SEQ ID NO: 39)	(SEQ ID NO: 779)
256	MLH1	TTGAGAAGTGTTAAGTAT
	(SEQ ID NO: 39)	(SEQ ID NO: 780)
257	MLH1	TTAGGTAGCGGGTAGTAG
	(SEQ ID NO: 39)	(SEQ ID NO: 781)
258	MLH1	TTAGGTAGTAGTAG
	(SEQ ID NO: 39)	(SEQ ID NO: 782)
259	MLH1	GTAGTAGTCGTTTTAGGG
	(SEQ ID NO: 39)	(SEQ ID NO: 783)
260	MLH1	GTAGTAGTTGTTTTAGGG
	(SEQ ID NO: 39)	(SEQ ID NO: 784)
261	MLH1	ATAGTTGTCGTTGAAGGG
	(SEQ ID NO: 39)	(SEQ ID NO: 785)
262	MLH1	ATAGTTGTTGAAGGG
	(SEQ ID NO: 39)	(SEQ ID NO: 786)
263	MLH1	TTGGATGGCGTAAGTTAT
	(SEQ ID NO: 39)	(SEQ ID NO: 787)

No:	Gene	Oligo:
264	MLH1	TTGGATGGTGTAAGTTAT
ļ	(SEQ ID NO: 39)	(SEQ ID NO: 788)
265	MOS	AGTAGTTTCGTAGGTAGT
	(SEQ ID NO: 40)	(SEQ ID NO: 789)
266	MOS	AGTÄGTTTTGTAGGTAGT
	(SEQ ID NO: 40)	(SEQ ID NO: 790)
267	MOS	GTAAGTCGTTTTGTATAT
	(SEQ ID NO: 40)	(SEQ ID NO: 791)
268	MOS	GTAAGTTGTTTTGTATAT
	(SEQ ID NO: 40)	(SEQ ID NO: 792)
269	MOS	ATGTTAGTCGGTTTTTGG
	(SEQ ID NO: 40)	(SEQ ID NO: 793)
270	MOS	ATGTTAGTTGGTTTTTGG
	(SEQ ID NO: 40)	(SEQ ID NO: 794)
271	MPL	GTGGTGGCGTTTGTAAT
	(SEQ ID NO: 42)	(SEQ ID NO: 795)
272	MPL	GTGGTGGTGTTTGTAAT
	(SEQ ID NO: 42)	(SEQ ID NO: 796)
273	MPL	TTTGAGGTCGGGAGTTTA
1	(SEQ ID NO: 42)	(SEQ ID NO: 797)
274	MPL	TTTGAGGTTGGGAGTTTA
	(SEQ ID NO: 42)	(SEQ ID NO: 798)
275	MPL	TGTAGTGAGTCGAGATTA
	(SEQ ID NO: 42)	(SEQ ID NO: 1227)
276	MPL	TGTAGTGAGTTA
	(SEQ ID NO: 42)	(SEQ ID NO: 1228)
277	MPL	GGTAATAGAGCGAGATTT
	(SEQ ID NO: 42)	(SEQ ID NO: 799)
278	MPL	GGTAATAGAGTGAGATTT
	(SEQ ID NO: 42)	(SEQ ID NO: 800)
279	MPL	TTAAGTAGCGGGTAAGGT
	(SEQ ID NO: 42)	(SEQ ID NO: 801)
280	MPL	TTAAGTAGTGGGTAAGGT
	(SEQ ID NO: 42)	(SEQ ID NO: 802)
281	MYC	TTAGAGTGTTCGGTTGTT
	(SEQ ID NO: 43)	(SEQ ID NO: 803)
282	MYC	TTAGAGTGTTTGGTT
	(SEQ ID NO: 43)	(SEQ ID NO: 804)

No:	Gene	Oligo:
283	MYC	TTATAATGCGAGGGTTTG
	(SEQ ID NO: 43)	(SEQ ID NO: 805)
284	MYC	TTATAATGTGAGGGTTTG
	(SEQ ID NO: 43)	(SEQ ID NO: 806)
285	MYC	AGGATTTTCGAGTTGTGT
	(SEQ ID NO: 43)	(SEQ ID NO: 807)
286	MYC	AGGATTTTTGAGTTGTGT
	(SEQ ID NO: 43)	(SEQ ID NO: 808)
287	MYC	AATTTTAGCGAGAGGTAG
	(SEQ ID NO: 43)	(SEQ ID NO: 809)
288	MYC	AATTTTAGTGAGAGGTAG
	(SEQ ID NO: 43)	(SEQ ID NO: 810)
289	MYC	TTGTGGGCGTTTTGGGAA
	(SEQ ID NO: 43)	(SEQ ID NO: 811)
290	MYC	TTGTGGGTGTTTTGGGAA
	(SEQ ID NO: 43)	(SEQ ID NO: 812)
291	N33	TTGGTTCGGGAAAGGTAA
	(SEQ ID NO: 46)	(SEQ ID NO: 1211)
292	N33	TTGGTTTGGGAAAGGTAA
···	(SEQ ID NO: 46)	(SEQ ID NO: 1212)
293	N33	TGTTATTTCGGAGGGTTT
	(SEQ ID NO: 46)	(SEQ ID NO: 1217)
294	N33	TGTTATTTTGGAGGGTTT
	(SEQ ID NO: 46)	(SEQ ID NO: 1218)
295	N33	GTTTAGTTAGCGGGTTTT
	(SEQ ID NO: 46)	(SEQ ID NO: 1219)
296	N33	GTTTAGTTAGTGGGTTTT
	(SEQ ID NO: 46)	(SEQ ID NO: 1220)
297	N33	ATTTAGTTCGGGGGAGGA
	(SEQ ID NO: 46)	(SEQ ID NO: 1215)
298	N33	ATTTAGTTTGGGGGAGGA
	(SEQ ID NO: 46)	(SEQ ID NO: 1216)
299	PMS2	AGATTATACGTTGTATGT
200	(SEQ ID NO: 49)	(SEQ ID NO: 813)
300	PMS2	AGATTATATGTTGTATGT
201	(SEQ ID NO: 49)	(SEQ ID NO: 814)
301	PMS2	TAGGAGTTCGAGATTAGT
	(SEQ ID NO: 49)	(SEQ ID NO: 815)

No:	Gene	Oligo:
302	PMS2	TAGGAGTTTGAGATTAGT
	(SEQ ID NO: 49)	(SEQ ID NO: 816)
303	PMS2	AAATTAGTCGGGTTTAGT
}	(SEQ ID NO: 49)	(SEQ ID NO: 817)
304	PMS2	AAATTAGTTGGGTTTAGT
	(SEQ ID NO: 49)	(SEQ ID NO: 818)
305	PMS2	GTGGATAGCGTTTGTAAT
	(SEQ ID NO: 49)	(SEQ ID NO: 819)
306	PMS2	GTGGATAGTGTTTGTAAT
	(SEQ ID NO: 49)	(SEQ ID NO: 820)
307	PTEN	GGATTTTGCGTTCGTATT
	(SEQ ID NO: 51)	(SEQ ID NO: 821)
308	PTEN	GGATTTTGTGTTTTT
	(SEQ ID NO: 51)	(SEQ ID NO: 822)
309	PTEN	AGAGTTATCGTTTTGTTT
	(SEQ ID NO: 51)	(SEQ ID NO: 823)
310	PTEN	AGAGTTATTGTTTTGTTT
	(SEQ ID NO: 51)	(SEQ ID NO: 824)
311	PTEN	TGATGTGGCGGGATTTTT
	(SEQ ID NO: 51)	(SEQ ID NO: 825)
312	PTEN	TGATGTGGTGGGATTTTT
÷	(SEQ ID NO: 51)	(SEQ ID NO: 826)
313	RBL2	ATTAGTGTCGTTGTTAAG
	(SEQ ID NO: 53)	(SEQ ID NO: 827)
314	RBL2	ATTAGTGTTGTTAAG
	(SEQ ID NO: 53)	(SEQ ID NO: 828)
315	RBL2	AGATTATACGGATAAGGG
	(SEQ ID NO: 53)	(SEQ ID NO: 829)
316	RBL2	AGATTATATGGATAAGGG
	(SEQ ID NO: 53)	(SEQ ID NO: 830)
317	SFN	ATAGAGTTCGGTATTGGT
	(SEQ ID NO: 55)	(SEQ ID NO: 831)
318	SFN	ATAGAGTTTGGTATTGGT
	(SEQ ID NO: 55)	(SEQ ID NO: 832)
319	SFN	GAGTAGGTCGAACGTTAT
	(SEQ ID NO: 55)	(SEQ ID NO: 833)
320	SFN	GAGTAGGTTGAATGTTAT
	(SEQ ID NO: 55)	(SEQ ID NO: 834)

No:	Gene	Oligo:
321	SFN	AAAAGTAACGAGGAGGGT
	(SEQ ID NO: 55)	(SEQ ID NO: 835)
322	SFN	AAAAGTAATGAGGAGGGT
	(SEQ ID NO: 55)	(SEQ ID NO: 836)
323	TGFBR2	ATTTGGAGCGAGGAATTT
	(SEQ ID NO: 57)	(SEQ ID NO: 837)
324	TGFBR2	ATTTGGAGTGAGGAATTT
	(SEQ ID NO: 57)	(SEQ ID NO: 838)
325	TGFBR2	TTGAAAGTCGGTTAAAGT
	(SEQ ID NO: 57)	(SEQ ID NO: 839)
326	TGFBR2	TTGAAAGTTGGTTAAAGT
	(SEQ ID NO: 57)	(SEQ ID NO: 840)
327	TGFBR2	AAAGTTTTCGGAGGGTT
	(SEQ ID NO: 57)	(SEQ ID NO: 841)
328	TGFBR2	AAAGTTTTTGGAGGGTT
'	(SEQ ID NO: 57)	(SEQ ID NO: 842)
329	TGFBR2	GGTAGTTACGAGAGTT
	(SEQ ID NO: 57)	(SEQ ID NO: 843)
330	TGFBR2	GGTAGTTATGAGAGAGTT
	(SEQ ID NO: 57)	(SEQ ID NO: 844)
331	TP73	TGATTTAGCGTAGGTTTG
	(SEQ ID NO: 58)	(SEQ ID NO: 845)
332	TP73	TGATTTAGTGTAGGTTTG
	(SEQ ID NO: 58)	(SEQ ID NO: 846)
333	TP73	TTTGGTGCGCGTAGAGAA
	(SEQ ID NO: 58)	(SEQ ID NO: 1229)
334	TP73	TTTGGTGTGTGAGAA
	(SEQ ID NO: 58)	(SEQ ID NO: 1230)
335	TP73	TTAGAGTTCGAGTTTATA
	(SEQ ID NO: 58)	(SEQ ID NO: 847)
336	TP73	TTAGAGTTTGAGTTTATA
	(SEQ ID NO: 58)	(SEQ ID NO: 848)
337	TP73	GAGTGTTCGCGTTTTGGG
	(SEQ ID NO: 58)	(SEQ ID NO: 849)
338	TP73	GAGTGTTTGTGTTTTGGG
	(SEQ ID NO: 58)	(SEQ ID NO: 850)
339	TP73	AAGATGTGCGAGTTAGTC
	(SEQ ID NO: 58)	(SEQ ID NO: 851)

No:	Gene	Oligo:
340	TP73	AAGATGTGTGAGTTAGTC
	(SEQ ID NO: 58)	(SEQ ID NO: 852)
341	TP73	AAGTTACGGGTTTTATTG
	(SEQ ID NO: 58)	(SEQ ID NO: 853)
342	TP73	AAGTTATGGGTTTTATTG
	(SEQ ID NO: 58)	(SEQ ID NO: 854)
343	TP73	TTAGATTACGGGTTTTAT
	(SEQ ID NO: 58)	(SEQ ID NO: 855)
344	TP73	TTAGATTATGGGTTTTAT
	(SEQ ID NO: 58)	(SEQ ID NO: 856)
345	TP73	TAAGTAGCGTCGTTATTG
	(SEQ ID NO: 58)	
346	TP73	(SEQ ID NO: 857) TAAGTAGTGTTGTTATTG
	(SEQ ID NO: 58)	
347	TP73	(SEQ ID NO: 858) GGAAGTTTCGATGGTTTA
J .,		
348	(SEQ ID NO: 58) TP73	(SEQ ID NO: 859) GGAAGTTTTGATGGTTTA
340		
349	(SEQ ID NO: 58) CDKN1C	(SEQ ID NO: 860)
349		ATGAAGAACGGTTAAGGG
350	(SEQ ID NO: 21) CDKN1C	(SEQ ID NO: 861)
330		ATGAAGAATGGTTAAGGG
251	(SEQ ID NO: 21)	(SEQ ID NO: 862)
351	CDKN1C	TTAAGTTACGGTTATTAG
0.70	(SEQ ID NO: 21)	(SEQ ID NO: 863)
352	CDKN1C	TTAAGTTATGGTTATTAG
2.72	(SEQ ID NO: 21)	(SEQ ID NO: 864)
353	CDKN1C	TTAGTGTTCGTTTGGAAT
	(SEQ ID NO: 21)	(SEQ ID NO: 865)
354	CDKN1C	TTAGTGTTTGGAAT
	(SEQ ID NO: 21)	(SEQ ID NO: 866)
355	CDKN1C	AGGAGTTGCGGTGGGAAT
	(SEQ ID NO: 21)	(SEQ ID NO: 867)
356	CDKN1C	AGGAGTTGTGGTGGAAT
	(SEQ ID NO: 21)	(SEQ ID NO: 868)
357	GSK3ß	GGGTAAAGCGCGGÁTATT
	(SEQ ID NO: 32)	(SEQ ID NO: 869)
358	GSK3ß	GGGTAAAGTGTGGATATT
	(SEQ ID NO: 32)	(SEQ ID NO: 870)

No:	Gene	Oligo:
359	GSK3ß	TATGTTTTCGGCGAATGG
	(SEQ ID NO: 32)	(SEQ ID NO: 871)
360	GSK3ß	TATGTTTTTGGTGAATGG
	(SEQ ID NO: 32)	(SEQ ID NO: 872)
361	GSK3ß	GGGGAATAGTCGAGGAGT
	(SEQ ID NO: 32)	(SEQ ID NO: 873)
362	GSK3ß	GGGGAATAGTTGAGGAGT
	(SEQ ID NO: 32)	(SEQ ID NO: 874)
363	GSK3ß	AGGAGTCGTTGTTTGGGG
	(SEQ ID NO: 32)	(SEQ ID NO: 875)
364	GSK3ß	AGGAGTTGTTTGGGG
	(SEQ ID NO: 32)	(SEQ ID NO: 876)
365	CDC2	TAGTTATTCGGGAAGGTA
	(SEQ ID NO: 14)	
366	CDC2	(SEQ ID NO: 877) TAGTTATTTGGGAAGGTA
367	(SEQ ID NO: 14) CDC2	(SEQ ID NO: 878) AGGTATTGCGGTAGTTGG
368	(SEQ ID NO: 14) CDC2	(SEQ ID NO: 879) AGGTATTGTGGTAGTTGG
		·
369	(SEQ ID NO: 14) CDC2	(SEQ ID NO: 880) GGGTTATTCGATTGGTGA
307		
370	(SEQ ID NO: 14) CDC2	(SEQ ID NO: 881) GGGTTATTTGATTGGTGA
370		
371	(SEQ ID NO: 14) CDC2	(SEQ ID NO: 882)
3/1]	AAATTGTTCGTATTTGGT
372	(SEQ ID NO: 14) CDC2	(SEQ ID NO: 883)
		AAATTGTTTGTATTTGGT
373	(SEQ ID NO: 14) ESR1	(SEQ ID NO: 884)
313	ļ	AGATATATCGGAGTTTGG
374	(SEQ ID NO: 28)	(SEQ ID NO: 885)
3/4	ESR1	AGATATATTGGAGTTTGG
275	(SEQ ID NO: 28)	(SEQ ID NO: 886)
375	ESR1	GTTTGGTACGGGGTATAT
276	(SEQ ID NO: 28)	(SEQ ID NO: 887)
376	ESR1	GTTTGGTATGGGGTATAT
077	(SEQ ID NO: 28)	(SEQ ID NO: 888)
377	ESR1	TTTTAAATCGAGTTGTGT
	(SEQ ID NO: 28)	(SEQ ID NO: 889)

No:	Gene	Oligo:
378	ESR1	TTTTAAATTGAGTTGTGT
	(SEQ ID NO: 28)	(SEQ ID NO: 890)
379	ESR1	TATGAGTTCGGGAGATTA
	(SEQ ID NO: 28)	(SEQ ID NO: 891)
. 380	ESR1	TATGAGTTTGGGAGATTA
-	(SEQ ID NO: 28)	(SEQ ID NO: 892)
381	ESR1	TGGAGGTTCGGGAGTTTA
	(SEQ ID NO: 28)	(SEQ ID NO: 893)
382	ESR1	TGGAGGTTTGGGAGTTTA
	(SEQ ID NO: 28)	(SEQ ID NO: 894)
383	APAF1	TTTGGTATCGTTTAGAGT
	(SEQ ID NO: 3)	(SEQ ID NO: 895)
384	APAF1	TTTGGTATTGTTTAGAGT
	(SEQ ID NO: 3)	(SEQ ID NO: 896)
385	APAF1	GTATGAGTCGTGGTAGGA
	(SEQ ID NO: 3)	(SEQ ID NO: 897)
386	APAF1	GTATGAGTTGTGGTAGGA
	(SEQ ID NO: 3)	(SEQ ID NO: 898)
387	APAF1	GTGGATTCGGCGGATTT
	(SEQ ID NO: 3)	(SEQ ID NO: 899)
388	APAF1	GTGGATTTGGTGGGATTT
	(SEQ ID NO: 3)	(SEQ ID NO: 900)
389	APAF1	TTTAGAGGCGGAGAAGAA
	(SEQ ID NO: 3)	(SEQ ID NO: 901)
390	APAF1	TTTAGAGGTGGAGAAGAA
	(SEQ ID NO: 3)	(SEQ ID NO: 902)
391	APAF1	GAAGAGGTAGCGAGTGGA
	(SEQ ID NO: 3)	(SEQ ID NO: 903)
392	APAF1	GAAGAGGTAGTGGA
	(SEQ ID NO: 3)	(SEQ ID NO: 904)
393	BAK1	TAGGTTGTCGGTTTGTGC
[(SEQ ID NO: 7)	(SEQ ID NO: 905)
394	BAK1	TAGGTTGTTGTGC
	(SEQ ID NO: 7)	(SEQ ID NO: 906)
395	BAK1	TTTGTATTCGGTGGTTAT
	(SEQ ID NO: 7)	(SEQ ID NO: 907)
396	BAK1	TTTGTATTTGGTGGTTAT
	(SEQ ID NO: 7)	(SEQ ID NO: 908)
	(3-2-1,0,7)	(BLQ ID 140. 300)

No:	Gene	Oligo:
397	BAK1	GGAGTTTCGCGGGTTTTT
	(SEQ ID NO: 7)	(SEQ ID NO: 909)
398	BAK1	GGAGTTTTGTGGGTTTTT
	(SEQ ID NO: 7)	(SEQ ID NO: 910)
399	BAK1	TAGGATTTCGGTAGGTAA
	(SEQ ID NO: 7)	(SEQ ID NO: 911)
400	BAK1	TAGGATTTTGGTAGGTAA
	(SEQ ID NO: 7)	(SEQ ID NO: 912)
401	BAX	AGTTTGGGCGTGGGTTAT
	(SEQ ID NO: 8)	(SEQ ID NO: 913)
402	BAX	AGTTTGGGTGTGGGTTAT
}	(SEQ ID NO: 8)	(SEQ ID NO: 914)
403	BAX	ATTAGAGTTGCGATTGGA
	(SEQ ID NO: 8)	(SEQ ID NO: 915)
404	BAX	ATTAGAGTTGTGATTGGA
	(SEQ ID NO: 8)	(SEQ ID NO: 916)
405	BAX	GTATTTATCGGGAGATGT
	(SEQ ID NO: 8)	(SEQ ID NO: 917)
406	BAX	GTATTTATTGGGAGATGT
	(SEQ ID NO: 8)	(SEQ ID NO: 918)
407	BAX	TTTAGAGGCGGGGTGAG
	(SEQ ID NO: 8)	(SEQ ID NO: 919)
408	BAX	TTTAGAGGTGGGGTGAG
	(SEQ ID NO: 8)	(SEQ ID NO: 920)
409	CASP10	GAAGTAAACGGTTTTAG
	(SEQ ID NO: 10)	(SEQ ID NO: 921)
410	CASP10 .	GAAGTAAATGGTTTTAG
	(SEQ ID NO: 10)	(SEQ ID NO: 922)
411	CASP10	TTATAATACGAAGTTATG
	(SEQ ID NO: 10)	(SEQ ID NO: 923)
412	CASP10	TTATAATATGAAGTTATG
	(SEQ ID NO: 10)	(SEQ ID NO: 924)
413	CASP10	GAGGTTTTCGGATTTTT
	(SEQ ID NO: 10)	(SEQ ID NO: 925)
414	CASP10	GAGGTTTTTGGATTTTT
	(SEQ ID NO: 10)	(SEQ ID NO: 926)
415	CASP10	TTTTGTTTCGGTAAAAGG
	(SEQ ID NO: 10)	(SEQ ID NO: 927)
	(22 (110. 10)	(SLQ ID 1(0. 321)

No:	Gene	Oligo:
416	CASP10	TTTTGTTTTGGTAAAAGG
	(SEQ ID NO: 10)	(SEQ ID NO: 928)
417	CASP8	GAATGAGTCGAGGAAGGT
1	(SEQ ID NO: 11)	(SEQ ID NO: 929)
418	CASP8	GAATGAGTTGAGGAAGGT
	(SEQ ID NO: 11)	(SEQ ID NO: 930)
419	CASP8	TATTGAGACGTTAAGTAA
	(SEQ ID NO: 11)	(SEQ ID NO: 931)
420	CASP8	TATTGAGATGTTAAGTAA
	(SEQ ID NO: 11)	(SEQ ID NO: 932)
421	CASP8	TAAGGTTACGTAGTTAGT
1	(SEQ ID NO: 11)	(SEQ ID NO: 933)
422	CASP8	TAAGGTTATGTAGT
[(SEQ ID NO: 11)	(SEQ ID NO: 934)
423	CASP8	GTTAATAGCGGGGATTTT
	(SEQ ID NO: 11)	(SEQ ID NO: 935)
424	CASP8	GTTAATAGTGGGGATTTT
	(SEQ ID NO: 11)	(SEQ ID NO: 936)
425	CASP9	TATATGATCGAGGATATT
	(SEQ ID NO: 12)	(SEQ ID NO: 937)
426	CASP9	TATATGATTGAGGATATT
ł	(SEQ ID NO: 12)	(SEQ ID NO: 938)
427	CASP9	ATTTTTAGTCGGGTGTT
j	(SEQ ID NO: 12)	(SEQ ID NO: 939)
428	CASP9	ATTTTTAGTTGGGTGTT
	(SEQ ID NO: 12)	(SEQ ID NO: 940)
429	CASP9	TTAGTTAGCGTATATTGT
	(SEQ ID NO: 12)	(SEQ ID NO: 941)
430	CASP9	TTAGTTAGTGTATATTGT
	(SEQ ID NO: 12)	(SEQ ID NO: 942)
431	CASP9	TGATAAAACGTTAGAGGT
	(SEQ ID NO: 12)	(SEQ ID NO: 943)
432	CASP9	TGATAAAATGTTAGAGGT
	(SEQ ID NO: 12)	(SEQ ID NO: 944)
433	TCL1A	TTAGAGTACGTTTTTTGG
	(SEQ ID NO: 56)	(SEQ ID NO: 945)
434	TCL1A	TTAGAGTATGTTTTTGG
	(SEQ ID NO: 56)	(SEQ ID NO: 946)

No:	Gene	Oligo:
435	TCL1A	TTTATAGTCGATATTAGG
	(SEQ ID NO: 56)	(SEQ ID NO: 947)
436	TCL1A	TTTATAGTTGATATTAGG
	(SEQ ID NO: 56)	(SEQ ID NO: 948)
437	TCL1A	TAGATTCGGTTTTAGTTG
	(SEQ ID NO: 56)	(SEQ ID NO: 949)
438	TCL1A	TAGATTTGGTTTTAGTTG
	(SEQ ID NO: 56)	(SEQ ID NO: 950)
439	TCL1A	TTTTGGAACGGTGGTAGT
	(SEQ ID NO: 56)	(SEQ ID NO: 951)
440	TCL1A	TTTTGGAATGGTGGTAGT
	(SEQ ID NO: 56)	(SEQ ID NO: 952)
441	PML	TTAAGAGTCGTAGATGTA
	(SEQ ID NO: 48)	(SEQ ID NO: 953)
442	PML	TTAAGAGTTGTAGATGTA
	(SEQ ID NO: 48)	(SEQ ID NO: 954)
443	PML	GAGTTGGGCGTATAAAGA
	(SEQ ID NO: 48)	(SEQ ID NO: 955)
444	PML	GAGTTGGGTGTATAAAGA
	(SEQ ID NO: 48)	(SEQ ID NO: 956)
445	PML	GTAAAGCGGGAGAGGTAG
	(SEQ ID NO: 48)	(SEQ ID NO: 957)
446	PML	GTAAAGTGGGAGAGGTAG
	(SEQ ID NO: 48)	(SEQ ID NO: 958)
447	PML	TTAAAGATACGTTTAGAG
	(SEQ ID NO: 48)	(SEQ ID NO: 959)
448	PML	TTAAAGATATGTTTAGAG
	(SEQ ID NO: 48)	(SEQ ID NO: 960)
449	HOXA5	AGTTAGTCGGGTTTTAAG
	(SEQ ID NO: 35)	(SEQ ID NO: 961)
450	HOXA5	AGTTAGTTGGGTTTTAAG "
	(SEQ ID NO: 35)	(SEQ ID NO: 962)
451	HOXA5	TTATAGGGTTCGGTTTTT
	(SEQ ID NO: 35)	(SEQ ID NO: 963)
452	HOXA5	TTATAGGGTTTGGTTTTT
	(SEQ ID NO: 35)	(SEQ ID NO: 964)
453	HOXA5	TTTTAAGGCGAGGTTAAA
	(SEQ ID NO: 35)	(SEQ ID NO: 965)

No:	Gene	Oligo:
454 -	HOXA5	TTTTAAGGTGAGGTTAAA
	(SEQ ID NO: 35)	(SEQ ID NO: 966)
455	HOXA5	ATGATAGGCGTTTATTAA
	(SEQ ID NO: 35)	(SEQ ID NO: 967)
456	HOXA5	ATGATAGGTGTTTATTAA
	(SEQ ID NO: 35)	(SEQ ID NO: 968)
457	SDC4	GGTTATACGATTTTTGTG
	(SEQ ID NO: 54)	(SEQ ID NO: 969)
458	SDC4	GGTTATATGATTTTTGTG
	(SEQ ID NO: 54)	(SEQ ID NO: 970)
459	SDC4	TAGGTTTTGCGGGTTATA
	(SEQ ID NO: 54)	(SEQ ID NO: 971)
460	SDC4	TAGGTTTTGTGGGTTATA
	(SEQ ID NO: 54)	(SEQ ID NO: 972)
461	SDC4	ATGTGAATCGTGTTAAGA
	(SEQ ID NO: 54)	(SEQ ID NO: 973)
462	SDC4	ATGTGAATTGTGTTAAGA
	(SEQ ID NO: 54)	(SEQ ID NO: 974)
463	SDC4	TTTTATATCGGGTGTGTT
	(SEQ ID NO: 54)	
464	SDC4	(SEQ ID NO: 1223) TTTTATATTGGGTGTGTT
	(SEQ ID NO: 54)	
465	PITX2	(SEQ ID NO: 1224) TTAAGAGACGGAATAAAG
	(SEQ ID NO: 47)	,
466	PITX2	(SEQ ID NO: 975) TTAAGAGATGGAATAAAG
467	(SEQ ID NO: 47) PITX2	(SEQ ID NO: 976) AATTAGGACGGGTTGGGA
468	(SEQ ID NO: 47) PITX2	(SEQ ID NO: 977) AATTAGGATGGGTTGGGA
100		
469	(SEQ ID NO: 47) PITX2	(SEQ ID NO: 978) TAGGAGTACGGATTTTA
105		·
470	(SEQ ID NO: 47) PITX2	(SEQ ID NO: 979) TAGGAGTATGGATTTTA
170		
471	(SEQ ID NO: 47) PITX2	(SEQ ID NO: 980) GAGTTAATCGAGGAGTAG
.,,		
472	(SEQ ID NO: 47) PITX2	(SEQ ID NO: 981)
712		GAGTTAATTGAGGAGTAG
	(SEQ ID NO: 47)	(SEQ ID NO: 982)

No:	Gene	Oligo:
473	PITX2	GGATTTGGCGGGGTATTT
	(SEQ ID NO: 47)	(SEQ ID NO: 983)
474	PITX2	GGATTTGGTGGGGTATTT
	(SEQ ID NO: 47)	(SEQ ID NO: 984)
475	GPR37	AGTTTAGACGAAGTTTTA
	(SEQ ID NO: 31)	(SEQ ID NO: 985)
476	GPR37	AGTTTAGATGAAGTTTTA
	(SEQ ID NO: 31)	(SEQ ID NO: 986)
477	GPR37	GTTATTTCGGTGGTTGA
	(SEQ ID NO: 31)	(SEQ ID NO: 987)
478	GPR37	GTTÀTTTTGGTGGTTGA
	(SEQ ID NO: 31)	(SEQ ID NO: 988)
479	GPR37	AGGGTGATAGCGGAAGAT
	(SEQ ID NO: 31)	(SEQ ID NO: 989)
480	GPR37	AGGTGATAGTGGAAGAT
	(SEQ ID NO: 31)	(SEQ ID NO: 990)
481	GPR37	GTTAGGATCGATATTTGT
	(SEQ ID NO: 31)	(SEQ ID NO: 991)
482	GPR37	GTTAGGATTGATATTTGT
	(SEQ ID NO: 31)	(SEQ ID NO: 992)
483	PRAME	GGGATTGTCGTATAAAAT
	(SEQ ID NO: 50)	(SEQ ID NO: 993)
484	PRAME	GGGATTGTTGTATAAAAT
	(SEQ ID NO: 50)	(SEQ ID NO: 994)
485	PRAME	TTTAGATTCGATTTGGGA
	(SEQ ID NO: 50)	(SEQ ID NO: 995)
486	PRAME	TTTAGATTTGGGA
	(SEQ ID NO: 50)	(SEQ ID NO: 996)
487	PRAME	TAGÀGTTGCGGGTTTTTA
	(SEQ ID NO: 50)	(SEQ ID NO: 997)
488	PRAME	TAGAGTTGTGGGTTTTTA
ļ	(SEQ ID NO: 50)	(SEQ ID NO: 998)
489	PRAME	AATAGGTCGCGTGGTGAG
	(SEQ ID NO: 50)	(SEQ ID NO: 999)
490	PRAME	AATAGGTTGTGTGAG
	(SEQ ID NO: 50)	(SEQ ID NO: 1000)
491	N-MYC	AGTTTTGGCGTTTTTTT
	(SEQ ID NO: 60)	(SEQ ID NO: 1001)

No:	Gene	Oligo:
492	N-MYC	AGTTTTGGTGTTTTTTT
	(SEQ ID NO: 60)	(SEQ ID NO: 1002)
493	N-MYC	GGTAAAGTCGTTTTTTT
	(SEQ ID NO: 60)	(SEQ ID NO: 1003)
494	N-MYC	GGTAAAGTTGTTTTTTT
	(SEQ ID NO: 60)	(SEQ ID NO: 1004)
495	N-MYC	TTTGTTTGCGTTATAGTT
	(SEQ ID NO: 60)	(SEQ ID NO: 1005)
496	N-MYC	TTTGTTTGTGTTATAGTT
	(SEQ ID NO: 60)	(SEQ ID NO: 1006)
497	N-MYC	ATATTTTCGAGTTTTAA
•	(SEQ ID NO: 60)	(SEQ ID NO: 1007)
498	N-MYC	ATATTTTTGAGTTTTAA
	(SEQ ID NO: 60)	(SEQ ID NO: 1008)
499	N-MYC	TTTTAAAGCGTAGGTTGT
}	(SEQ ID NO: 60)	(SEQ ID NO: 1009)
500	N-MYC	TTTTAAAGTGTAGGTTGT
	(SEQ ID NO: 60)	(SEQ ID NO: 1010)
501	N-MYC	GGTGGATGCGGGGGTTT
	(SEQ ID NO: 60)	(SEQ ID NO: 1011)
502	N-MYC	GGTGGATGTGGGGGTTT
	(SEQ ID NO: 60)	(SEQ ID NO: 1012)
503	L-MYC	GGGTGGGGCGGGAGTAG
	(SEQ ID NO: 61)	
504	L-MYC	(SEQ ID NO: 1013) GGGTGGGGTGGGGAGTAG
	(SEQ ID NO: 61)	(SEQ ID NO: 1014)
505	L-MYC	TTTTAGTTCGGAGTGGGT
	(SEQ ID NO: 61)	
506	L-MYC	(SEQ ID NO: 1015) TTTTAGTTTGGAGTGGGT
	(SEQ ID NO: 61)	(SEQ ID NO: 1016)
507	L-MYC	AAGGAGTACGGTTTAGTT
	(SEQ ID NO: 61)	(SEQ ID NO: 1017)
508	L-MYC	AAGGAGTATGGTTTAGTT
	(SEQ ID NO: 61)	(SEQ ID NO: 1018)
509	L-MYC	AGTTTAGTCGGTTGGTAT
	(SEQ ID NO: 61)	(SEQ ID NO: 1019)
510	L-MYC	AGTTTAGTTGGTTAT
	(SEQ ID NO: 61)	(SEQ ID NO: 1020)
	(~2 < 10.01)	(υτιζ τη 140, 1020)

No:	Gene	Oligo:
511	L-MYC	TTTGTTTGCGTTTTTTAG
	(SEQ ID NO: 61)	(SEQ ID NO: 1021)
512	L-MYC	TTTGTTTGTGTTTTTAG
	(SEQ ID NO: 61)	(SEQ ID NO: 1022)
513	L-MYC	TTTTTAGCGTGTAATTT
	(SEQ ID NO: 61)	(SEQ ID NO: 1023)
514	L-MYC	TTTTTAGTGTGTAATTT
	(SEQ ID NO: 61)	(SEQ ID NO: 1024)
515	L-MYC	GGGGTTATCGGGGATTGA
	(SEQ ID NO: 61)	(SEQ ID NO: 1025)
516	L-MYC	GGGGTTATTGGGGATTGA
	(SEQ ID NO: 61)	(SEQ ID NO: 1026)
517	L-MYC	GGTGGGGTCGGTTTTATT
	(SEQ ID NO: 61)	(SEQ ID NO: 1027)
518	L-MYC	GGTGGGGTTGGTTTATT
	(SEQ ID NO: 61)	(SEQ ID NO: 1028)
519	L-MYC	GGAGAGGGCGGGTTTGT
	(SEQ ID NO: 61)	(SEQ ID NO: 1029)
520	L-MYC	GGAGAGGGTGGGGTTTGT
	(SEQ ID NO: 61)	(SEQ ID NO: 1030)
521	L-MYC	TTGAGGGTCGTTAGGTGG
	(SEQ ID NO: 61)	(SEQ ID NO: 1031)
522	L-MYC	TTGAGGGTTGTTAGGTGG
	(SEQ ID NO: 61)	(SEQ ID NO: 1032)
523	C-ABL	ATTAGAGTCGGGAGGGC
	(SEQ ID NO: 62)	(SEQ ID NO: 1033)
524	C-ABL	ATTAGAGTTGGGAGGGC
	(SEQ ID NO: 62)	(SEQ ID NO: 1034)
525	C-ABL	GGGGTATTCGGGTTTTTT
	(SEQ ID NO: 62)	(SEQ ID NO: 1035)
526	C-ABL	GGGGTATTTGGGTTTTTT
	(SEQ ID NO: 62)	(SEQ ID NO: 1036)
527	C-ABL	TTTGGTTTCGTTGTGGAG
	(SEQ ID NO: 62)	(SEQ ID NO: 1037)
528	C-ABL	TTTGGTTTTGTTGGAG
	(SEQ ID NO: 62)	(SEQ ID NO: 1038)
529	C-ABL	TTGAGTAGCGTTGGAGTC
	(SEQ ID NO: 62)	(SEQ ID NO: 1039)

No:	Gene	Oligo:
530	C-ABL	TTGAGTAGTGTTGGAGTC
	(SEQ ID NO: 62)	(SEQ ID NO: 1040)
531	C-ABL	GGGTTTTGCGGTTGÁGGA
	(SEQ ID NO: 62)	(SEQ ID NO: 1041)
532	C-ABL	GGGTTTTGTGGTTGÁGGA
	(SEQ ID NO: 62)	(SEQ ID NO: 1042)
533	C-ABL	AGTTTTTCGTTGTTTAC
	(SEQ ID NO: 62)	(SEQ ID NO: 1043)
534	C-ABL	AGTTTTTTGTTGTTTAC
	(SEQ ID NO: 62)	(SEQ ID NO: 1044)
535	C-ABL	GTTGTTTACGGTTTTTAT
	(SEQ ID NO: 62)	(SEQ ID NO: 1045)
536	C-ABL	GTTGTTTATGGTTTTTAT
	(SEQ ID NO: 62)	(SEQ ID NO: 1046)
537	C-ABL	TTTGAGGGCGTGGTGGT
	(SEQ ID NO: 62)	(SEQ ID NO: 1047)
538	C-ABL	TTTGAGGGTGGTGGT
	(SEQ ID NO: 62)	(SEQ ID NO: 1048)
539	C-ABL	GGGTGTTTCGGGTAGAGA
	(SEQ ID NO: 62)	(SEQ ID NO: 1049)
540	C-ABL	GGGTGTTTTGGGTAGAGA
	(SEQ ID NO: 62)	(SEQ ID NO: 1050)
541	C-ABL	GAGATTTTCGGGTTTGGG
	(SEQ ID NO: 62)	(SEQ ID NO: 1051)
542	C-ABL	GAGATTTTTGGGTTTGGG
	(SEQ ID NO: 62)	(SEQ ID NO: 1052)
543	ELK1	TTTGTTTTCGTTGAGTAG
	(SEQ ID NO: 63)	(SEQ ID NO: 1053)
544	ELK1	TTTGTTTTGAGTAG
	(SEQ ID NO: 63)	(SEQ ID NO: 1054)
545	ELK1	TTTATTTTCGTTTTTGGG
	(SEQ ID NO: 63)	(SEQ ID NO: 1241)
546	ELK1	TTTATTTTTTTTTGGG
	(SEQ ID NO: 63)	(SEQ ID NO: 1242)
547	ELK1	ATAGGTTTCGAGTTTTTT
	(SEQ ID NO: 63)	(SEQ ID NO: 1055)
548	ELK1	ATAGGTTTTGAGTTTTTT
	(SEQ ID NO: 63)	(SEQ ID NO: 1056)

No:	Gene	Oligo:
549	ELK1	TTATAGTTCGTTAGTTTT
	(SEQ ID NO: 63)	(SEQ ID NO: 1057)
550	ELK1	TTATAGTTTGTTAGTTTT
	(SEQ ID NO: 63)	(SEQ ID NO: 1058)
551	ELK1	TGGAGTATCGTGAAGGTC
	(SEQ ID NO: 63)	(SEQ ID NO: 1059)
552	ELK1	TGGAGTATTGTGAAGGTC
	(SEQ ID NO: 63)	(SEQ ID NO: 1060)
553	ELK1	TTTTTTACGTTAATTAC
	(SEQ ID NO: 63)	(SEQ ID NO: 1061)
554	ELK1	TTTTTTATGTTAATTAC
	(SEQ ID NO: 63)	(SEQ ID NO: 1062)
555	ELK1	GTTÄATTACGGTATÄGTT
	(SEQ ID NO: 63)	(SEQ ID NO: 1063)
556	ELK1	GTTAATTATGGTATAGTT
·	(SEQ ID NO: 63)	(SEQ ID NO: 1064)
557	ELK1	GAAGGGTTCGTTTTTAA
}	(SEQ ID NO: 63)	(SEQ ID NO: 1065)
558	ELK1	GAAGGGTTTGTTTTTAA
	(SEQ ID NO: 63)	(SEQ ID NO: 1066)
559	ELK1	TAATTTTCGAGGTTTTG
	(SEQ ID NO: 63)	(SEQ ID NO: 1067)
560	ELK1	TAATTTTTGAGGTTTTG
	(SEQ ID NO: 63)	(SEQ ID NO: 1068)
561	ELK1	ATTAATAGCGTTTTGGTT
	(SEQ ID NO: 63)	(SEQ ID NO: 1069)
562	ELK1	ATTAATAGTGTTTTGGTT
	(SEQ ID NO: 63)	(SEQ ID NO: 1070)
563	ELK1	GTTTAGTTCGATTTATTT
	(SEQ ID NO: 63)	(SEQ ID NO: 1071)
564	ELK1	GTTTAGTTTGATTTATTT
	(SEQ ID NO: 63)	(SEQ ID NO: 1072)
565	ELK1	GGTTGTTACGTTTATTTT
	(SEQ ID NO: 63)	(SEQ ID NO: 1073)
566	ELK1	GGTTGTTATGTTTATTTT
	(SEQ ID NO: 63)	(SEQ ID NO: 1074)
567	Tubulin	TACTTCAACGTCCTAATA
}	(SEQ ID NO: 64)	(SEQ ID NO: 1239)

No:	Gene	Oligo:
568	Tubulin	TACTTCAACATCCTAATA
	(SEQ ID NO: 64)	(SEQ ID NO: 1240)
569	Tubulin	CAATACTCCGTCCATTAA
	(SEQ ID NO: 64)	(SEQ ID NO: 1243)
570	Tubulin	CAATACTCCATCCATTAA
ī	(SEQ ID NO: 64)	(SEQ ID NO: 1244)
571	Tubulin	TAATCAACCGCAACCTCC
	(SEQ ID NO: 64)	(SEQ ID NO: 1075)
572	Tubulin	TAATCAACCACAACCTCC
	(SEQ ID NO: 64)	(SEQ ID NO: 1076)
573	Tubulin	CTTACCCCCGTCAATCAA
	(SEQ ID NO: 64)	(SEQ ID NO: 1077)
574	Tubulin	CTTACCCCCATCAATCAA
}	(SEQ ID NO: 64)	(SEQ ID NO: 1078)
575	Tubulin	GCTAATCTCGAAATÁAAC
	(SEQ ID NO: 64)	(SEQ ID NO: 1079)
576	Tubulin	GCTAATCTCAAAATÁAAC
	(SEQ ID NO: 64)	(SEQ ID NO: 1080)
577	CSF1	TTATAGAGCGTTAGTATT
	(SEQ ID NO: 65)	(SEQ ID NO: 1081)
578	CSF1	TTATAGAGTGTTAGTATT
	(SEQ ID NO: 65)	(SEQ ID NO: 1082)
579	CSF1	AAGTGTAGCGTAGAAGAT
	(SEQ ID NO: 65)	(SEQ ID NO: 1083)
580	CSF1	AAGTGTAGTGTAGAAGAT
	(SEQ ID NO: 65)	(SEQ ID NO: 1084)
581	CSF1	AGGGGAGGCGGGAAGG
	(SEQ ID NO: 65)	(SEQ ID NO: 1085)
582	CSF1	AGGGGAGGTGGGGAAGG
	(SEQ ID NO: 65)	(SEQ ID NO: 1086)
583	CSF1	GGGGAAGGCGGTTGAGTG
	(SEQ ID NO: 65)	(SEQ ID NO: 1087)
584	CSF1	GGGGAAGGTGGTTGAGTG
	(SEQ ID NO: 65)	(SEQ ID NO: 1088)
585	CSF1	GTGTTTGGCGTTA
	(SEQ ID NO: 65)	(SEQ ID NO: 1089)
586	CSF1	GTGTTTGGTGTTA
	(SEQ ID NO: 65)	(SEQ ID NO: 1090)

WO 02/077272 PCT/EP02/03401

72 .

S87 CD1R3	No:	Gene	Oligo:
588 CD1R3 TAACTATACAAACTACAA (SEQ ID NO: 66) (SEQ ID NO: 1092) 589 CD1R3 ACCCTCCTCGCTAAACTT (SEQ ID NO: 66) (SEQ ID NO: 1093) 590 CD1R3 ACCCTCCTCACTAAACTT (SEQ ID NO: 66) (SEQ ID NO: 1094) 591 CD1R3 CCAAACCCGACAACCCT (SEQ ID NO: 66) (SEQ ID NO: 1095) 592 CD1R3 CCAAACCCCAACACCCT (SEQ ID NO: 66) (SEQ ID NO: 1096) 593 CD1R3 TCCCTTTCCGCTAAATAC (SEQ ID NO: 66) (SEQ ID NO: 1097) 594 CD1R3 TCCCTTTCCACTAAATAC (SEQ ID NO: 66) (SEQ ID NO: 1098) 595 CD1R3 TCAACTCACGTCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1099) 596 CD1R3 TCAACTCACATCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1099) 597 CSNK2B AGGAGTTTTGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1236) 599 CSNK2B AGGAGTTTTGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1237) <td>587</td> <td>CD1R3</td> <td> </td>	587	CD1R3	
588 CD1R3 TAACTATACAAACTACAA (SEQ ID NO: 66) (SEQ ID NO: 1092) 589 CD1R3 ACCCTCCTCGCTAAACTT (SEQ ID NO: 66) (SEQ ID NO: 1093) 590 CD1R3 ACCCTCCTCACTAAACTT (SEQ ID NO: 66) (SEQ ID NO: 1094) 591 CD1R3 CCAAACCCGACAACCCT (SEQ ID NO: 66) (SEQ ID NO: 1095) 592 CD1R3 CCAAACCCCAACACCCT (SEQ ID NO: 66) (SEQ ID NO: 1096) 593 CD1R3 TCCCTTTCCGCTAAATAC (SEQ ID NO: 66) (SEQ ID NO: 1097) 594 CD1R3 TCCCTTTCCACTAAATAC (SEQ ID NO: 66) (SEQ ID NO: 1098) 595 CD1R3 TCAACTCACGTCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1099) 596 CD1R3 TCAACTCACATCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1099) 597 CSNK2B AGGAGTTTTGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1236) 599 CSNK2B AGGAGTTTTGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1237) <td></td> <td>(SEQ ID NO: 66)</td> <td>(SEO ID NO: 1091)</td>		(SEQ ID NO: 66)	(SEO ID NO: 1091)
589 CDIR3 ACCCTCCTCGCTAAACTT (SEQ ID NO: 66) (SEQ ID NO: 1093) 590 CDIR3 ACCCTCCTCACTAAACTT (SEQ ID NO: 66) (SEQ ID NO: 1094) 591 CDIR3 CCAAACCCCGACAACCCT (SEQ ID NO: 66) (SEQ ID NO: 1095) 592 CDIR3 CCAAACCCCAACAACCCT (SEQ ID NO: 66) (SEQ ID NO: 1096) 593 CDIR3 TCCCTTTCCGCTAAATAC (SEQ ID NO: 66) (SEQ ID NO: 1097) 594 CDIR3 TCCCTTTCCACTAAATAC (SEQ ID NO: 66) (SEQ ID NO: 1098) 595 CDIR3 TCAACTCACGTCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1099) 596 CDIR3 TCAACTCACATCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1100) 597 CSNK2B AGGAGTTTCGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1236) 599 CSNK2B AGGAGTTTTGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1237) 600 CSNK2B ATTTCCTCCGAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1238) <	588		
589 CDIR3 ACCCTCCTCGCTAAACTT (SEQ ID NO: 66) (SEQ ID NO: 1093) 590 CDIR3 ACCCTCCTCACTAAACTT (SEQ ID NO: 66) (SEQ ID NO: 1094) 591 CDIR3 CCAAACCCCGACAACCCT (SEQ ID NO: 66) (SEQ ID NO: 1095) 592 CDIR3 CCAAACCCCAACAACCCT (SEQ ID NO: 66) (SEQ ID NO: 1096) 593 CDIR3 TCCCTTTCCGCTAAATAC (SEQ ID NO: 66) (SEQ ID NO: 1097) 594 CDIR3 TCCCTTTCCACTAAATAC (SEQ ID NO: 66) (SEQ ID NO: 1098) 595 CDIR3 TCAACTCACGTCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1099) 596 CDIR3 TCAACTCACATCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1100) 597 CSNK2B AGGAGTTTCGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1236) 599 CSNK2B AGGAGTTTTGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1237) 600 CSNK2B ATTTCCTCCGAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1238) <		(SEO ID NO: 66)	(SEO ID NO: 1092)
590 CD1R3 ACCCTCCTCACTAAACTT (SEQ ID NO: 66) (SEQ ID NO: 1094) 591 CD1R3 CCAAACCCCGACAACCCT (SEQ ID NO: 66) (SEQ ID NO: 1095) 592 CD1R3 CCAAACCCCAACAACCCT (SEQ ID NO: 66) (SEQ ID NO: 1096) 593 CD1R3 TCCCTTTCCGCTAAATAC (SEQ ID NO: 66) (SEQ ID NO: 1097) 594 CD1R3 TCCCTTTCCACTAAATAC (SEQ ID NO: 1098) (SEQ ID NO: 1098) 595 CD1R3 TCAACTCACGTCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1099) 596 CD1R3 TCAACTCACATCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1100) 597 CSNK2B AGGAGTTTCGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1235) 598 CSNK2B AGGAGTTTTGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1236) 599 CSNK2B ATTTCCTCCGAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1237) 600 CSNK2B ATTTCCTCCAAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1238)	589		
590 CD1R3 ACCCTCCTCACTAAACTT (SEQ ID NO: 66) (SEQ ID NO: 1094) 591 CD1R3 CCAAACCCCGACAACCCT (SEQ ID NO: 66) (SEQ ID NO: 1095) 592 CD1R3 CCAAACCCCAACAACCCT (SEQ ID NO: 66) (SEQ ID NO: 1096) 593 CD1R3 TCCCTTTCCGCTAAATAC (SEQ ID NO: 66) (SEQ ID NO: 1097) 594 CD1R3 TCCCTTTCCACTAAATAC (SEQ ID NO: 1098) (SEQ ID NO: 1098) 595 CD1R3 TCAACTCACGTCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1099) 596 CD1R3 TCAACTCACATCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1100) 597 CSNK2B AGGAGTTTCGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1235) 598 CSNK2B AGGAGTTTTGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1236) 599 CSNK2B ATTTCCTCCGAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1237) 600 CSNK2B ATTTCCTCCAAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1238)	i	(SEO ID NO: 66)	(SEO ID NO: 1093)
591 CD1R3 CCAAACCCGGACAACCCT (SEQ ID NO: 66) (SEQ ID NO: 1095) 592 CD1R3 CCAAACCCCAACAACCCT (SEQ ID NO: 66) (SEQ ID NO: 1096) 593 CD1R3 TCCCTTTCCGCTAAATAC (SEQ ID NO: 66) (SEQ ID NO: 1097) 594 CD1R3 TCCCTTTCCACTAAATAC (SEQ ID NO: 66) (SEQ ID NO: 1098) 595 CD1R3 TCAACTCACGTCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1099) 596 CD1R3 TCAACTCACATCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1100) 597 CSNK2B AGGAGTTTCGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1235) 598 CSNK2B AGGAGTTTTGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1236) 599 CSNK2B ATTTCCTCCGAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1237) 600 CSNK2B ATTTCCTCCAAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1238) 601 CSNK2B GATTGATTCGTATTAGGG (SEQ ID NO: 67) (SEQ ID NO: 1101)	590		
591 CD1R3 CCAAACCCGGACAACCCT (SEQ ID NO: 66) (SEQ ID NO: 1095) 592 CD1R3 CCAAACCCCAACAACCCT (SEQ ID NO: 66) (SEQ ID NO: 1096) 593 CD1R3 TCCCTTTCCGCTAAATAC (SEQ ID NO: 66) (SEQ ID NO: 1097) 594 CD1R3 TCCCTTTCCACTAAATAC (SEQ ID NO: 66) (SEQ ID NO: 1098) 595 CD1R3 TCAACTCACGTCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1099) 596 CD1R3 TCAACTCACATCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1100) 597 CSNK2B AGGAGTTTCGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1235) 598 CSNK2B AGGAGTTTTGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1236) 599 CSNK2B ATTTCCTCCGAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1237) 600 CSNK2B ATTTCCTCCAAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1238) 601 CSNK2B GATTGATTCGTATTAGGG (SEQ ID NO: 67) (SEQ ID NO: 1101)		(SEO ID NO: 66)	(SEO ID NO: 1094)
592 CD1R3 CCAAACCCCAACAACCCT (SEQ ID NO: 66) (SEQ ID NO: 1096) 593 CD1R3 TCCCTTTCCGCTAAATAC (SEQ ID NO: 66) (SEQ ID NO: 1097) 594 CD1R3 TCCCTTTCCACTAAATAC (SEQ ID NO: 66) (SEQ ID NO: 1098) 595 CD1R3 TCAACTCACGTCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1099) 596 CD1R3 TCAACTCACATCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1100) 597 CSNK2B AGGAGTTTCGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1235) 598 CSNK2B AGGAGTTTTGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1236) 599 CSNK2B ATTTCCTCCGAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1237) 600 CSNK2B ATTTCCTCCAAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1238) 601 CSNK2B GATTGATTCGTATTAGGG (SEQ ID NO: 67) (SEQ ID NO: 1101)	591		
592 CD1R3 CCAAACCCCAACAACCCT (SEQ ID NO: 66) (SEQ ID NO: 1096) 593 CD1R3 TCCCTTTCCGCTAAATAC (SEQ ID NO: 66) (SEQ ID NO: 1097) 594 CD1R3 TCCCTTTCCACTAAATAC (SEQ ID NO: 66) (SEQ ID NO: 1098) 595 CD1R3 TCAACTCACGTCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1099) 596 CD1R3 TCAACTCACATCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1100) 597 CSNK2B AGGAGTTTCGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1235) 598 CSNK2B AGGAGTTTTGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1236) 599 CSNK2B ATTTCCTCCGAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1237) 600 CSNK2B ATTTCCTCCAAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1238) 601 CSNK2B GATTGATTCGTATTAGGG (SEQ ID NO: 67) (SEQ ID NO: 1101)		(SEO ID NO: 66)	(SEO ID NO: 1095)
593 CD1R3 TCCCTTTCCGCTAAATAC (SEQ ID NO: 66) (SEQ ID NO: 1097) 594 CD1R3 TCCCTTTCCACTAAATAC (SEQ ID NO: 66) (SEQ ID NO: 1098) 595 CD1R3 TCAACTCACGTCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1099) 596 CD1R3 TCAACTCACATCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1100) 597 CSNK2B AGGAGTTTCGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1235) 598 CSNK2B AGGAGTTTTGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1236) 599 CSNK2B ATTTCCTCCGAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1237) 600 CSNK2B ATTTCCTCCAAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1238) 601 CSNK2B GATTGATTCGTATTAGGG (SEQ ID NO: 67) (SEQ ID NO: 1101)	592		
593 CD1R3 TCCCTTTCCGCTAAATAC (SEQ ID NO: 66) (SEQ ID NO: 1097) 594 CD1R3 TCCCTTTCCACTAAATAC (SEQ ID NO: 66) (SEQ ID NO: 1098) 595 CD1R3 TCAACTCACGTCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1099) 596 CD1R3 TCAACTCACATCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1100) 597 CSNK2B AGGAGTTTCGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1235) 598 CSNK2B AGGAGTTTTGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1236) 599 CSNK2B ATTTCCTCCGAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1237) 600 CSNK2B ATTTCCTCCAAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1238) 601 CSNK2B GATTGATTCGTATTAGGG (SEQ ID NO: 67) (SEQ ID NO: 1101)		(SEO ID NO: 66)	(SEO ID NO: 1096)
594 CD1R3 TCCCTTTCCACTAAATAC (SEQ ID NO: 66) (SEQ ID NO: 1098) 595 CD1R3 TCAACTCACGTCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1099) 596 CD1R3 TCAACTCACATCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1100) 597 CSNK2B AGGAGTTTCGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1235) 598 CSNK2B AGGAGTTTTGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1236) 599 CSNK2B ATTTCCTCCGAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1237) 600 CSNK2B ATTTCCTCCAAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1238) 601 CSNK2B GATTGATTCGTATTAGGG (SEQ ID NO: 67) (SEQ ID NO: 1101)	593		
594 CD1R3 TCCCTTTCCACTAAATAC (SEQ ID NO: 66) (SEQ ID NO: 1098) 595 CD1R3 TCAACTCACGTCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1099) 596 CD1R3 TCAACTCACATCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1100) 597 CSNK2B AGGAGTTTCGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1235) 598 CSNK2B AGGAGTTTTGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1236) 599 CSNK2B ATTTCCTCCGAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1237) 600 CSNK2B ATTTCCTCCAAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1238) 601 CSNK2B GATTGATTCGTATTAGGG (SEQ ID NO: 67) (SEQ ID NO: 1101)		(SEO ID NO: 66)	(SEO ID NO: 1097)
595 CD1R3 TCAACTCACGTCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1099) 596 CD1R3 TCAACTCACATCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1100) 597 CSNK2B AGGAGTTTCGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1235) 598 CSNK2B AGGAGTTTTGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1236) 599 CSNK2B ATTTCCTCCGAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1237) 600 CSNK2B ATTTCCTCCAAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1238) 601 CSNK2B GATTGATTCGTATTAGGG (SEQ ID NO: 67) (SEQ ID NO: 1101)	594		
595 CD1R3 TCAACTCACGTCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1099) 596 CD1R3 TCAACTCACATCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1100) 597 CSNK2B AGGAGTTTCGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1235) 598 CSNK2B AGGAGTTTTGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1236) 599 CSNK2B ATTTCCTCCGAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1237) 600 CSNK2B ATTTCCTCCAAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1238) 601 CSNK2B GATTGATTCGTATTAGGG (SEQ ID NO: 67) (SEQ ID NO: 1101)		(SEO ID NO: 66)	(SEO ID NO: 1098)
596 CD1R3 TCAACTCACATCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1100) 597 CSNK2B AGGAGTTTCGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1235) 598 CSNK2B AGGAGTTTTGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1236) 599 CSNK2B ATTTCCTCCGAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1237) 600 CSNK2B ATTTCCTCCAAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1238) 601 CSNK2B GATTGATTCGTATTAGGG (SEQ ID NO: 67) (SEQ ID NO: 1101)	595		
596 CD1R3 TCAACTCACATCCCTTTC (SEQ ID NO: 66) (SEQ ID NO: 1100) 597 CSNK2B AGGAGTTTCGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1235) 598 CSNK2B AGGAGTTTTGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1236) 599 CSNK2B ATTTCCTCCGAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1237) 600 CSNK2B ATTTCCTCCAAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1238) 601 CSNK2B GATTGATTCGTATTAGGG (SEQ ID NO: 67) (SEQ ID NO: 1101)		(SEO ID NO: 66)	(SEO ID NO: 1099)
597 CSNK2B AGGAGTTTCGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1235) 598 CSNK2B AGGAGTTTTGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1236) 599 CSNK2B ATTTCCTCCGAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1237) 600 CSNK2B ATTTCCTCCAAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1238) 601 CSNK2B GATTGATTCGTATTAGGG (SEQ ID NO: 67) (SEQ ID NO: 1101)	596		
597 CSNK2B AGGAGTTTCGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1235) 598 CSNK2B AGGAGTTTTGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1236) 599 CSNK2B ATTTCCTCCGAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1237) 600 CSNK2B ATTTCCTCCAAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1238) 601 CSNK2B GATTGATTCGTATTAGGG (SEQ ID NO: 67) (SEQ ID NO: 1101)		(SEO ID NO: 66)	(SEO ID NO: 1100)
598 CSNK2B AGGAGTTTTGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1236) 599 CSNK2B ATTTCCTCCGAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1237) 600 CSNK2B ATTTCCTCCAAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1238) 601 CSNK2B GATTGATTCGTATTAGGG (SEQ ID NO: 67) (SEQ ID NO: 1101)	597		
598 CSNK2B AGGAGTTTTGGAGGAAAT (SEQ ID NO: 67) (SEQ ID NO: 1236) 599 CSNK2B ATTTCCTCCGAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1237) 600 CSNK2B ATTTCCTCCAAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1238) 601 CSNK2B GATTGATTCGTATTAGGG (SEQ ID NO: 67) (SEQ ID NO: 1101)		(SEO ID NO: 67)	(SEO ID NO: 1235)
599 CSNK2B ATTTCCTCCGAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1237) 600 CSNK2B ATTTCCTCCAAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1238) 601 CSNK2B GATTGATTCGTATTAGGG (SEQ ID NO: 67) (SEQ ID NO: 1101)	598	 	
599 CSNK2B ATTTCCTCCGAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1237) 600 CSNK2B ATTTCCTCCAAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1238) 601 CSNK2B GATTGATTCGTATTAGGG (SEQ ID NO: 67) (SEQ ID NO: 1101)		(SEQ ID NO: 67)	(SEO ID NO: 1236)
600 CSNK2B ATTTCCTCCAAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1238) 601 CSNK2B GATTGATTCGTATTAGGG (SEQ ID NO: 67) (SEQ ID NO: 1101)	599		
600 CSNK2B ATTTCCTCCAAAACTCCT (SEQ ID NO: 67) (SEQ ID NO: 1238) 601 CSNK2B GATTGATTCGTATTAGGG (SEQ ID NO: 67) (SEQ ID NO: 1101)		(SEQ ID NO: 67)	(SEO ID NO: 1237)
601 CSNK2B GATTGATTCGTATTAGGG (SEQ ID NO: 67) (SEQ ID NO: 1101)	600	CSNK2B	
(SEQ ID NO: 67) (SEQ ID NO: 1101)		(SEQ ID NO: 67)	(SEQ ID NO: 1238)
	601	CSNK2B	GATTGATTCGTATTÁGGG
602 CSNK2B GATTGATTTGTATTAGGG		(SEQ ID NO: 67)	(SEQ ID NO: 1101)
	602	CSNK2B	GATTGATTTGTATTAGGG
(SEQ ID NO: 67) (SEQ ID NO: 1102)		(SEQ ID NO: 67)	(SEQ ID NO: 1102)
603 CSNK2B CCCTAATACGAATCAATC	603	CSNK2B	
(SEQ ID NO: 67) (SEQ ID NO: 1103)			(SEQ ID NO: 1103)
604 CSNK2B CCCTAATACAAATCAATC	604	CSNK2B	CCCTAATACAAATCAATC
(SEQ ID NO: 67) (SEQ ID NO: 1104)		(SEQ ID NO: 67)	(SEQ ID NO: 1104)
605 CSNK2B GGATTTAGCGGATTGGTC	605	CSNK2B	
(SEQ ID NO: 67) (SEQ ID NO: 1105)		(SEQ ID NO: 67)	(SEQ ID NO: 1105)

No:	Gene	Oligo:
606	CSNK2B	GGATTTAGTGGATTGGTC
	(SEQ ID NO: 67)	(SEQ ID NO: 1106)
607	CSNK2B	GACCAATCCGCTAAATCC
	(SEQ ID NO: 67)	(SEQ ID NO: 1107)
608	CSNK2B	GACCAATCCACTAAATCC
	(SEQ ID NO: 67)	(SEQ ID NO: 1108)
609	CSNK2B	GGATTGGTCGAGGATAGG
	(SEQ ID NO: 67)	(SEQ ID NO: 1109)
610	CSNK2B	GGATTGGTTGAGGATAGG
	(SEQ ID NO: 67)	(SEQ ID NO: 1110)
611	CSNK2B	CCTATCCTCGACCAATCC
	(SEQ ID NO: 67)	(SEQ ID NO: 1111)
612	CSNK2B	CCTATCCTCAACCAATCC
	(SEQ ID NO: 67)	(SEQ ID NO: 1112)
613	CSNK2B	GAGGGGAACGTGAGGAGA
	(SEQ ID NO: 67)	(SEQ ID NO: 1113)
614	CSNK2B	GAGGGAATGTGAGGAGA
	(SEQ ID NO: 67)	(SEQ ID NO: 1114)
615	CSNK2B	TCTCCTCACGTTCCCCTC
	(SEQ ID NO: 67)	(SEQ ID NO: 1115)
616	CSNK2B	TCTCCTCACATTCCCCTC
	(SEQ ID NO: 67)	(SEQ ID NO: 1116)
617	CSNK2B	TAGGTTAGCGTATTGGGA
	(SEQ ID NO: 67)	(SEQ ID NO: 1117)
618	CSNK2B	TAGGTTAGTGTATTGGGA
	(SEQ ID NO: 67)	(SEQ ID NO: 1118)
619	CSNK2B	TCCCAATACGCTAACCTA
	(SEQ ID NO: 67)	(SEQ ID NO: 1119)
620	CSNK2B	TCCCAATACACTAACCTA
	(SEQ ID NO: 67)	(SEQ ID NO: 1120)
621	CSNK2B	TATTGGGACGTTTTAGTT
	(SEQ ID NO: 67)	(SEQ ID NO: 1121)
622	CSNK2B	TATTGGGATGTTTTAGTT
	(SEQ ID NO: 67)	(SEQ ID NO: 1122)
623	CSNK2B	AACTAAAACGTCCCAATA
	(SEQ ID NO: 67)	(SEQ ID NO: 1123)
624	CSNK2B	AACTAAAACATCCCAATA
	(SEQ ID NO: 67)	(SEQ ID NO: 1124)
L	(BEQ 11 110.07)	(0000 10.1127)

No:	Gene	Oligo:
625	Me491/TD63	TTGTTATTCGGTTTTAGT
	(SEQ ID NO: 68)	(SEQ ID NO: 1125)
626	Me491/TD63	TTGTTATTTGGTTTTAGT
	(SEQ ID NO: 68)	(SEQ ID NO: 1126)
627	Me491/TD63	AGTTGGTACGGAGGTATA
	(SEQ ID NO: 68)	(SEQ ID NO: 1127)
628	Me491/TD63	AGTTGGTATGGAGGTATA
	(SEQ ID NO: 68)	(SEQ ID NO: 1128)
629	Me491/TD63	AGGTATAGCGTTAGAGGG
	(SEQ ID NO: 68)	(SEQ ID NO: 1129)
630	Me491/TD63	AGGTATAGTGTTAGAGGG
	(SEQ ID NO: 68)	(SEQ ID NO: 1130)
631	Me491/TD63	AGGGGTTACGTGGTTTTG
	(SEQ ID NO: 68)	(SEQ ID NO: 1131)
632	Me491/TD63	AGGGGTTATGTGGTTTTG
	(SEQ ID NO: 68)	(SEQ ID NO: 1132)
633	Me491/TD63	TGTATAGTCGGATGGGGA
	(SEQ ID NO: 68)	(SEQ ID NO: 1133)
634	Me491/TD63	TGTATAGTTGGATGGGGA
	(SEQ ID NO: 68)	(SEQ ID NO: 1134)
635	AR	ACCTCTCTCGAAATAACA
	(SEQ ID NO: 69)	(SEQ ID NO: 1135)
636	AR	ACCTCTCAAAATAACA
	1	
637	(SEQ ID NO: 69) AR	(SEQ ID NO: 1136) CTCTAAAACGCAACCTCT
057]	
638	(SEQ ID NO: 69)	(SEQ ID NO: 1137) CTCTAAAACACAACCTCT
639	(SEQ ID NO: 69)	(SEQ ID NO: 1138) AACTACTACGACAACCCC
027		
640	(SEQ ID NO: 69) AR	(SEQ ID NO: 1139) AACTACTACAACACCCC
010		
641	(SEQ ID NO: 69) AR	(SEQ ID NO: 1140) ACTAACCTCGCTCAAAAT
O II		
642	(SEQ ID NO: 69) AR	(SEQ ID NO: 1141) ACTAACCTCACTCAAAAT
U-T4		
643	(SEQ ID NO: 69) AR	(SEQ ID NO: 1142) ACTACCTTCGAATACTAC
U+3		
	(SEQ ID NO: 69)	(SEQ ID NO: 1257)

No:	Gene	Oligo:
644	AR	ACTACCTTCAAATACTAC
	(SEQ ID NO: 69)	(SEQ ID NO: 1258)
645	CDK 4	TATGGGGTCGTAGGAATC
	(SEQ ID NO: 70)	(SEQ ID NO: 1143)
646	CDK 4	TATGGGGTTGTAGGAATC
	(SEQ ID NO: 70)	(SEQ ID NO: 1144)
647	CDK 4	GATTCCTACGACCCCATA
	(SEQ ID NO: 70)	(SEQ ID NO: 1145)
648	CDK 4	GATTCCTACAACCCCATA
	(SEQ ID NO: 70)	(SEQ ID NO: 1146)
649	CDK 4	GGAAGGGTCGTTTAAGGG
		(SEQ ID NO: 1147)
650	(SEQ ID NO: 70) CDK 4	GGAAGGGTTGTTTAAGGG
		(SEQ ID NO: 1148)
651	(SEQ ID NO: 70) CDK 4	CCCTTAAACGACCCTTCC
031		(SEQ ID NO: 1149)
652	(SEQ ID NO: 70) CDK 4	CCCTTAAACAACCCTTCC
032		
653	(SEQ ID NO: 70) CDK 4	(SEQ ID NO: 1150) TTTAAGGGCGGGAAGTGG
033		
654	(SEQ ID NO: 70) CDK 4	(SEQ ID NO: 1247) TTTAAGGGTGGGAAGTGG
054		
655	(SEQ ID NO: 70) CDK 4	(SEQ ID NO: 1248) CCACTTCCCGCCCTTAAA
033		
656	(SEQ ID NO: 70)	(SEQ ID NO: 1249) CCACTTCCCACCCTTAAA
050		
657	(SEQ ID NO: 70) CDK 4	(SEQ ID NO: 1250) AGGATTTTCGATGTAAGG
05/		
658	(SEQ ID NO: 70) CDK 4	(SEQ ID NO: 1151) AGGATTTTTGATGTAAGG
860		
650	(SEQ ID NO: 70) CDK 4	(SEQ ID NO: 1152) CCTTACATCGAAAATCCT
659		•
660	(SEQ ID NO: 70) CDK 4	(SEQ ID NO: 1153) CCTTACATCAAAAATCCT
660		
661	(SEQ ID NO: 70) CDK 4	(SEQ ID NO: 1154) GGGTTTTACGTGGTTGGA
661		
660	(SEQ ID NO: 70)	(SEQ ID NO: 1231) GGGTTTTATGTGGTTGGA
662	CDK 4	·
	(SEQ ID NO: 70)	(SEQ ID NO: 1232)

No:	Gene	Oligo:
663	CDK 4	TCCAACCACGTAAAACCC
	(SEQ ID NO: 70)	(SEQ ID NO: 1233)
664	CDK 4	TCCAACCACATAAAACCC
	(SEQ ID NO: 70)	(SEQ ID NO: 1234)
665	Humos	CCTTACTACGTTAAACTC
	(SEQ ID NO: 71)	(SEQ ID NO: 1155)
666	Humos	CCTTACTACATTAAACTC
	(SEQ ID NO: 71)	(SEQ ID NO: 1156)
667	Humos	GTTACCACCGAACTCCAT
	(SEQ ID NO: 71)	(SEQ ID NO: 1245)
668	Humos	GTTACCACCAAACTCCAT
	(SEQ ID NO: 71)	(SEQ ID NO: 1246)
669	Humos	CTCCCTACGTCCCCTC
	(SEQ ID NO: 71)	(SEQ ID NO: 1253)
670	Humos	CTCCCCTACATCCCCCTC
	(SEQ ID NO: 71)	(SEQ ID NO: 1254)
671	Humos	CTCCAATACGACAATAAA
	(SEQ ID NO: 71)	(SEQ ID NO: 1157)
672	Humos	CTCCAATACAACAATAAA
0,2	(SEQ ID NO: 71)	(SEQ ID NO: 1158)
673	Humos	AAACAAACCGTTCACAAC
0,5		(SEQ ID NO: 1251)
674	(SEQ ID NO: 71) Humos	AAACAAACCATTCACAAC
071		(SEQ ID NO: 1252)
675	(SEQ ID NO: 71) CDC25A	GTGTAGGTCGGTTTGGTT
075		(SEQ ID NO: 1159)
676	(SEQ ID NO: 72) CDC25A	GEQ ID NO. 1139) GTGTAGGTTGGTT
070		
677	(SEQ ID NO: 72) CDC25A	(SEQ ID NO: 1160) AACCAAACCGACCTACAC
077		
678	(SEQ ID NO: 72) CDC25A	(SEQ ID NO: 1161) AACCAAACCAACCTACAC
078		
679	(SEQ ID NO: 72) CDC25A	(SEQ ID NO: 1162) TTGTTATTCGGAGTTGGG
013		
680	(SEQ ID NO: 72) CDC25A	(SEQ ID NO: 1163) TTGTTATTTGGAGTTGGG
OQU		
(01	(SEQ ID NO: 72)	(SEQ ID NO: 1164) CCCAACTCCGAATAACAA
681	CDC25A	
	(SEQ ID NO: 72)	(SEQ ID NO: 1165)

No:	Gene	Oligo:
682	CDC25A	CCCAACTCCAAATAACAA
	(SEQ ID NO: 72)	(SEQ ID NO: 1166)
683	CDC25A	TGGGTAAGCGGGTGGGAG
	(SEQ ID NO: 72)	(SEQ ID NO: 1167)
684	CDC25A	TGGGTAAGTGGGTGGGAG
	(SEQ ID NO: 72)	(SEQ ID NO: 1168)
685	CDC25A	CTCCCACCCGCTTACCCA
ĺ	(SEQ ID NO: 72)	(SEQ ID NO: 1169)
686	CDC25A	CTCCCACCCACTTACCCA
ļ	(SEQ ID NO: 72)	(SEQ ID NO: 1170)
687	CDC25A	GAGAATAGCGAAGATAGC
	(SEQ ID NO: 72)	(SEQ ID NO: 1171)
688	CDC25A	GAGAATAGTGAAGATAGC
	(SEQ ID NO: 72)	
689	CDC25A	(SEQ ID NO: 1172) GCTATCTTCGCTATTCTC
003		
690	(SEQ ID NO: 72) CDC25A	(SEQ ID NO: 1173) GCTATCTTCACTATTCTC
050		
691	(SEQ ID NO: 72) CDC25A	(SEQ ID NO: 1174) TATTGAGTCGTTATTATC
091		
692	(SEQ ID NO: 72) CDC25A	(SEQ ID NO: 1175)
092		TATTGAGTTGTTATTATC
602	(SEQ ID NO: 72)	(SEQ ID NO: 1176)
693	CDC25A	GATAATAACGACTCAATA
604	(SEQ ID NO: 72)	(SEQ ID NO: 1177)
694	CDC25A	GATAATAACAACTCAATA
	(SEQ ID NO: 72)	(SEQ ID NO: 1178)
695	CMYCex3	GAAGAAATCGATGTTGTT
	(SEQ ID NO: 73)	(SEQ ID NO: 1179)
696	CMYCex3	GAAGAAATTGATGTT
	(SEQ ID NO: 73)	(SEQ ID NO: 1180)
697	CMYCex3	AGGTGTTACGTTTTTATA
	(SEQ ID NO: 73)	(SEQ ID NO: 1181)
698	CMYCex3	AGGTGTTATGTTTTATA
	(SEQ ID NO: 73)	(SEQ ID NO: 1182)
699	CMYCex3	TTTTATTCGGAAGGATT
	(SEQ ID NO: 73)	(SEQ ID NO: 1183)
700	CMYCex3	TTTTATTTGGAAGGATT
	(SEQ ID NO: 73)	(SEQ ID NO: 1184)

No:	Gene	Oligo:
701	CMYCex3	GTAATAATCGAAAATGTA
	(SEQ ID NO: 73)	(SEQ ID NO: 1185)
702	CMYCex3	GTAATAATTGAAAATGTA
	(SEQ ID NO: 73)	(SEQ ID NO: 1186)
703	CMYCex3	TTAAGAGGCGAATATATA
	(SEQ ID NO: 73)	(SEQ ID NO: 1187)
704	CMYCex3	TTAAGAGGTGAATATATA
	(SEQ ID NO: 73)	(SEQ ID NO: 1188)
705	CMYCex3	ATATATAACGTTTTGGAG
}	(SEQ ID NO: 73)	(SEQ ID NO: 1189)
706	CMYCex3	ATATATAATGTTTTGGAG
	(SEQ ID NO: 73)	(SEQ ID NO: 1190)
707	CMYCex3	TTTTGGAGCGTTAGAGGA
	(SEQ ID NO: 73)	(SEQ ID NO: 1191)
708	CMYCex3	TTTTGGAGTGTTAGAGGA
	(SEQ ID NO: 73)	(SEQ ID NO: 1192)
709	CMYCex3	AGGAGGAACGAGTTAAAA
	(SEQ ID NO: 73)	(SEQ ID NO: 1193)
710	CMYCex3	AGGAGGAATGAGTTAAAA
	(SEQ ID NO: 73)	(SEQ ID NO: 1194)
711	CMYCex3	AGTTAAAACGGAGTTTTT
	(SEQ ID NO: 73)	(SEQ ID NO: 1195)
712	CMYCex3	AGTTAAAATGGAGTTTTT
	(SEQ ID NO: 73)	(SEQ ID NO: 1196)
713	CMYCex3	TTGTTTTGCGTGATTAGA
	(SEQ ID NO: 73)	(SEQ ID NO: 1197)
714	CMYCex3	TTGTTTTGTGTGATTAGA
	(SEQ ID NO: 73)	(SEQ ID NO: 1198)
715	CMYCex3	TTAGATTTCGGAGTTGGA
	(SEQ ID NO: 73)	(SEQ ID NO: 1199)
716	CMYCex3	TTAGATTTTGGAGTTGGA
	(SEQ ID NO: 73)	(SEQ ID NO: 1200)
717	CMYCex3	ATTITGTTCGTTTAAGTA
	(SEQ ID NO: 73)	(SEQ ID NO: 1201)
718	CMYCex3	ATTTTGTTTGTTTAAGTA
	(SEQ ID NO: 73)	(SEQ ID NO: 1202)
719	CMYCex3	AATAGTTACGGAATTTTT
	(SEQ ID NO: 73)	(SEQ ID NO: 1203)
L	(900 10.73)	(ppd m 140, 1503)

(SEQ ID NO: 1210)

79

No:	Gene	Oligo:
720	CMYCex3	AATAGTTATGGAATTTTT
	(SEQ ID NO: 73)	(SEQ ID NO: 1204)
721	CMYCex3	TTTTTGTGCGTAAGGAAA
	(SEQ ID NO: 73)	(SEQ ID NO: 1205)
722	CMYCex3	TTTTTGTGTGTAAGGAAA
	(SEQ ID NO: 73)	(SEQ ID NO: 1206)
723	CMYCex3	AAGGAAAACGATTTTTT
	(SEQ ID NO: 73)	(SEQ ID NO: 1207)
724	CMYCex3	AAGGAAAATGATTTTTT
	(SEQ ID NO: 73)	(SEQ ID NO: 1208)
725	CMYCex3	TGAGTAATCGTTTATGAA
	(SEQ ID NO: 73)	(SEQ ID NO: 1209)
726	CMYCex3	TGAGTAATTGTTTATGAA

Table 4: Hybridisation oligonucleotides used in the differentiation between healthy and leukemia cells (Figure 3)

(SEQ ID NO: 73)

No:	Gene	Oligo:
1	Humos	GTTACCACCGAACTCCAT
	(SEQ ID NO: 71)	(SEQ ID NO: 1245)
2	Humos	GTTACCACCAAACTCCAT .
	(SEQ ID NO: 71)	(SEQ ID NO: 1246)
3	Humos	GTTACCACCGAACTCCAT
	(SEQ ID NO: 71)	(SEQ ID NO: 1245)
4	Humos	GTTACCACCAAACTCCAT
	(SEQ ID NO: 71)	(SEQ ID NO: 1246)
5	CDK 4	TTTAAGGGCGGAAGTGG
	(SEQ ID NO: 70)	(SEQ ID NO: 1247)
6	CDK 4	TTTAAGGGTGGGAAGTGG
	(SEQ ID NO: 70)	(SEQ ID NO: 1248)
7	CDK 4	CCACTTCCCGCCCTTAAA
	(SEQ ID NO: 70)	(SEQ ID NO: 1249)
8	CDK 4	CCACTTCCCACCCTTAAA
	(SEQ ID NO: 70)	(SEQ ID NO: 1250)

No:	Gene	Oligo:
9	CDK 4	TTTAAGGGCGGAAGTGG
	(SEQ ID NO: 70)	(SEQ ID NO: 1247)
10	CDK 4	TTTAAGGGTGGGAAGTGG
	(SEQ ID NO: 70)	(SEQ ID NO: 1248)
11	CDK 4	CCACTTCCCGCCCTTAAA
	(SEQ ID NO: 70)	(SEQ ID NO: 1249)
12	CDK 4	CCACTTCCCACCCTTAAA
	(SEQ ID NO: 70)	(SEQ ID NO: 1250)
13	Humos	AAACAAACCGTTCACAAC
	(SEQ ID NO: 71)	(SEQ ID NO: 1251)
14	Humos	AAACAAACCATTCACAAC
	(SEQ ID NO: 71)	(SEQ ID NO: 1252)
15	Humos	AAACAAACCGTTCACAAC
	(SEQ ID NO: 71)	(SEQ ID NO: 1251)
16	Humos	AAACAAACCATTCACAAC
	(SEQ ID NO: 71)	(SEQ ID NO: 1252)
17	Humos	CTCCCTACGTCCCCTC
	(SEQ ID NO: 71)	(SEQ ID NO: 1253)
18	Humos	CTCCCTACATCCCCCTC
	(SEQ ID NO: 71)	(SEQ ID NO: 1254)
1.9	Humos	CTCCCTACGTCCCCTC
	(SEQ ID NO: 71)	(SEQ ID NO: 1253)
20	Humos	CTCCCTACATCCCCCTC
	(SEQ ID NO: 71)	(SEQ ID NO: 1254)
21	AR	GTAGTATTCGAAGGTAGT
		(SEQ ID NO: 1255)
22	(SEQ ID NO: 69) AR	GTAGTATTTGAAGGTAGT
	(SEQ ID NO: 69)	(SEQ ID NO: 1256)
23	AR	ACTACCTTCGAATACTAC
	(SEQ ID NO: 69)	(SEQ ID NO: 1257)
24	AR	ACTACCTTCAAATACTAC
-	(SEQ ID NO: 69)	(SEQ ID NO: 1258)
25	AR	GTAGTATTCGAAGGTAGT
	(SEQ ID NO: 69)	(SEQ ID NO: 1255)
26	AR	GTAGTATTTGAAGGTAGT
	(SEQ ID NO: 69)	(SEQ ID NO: 1256)
27	AR	ACTACCTTCGAATACTAC
- ·		(SEQ ID NO: 1257)
	(SEQ ID NO: 69)	(PDG ID 140, 1731)

No:	Gene	Oligo:
28	AR	ACTACCTTCAAATACTAC
	(SEQ ID NO: 69)	(SEQ ID NO: 1258)

Table 3: Hybridisation oligonucleotides used in differentiation between acute lymphocytic leukemia and acute myelogenous leukemia (Figures 5 and 6).

No:	Gene	Oligo:
1	N33	TTGGTTCGGGAAAGGTAA
	(SEQ ID NO: 46)	(SEQ ID NO: 1211)
2	N33	TTGGTTTGGGAAAGGTAA
	(SEQ ID NO: 46)	(SEQ ID NO: 1212)
3	EGR4	GGTGGGAAGCGTATTTAT
	(SEQ ID NO: 26)	(SEQ ID NO: 1213)
4	EGR4	GGTGGGAAGTGTATTTAT
	(SEQ ID NO: 26)	(SEQ ID NO: 1214)
5	N33	ATTTAGTTCGGGGGAGGA
	(SEQ ID NO: 46)	(SEQ ID NO: 1215)
6	N33	ATTTAGTTTGGGGGAGGA
	(SEQ ID NO: 46)	(SEQ ID NO: 1216)
7	N33	TGTTATTTCGGAGGGTTT
ł	(SEQ ID NO: 46)	(SEQ ID NO: 1217)
8	N33	TGTTATTTTGGAGGGTTT
	(SEQ ID NO: 46)	(SEQ ID NO: 1218)
9	N33	GTTTAGTTAGCGGGTTTT
	(SEQ ID NO: 46)	(SEQ ID NO: 1219)
10	N33	GTTTAGTTAGTGGGTTTT
	(SEQ ID NO: 46)	(SEQ ID NO: 1220)
11	Humos	GAGTTTAACGTAGTAAGG
	(SEQ ID NO: 71)	(SEQ ID NO: 1221)
12	Humos	GAGTTTAATGTAGTAAGG
	(SEQ ID NO: 71)	(SEQ ID NO: 1222)
13	SDC4	TTTTATATCGGGTGTGTT
	(SEQ ID NO: 54)	(SEQ ID NO: 1223)
14	SDC4	TTTTATATTGGGTGTGTT
	(SEQ ID NO: 54)	(SEQ ID NO: 1224)

Gene	Oligo:
EGR4	TTATAGTTCGAGTTTTTT
(SEQ ID NO: 26)	(SEQ ID NO: 1225)
EGR4	TTATAGTTTGAGTTTTT
(SEO ID NO: 26)	(SEQ ID NO: 1226)
MPL	TGTAGTGAGTCGAGATTA
(SEO ID NO: 42)	(SEQ ID NO: 1227)
MPL	TGTAGTGAGTTGAGATTA
(SEO ID NO: 42)	(SEQ ID NO: 1228)
TP73	TTTGGTGCGCGTAGAGAA
(SEO ID NO: 58)	(SEQ ID NO: 1229)
TP73	TTTGGTGTGTGTAGAGAA
(SEO ID NO: 58)	(SEQ ID NO: 1230)
	GGGTTTTACGTGGTTGGA
	(SEQ ID NO: 1231)
	GGGTTTTATGTGGTTGGA
	,
	(SEQ ID NO: 1232) TCCAACCACGTAAAACCC
	(SEQ ID NO: 1233) TCCAACCACATAAAACCC
	(SEQ ID NO: 1234) GGGTTTTACGTGGTTGGA
	(SEQ ID NO: 1231) GGGTTTTATGTGGTTGGA
	(SEQ ID NO: 1232)
CDK 4	TCCAACCACGTAAAACCC
(SEQ ID NO: 70)	(SEQ ID NO: 1233)
CDK 4	TCCAACCACATAAAACCC
(SEQ ID NO: 70)	(SEQ ID NO: 1234)
CSNK2B	AGGAGTTTCGGAGGAAAT
(SEQ ID NO: 67)	(SEQ ID NO: 1235)
CSNK2B	AGGAGTTTTGGAGGAAAT
(SEQ ID NO: 67)	(SEQ ID NO: 1236)
CSNK2B	ATTTCCTCCGAAACTCCT
(SEQ ID NO: 67)	(SEQ ID NO: 1237)
CSNK2B	ATTTCCTCCAAAACTCCT
(SEQ ID NO: 67)	(SEQ ID NO: 1238)
CSNK2B	AGGAGTTTCGGAGGAAAT
(SEQ ID NO: 67)	(SEQ ID NO: 1235)
	EGR4 (SEQ ID NO: 26) EGR4 (SEQ ID NO: 26) MPL (SEQ ID NO: 42) MPL (SEQ ID NO: 42) TP73 (SEQ ID NO: 58) TP73 (SEQ ID NO: 58) CDK 4 (SEQ ID NO: 70) CSNK2B (SEQ ID NO: 67) CSNK2B

No:	Gene	Oligo:
34	CSNK2B	AGGAGTTTTGGAGGAAAT
	(SEQ ID NO: 67)	(SEQ ID NO: 1236)
35	CSNK2B	ATTTCCTCCGAAACTCCT
	(SEQ ID NO: 67)	(SEQ ID NO: 1237)
36	CSNK2B	ATTTCCTCCAAAACTCCT
	(SEQ ID NO: 67)	(SEQ ID NO: 1238)
37	Tubulin	TACTTCAACGTCCTAATA
	(SEQ ID NO: 64)	(SEQ ID NO: 1239)
38	Tubulin	TACTTCAACATCCTAATA
	(SEQ ID NO: 64)	(SEQ ID NO: 1240)
39	Tubulin	TACTTCAACGTCCTAATA
	(SEQ ID NO: 64)	(SEQ ID NO: 1239)
40	Tubulin	TACTTCAACATCCTAATA
	(SEQ ID NO: 64)	(SEQ ID NO: 1240)
41	ELK1	TTTATTTTCGTTTTTGGG
	(SEQ ID NO: 63)	(SEQ ID NO: 1241)
42	ELK1	TTTATTTTTGTTTTTGGG
	(SEQ ID NO: 63)	(SEQ ID NO: 1242)
43	ELK1	TTTATTTTCGTTTTTGGG
	(SEQ ID NO: 63)	(SEQ ID NO: 1241)
44	ELK1	TTTATTTTTGTTTTTGGG
	(SEQ ID NO: 63)	(SEQ ID NO: 1242)
. 45	Tubulin	CAATACTCCGTCCATTAA
	(SEQ ID NO: 64)	(SEQ ID NO: 1243)
46	Tubulin	CAATACTCCATCCATTAA
	(SEQ ID NO: 64)	(SEQ ID NO: 1244)
47	Tubulin	CAATACTCCGTCCATTAA
	(SEQ ID NO: 64)	(SEQ ID NO: 1243)
48	Tubulin	CAATACTCCATCCATTAA
	(SEQ ID NO: 64)	(SEQ ID NO: 1244)

Description of the figures:

Figure 1

Methylation analysis and quantification of two CpG dinucleotides in exon 14 of the human factor VIII gene. For calibration purpose a series of hybridisations was performed with mixtures of artificially up- and down-methylated DNA fragments of the factor VIII exon 14 gene. Down- and up-methylated DNA fragments were mixed at ratios: 0:3, 1:2, 2:1, 3:0, representing a methylation status of 100 %, 66 %, 33 % and 0 %, respectively. A, Methylation detection by oligonucleotide microarray hybridisation. The fluorescence signals of the CG and TG version of the factor VIII exon 14 oligonucleotides F8-5 (TTATTAACGGGAAATAAT, TTATTAATGGGAAATAAT) and F8-3 (AATAAGTTCGAAATAGAA, AATAAGTTTGAAATAGAA) are shown which are generated by samples reflecting the methylation status 0%, 33 %, 66 % and 100 %. The hybridisation signals are shown in a false colour image with the colours blue, green and yellow indicating fluorescence signal ranges at 635nm of 200 to 800, 800 to 2000 and 2000 to 8000, respectively. B. Distribution and probability of measurement of different methylated DNA for two CpG methylation. For each CpG position two kinds of detection oligomers were used. Oligomers that hybridise if the CpG was methylated are referred to as CG oligomer and the oligomers that hybridise if the CpG was not methylated are referred to as TG oligos. For the 4 kinds of compounds 59, 36, 40, 63 identical slides were made. The log-ratio of the CG and the TG detection oligomer hybridisation intensity was calculated and then averaged for experimental subgroups each containing 3 identical experiments. The distribution function of the CG:TG ratios shows that measurement values of the different mixtures are well separated and therefore allow a high resolution detection of the methylation level of a single CpG. This is an essential prerequisite for methylation dependent class prediction or class discovery. Taking into account only the 100% and 0% methylated DNA, and averaging for the 22 CpG sites investigated in the calibration experiments, the average error for methylation detection is 4%. The log-ratios are not grouped symmetrically around zero but shifted towards negative values. We assume that the energetically different effects of G-T and A-C mismatches allow hybridisations of the methylated allele to the oligonucleotide representing the unmethylated more easily than vice versa.

Figure 2

Gender separation and CpG site prediction of ALL/AML and healthy blood samples. High probability of methylation corresponds to red, uncertainty to black and low probability to green. The labels on the left side of the plot are gene and CpG identifiers. The labels on the right side give the significance of the difference between the means of the two groups. Each row corresponds to a single CpG and each column to the methylation levels of one sample.

Figure 2 A, Gender separation. The 20 CpG sites with the most significant difference between female and male samples are shown. Only non cell lines were used. As expected the significant CpG dinucleotides come from the two X-chromosome genes (ELK1, AR).

Figure 2 B, differentiation between healthy samples and ALL samples. The 54 CpG sites with the most significant difference between healthy and ALL samples are shown. High probability of methylation corresponds to red, uncertainty to black and low probability to green. The labels on the left side of the plot are gene and CpG identifiers. The labels on the right side give the significance of the difference between the means of the two groups. Each row corresponds to a single CpG and each column to the methylation levels of one sample.

Figure 2C shows Figure 2B reproduced in greyscale. High probability of methylation corresponds to black, uncertainty to grey and low probability to white. The labels on the left side of the plot are gene and CpG identifiers. The labels on the right side give the significance of the difference between the means of the two groups. Each row corresponds to a single CpG and each column to the methylation levels of one sample.

Legend Fig. 2A:

- 1 Female
- 2 Male

Legend Fig. 2B:

1 – Healthy

2 – ALL

Figure 3

Class prediction and class discovery of blood samples. A, The plot shows a SVM trained on the two most significant CpG sites for the healthy-ALL discrimination using all available healthy and ALL samples as training data. The red points are healthy, the yellow ones ALL samples. Circled points are the support vectors defining the white borderline between the green area of healthy prediction and the blue area of an ALL prediction. The colour intensity corresponds to the prediction strength. B, Support Vector Machine for AML-ALL class prediction. The red points are ALL, the yellow ones AML samples. Circled points are the support vectors defining the white borderline between the green area of ALL prediction and the blue area of an AML prediction. The colour intensity corresponds to the prediction strength.C, Class discovery. The figure shows a hierarchical clustering of all available samples. Healthy individuals are coloured green, patients with ALL red and patients with AML blue. Asterisks indicate cell line samples. The feature space consisted of all CpG sites except those from the two X-chromosomal genes. The diagnosis was unknown to the algorithm.

Legend Fig. 3A: t-Test ranking (Healty vs. ALL)

1 – MOS CpG5

2 - CDK4 CpG2

Legend Fig. 3B: t-Test ranking (AML vs. ALL)

1 - CSNK2B CpG2

2 - CDK4 CpG10

Figure 4: Feature selection methods.

Figure 4A shows principle component analysis. The entire data set was projected onto its first 2

principle components. Circles represent cell lines, triangles primary patient tissue. Filled circles or

triangles are AML, empty ones ALL samples.

Figure 4B: Fisher criterion. The 20 highest ranking CpG sites according to the Fisher criterion are

shown. The highest ranking features are on the bottom of the plot. High probability of methylation

corresponds to black, uncertainty to grey and low probability to white.

Figure 4C: Two sample t-test.

Figure 4D: Backward elimination.

Figure 5

Support Vector Machine on two best features of the Fisher criterion. The plot shows a SVM trained

on the two highest ranking CpG sites according to the Fisher criterion with all ALL and AML

samples used as training data. The black points are AML, the grey ones ALL samples. Circled

points are the support vectors defining the white borderline between

the areas of AML and ALL prediction. The grey value of the background corresponds to the

prediction strength.

Legend Fig. 5:

1 – CSNK2B CpG2

2 – CDK4 CpG3

Figure 6

Dimension dependence of feature selection performance. The plot shows the generalisation performance of a linear SVM with four different feature selection methods against the number of selected features. The x-axis is scaled logarithmically and gives the number of input features for the SVM, starting with two. The y-axis gives the achieved generalisation performance. Note that the maximum number of principle components corresponds to the number of available samples.

Legend Fig. 6:

- 1 Feature Number
- 2 Test Error
- 3 Fisher Criterion
- 4 t-Test
- 5 Backward Elimination
- 6 Principle Component Analysis

Figure 7A

Differentiation between healthy and acute lymphocytic leukemia samples. High probability of methylation corresponds to red, uncertainty to black and low probability to green. The labels on the left side of the plot are gene and CpG identifiers. The labels on the right side give the significance of the difference between the means of the two groups. Each row corresponds to a single CpG and each column to the methylation levels of one sample.

Figure 7B

Differentiation between acute lymphocytic leukemia and acute myelogenous leukemia, this shows Figure 7A reproduced in greyscale. High probability of methylation

corresponds to black, uncertainty to grey and low probability to white. The labels on the left side of the plot are gene and CpG identifiers. The labels on the right side give the significance of the difference between the means of the two groups. Each row corresponds to a single CpG and each column to the methylation levels of one sample.

Figure 8A

Differentiation between acute lymphocytic leukemia and acute myelogenous leukemia samples. High probability of methylation corresponds to red, uncertainty to black and low probability to green. The labels on the left side of the plot are gene and CpG identifiers. The labels on the right side give the significance of the difference between the means of the two groups. Each row corresponds to a single CpG and each column to the methylation levels of one sample.

Figure 8B

Differentiation between acute lymphocytic leukemia and acute myelogenous leukemia, this shows Figure 8A reproduced in greyscale. High probability of methylation corresponds to red, uncertainty to black and low probability to green. The labels on the left side of the plot are gene and CpG identifiers. The labels on the right side give the significance of the difference between the means of the two groups. Each row corresponds to a single CpG and each column to the methylation levels of one sample.

Patent Claims

- 1. A method for detecting and differentiating between hematopoietic cell proliferative disorders associated with at least one gene and/or their regulatory regions from the group comprising ABL1, ABL1, APAF1, APC, AR, ARHI, BAK1, BAX, BCL2, CASP10, CASP8, CASP9, CCND2, CDC2, CDC25A, CDH1, CDH3, CDK 4, CDKN1A, CDKN1B (p27 Kip1), CDKN1C, CDKN2a, CDKN2B, CSNK2B, DAPK1, EGR4, ELK1, ESR1, FOS, GPIb beta, GPR37, GSK3ß, GSTP1, HIC-1, HOXA5, IGF2, MDR1, MGMT, MLH1, MOS, Humos, MPL, MYC, MYCL1, MYOD1, N33, PITX2, PML, PMS2, PRAME, PTEN, RB1, RBL2, SDC4, SFN, TCL1A, TGFBR2, TP73, WT1, N-MYC, L-MYC, C-ABL, ELK1, Tubulin, CSF1, CD1R3, CSNK2B, Me491/TD63, AR, CDK 4, Humos, CDC25A, CMYCex3 in a subject, said method comprising contacting a target nucleic acid in a biological sample obtained from said subject with at least one reagent or a series of reagents, wherein said reagent or series of reagents, distinguishes between methylated and non methylated CpG dinucleotides within the target nucleic acid.
- 2. A method according to claim 1 wherein,

said method differentiates between normal hematopoietic cells and proliferative disorder hematopoietic cells.

3. A method according to claim 1 wherein,

said method differentiates between acute lymphocytic leukemia and acute myelogenous leukemia.

4. Use of methods according to claims 1 and 2 wherein,

CDKN1A, CDK 4, AR are used to differentiate between healthy hematopoietic cells and proliferative disorder hematopoietic cells.

5. Use of methods according to claims 1 and 3 wherein,

N33,EGR4, CDKN1A, SDC4, MPL, TP 73, BAK 1, CSNK2B, ARHI, CASP 10 are used to

differentiate between acute lymphocytic leukemia and acute myelogenous leukemia.

- 6. A method according to any one of Claims 1 to 5 comprising the following steps:
 - obtaining a biological sample containing genomic DNA
 - extracting the genomic DNA
 - converting cytosine bases in the genomic DNA sample which are unmethylated at the 5-position, by treatment, to uracil or another base which is dissimilar to cytosine in terms of base pairing behaviour;
 - fragments of the pretreated genomic DNA are amplified
 - identification of the methylation status of one or more cytosine positions
- 7. The method according to claim 6, characterized in that the reagent is a solution of bisulfite, hydrogen sulfite or disulfite.
- 8. The method as recited in Claims 6 and 7, characterized in that the amplification is carried out by means of the polymerase chain reaction (PCR).
- 9. The method as recited in one of the Claims 6 through 8, characterized in that the amplification is carried out by means of a heat-resistant DNA polymerase.
- 10. The method as recited in one of the Claims 6 through 9, characterized in that more than ten different fragments having a length of 100 2000 base pairs are amplified.
- 11. The method as recited in one of claims 6 through 10,

wherein the amplification step is carried out using a set of primer oligonucleotides

comprising SEO ID NO: 387 through SEO ID NO: 534.

- 12. The method as recited in one of the Claims 6 through 11, characterized in that the amplification of several DNA segments is carried out in one reaction vessel.
- 13. The method as recited in one of Claims 6 through 12, characterized in that the amplification step preferentially amplifies DNA which is of particular interest in healthy and/or diseased hematopoietic cells, based on the specific genomic methylation status of hematopoietic cells, as opposed to background DNA.
- 14. The method according to one of Claims 6 through 13,

characterized in that the methylation status within at least one gene and/or their regulatory regions from the group comprising ABL1, ABL1, APAF1, APC, AR, ARHI, BAK1, BAX, BCL2, CASP10, CASP8, CASP9, CCND2, CDC2, CDC25A, CDH1, CDH3, CDK 4, CDKN1A, CDKN1B (p27 Kip1), CDKN1C, CDKN2a, CDKN2B, CSNK2B, DAPK1, EGR4, ELK1. ESR1. FOS. GPIb beta, GPR37, GSK3B, GSTP1, HIC-1, HOXA5, IGF2, MDR1, MGMT, MLH1, MOS, Humos, MPL, MYC, MYCL1, MYOD1, N33, PITX2, PML, PMS2, PRAME, PTEN, RB1, RBL2, SDC4, SFN, TCL1A, TGFBR2, TP73, WT1, N-MYC, L-MYC, C-ABL, ELK1, Tubulin, CSF1, CD1R3, CSNK2B, Me491/TD63, AR, CDK 4, Humos, CDC25A, CMYCex3 is detected by hybridization of each amplificate to an oligonucleotide or peptide nucleic acid (PNA)-oligomer.

15.A method according to Claim 14,

characterized in that the olignonucleotide or peptide nucleic acid (PNA)-oligomer is taken from the group comprising SEQ ID NO: 535 to SEQ ID NO: 1258.

16. The method according to Claims 6 through 15,

characterized in that the amplificates are labelled.

- 17. The method as recited in Claim 16, characterized in that the labels of the amplificates are fluorescence labels.
- 18. The method as recited in Claim 16, characterized in that the labels of the amplificates are radionuclides.
- 19. The method as recited in Claims 16, characterized in that the labels of the amplificates are detachable molecule fragments having a typical mass which are detected in a mass spectrometer.
- 20. The method as recited in one of the Claims 6 through 19, characterized in that the amplificates or fragments of the amplificates are detected in the mass spectrometer.
- 21. The method as recited in one of the Claims 19 and 20, characterized in that the produced fragments have a single positive or negative net charge.
- 22. The method as recited in one of the Claims 19 through 21, characterized in that detection is carried out and visualized by means of matrix assisted laser desorption/ionization mass spectrometry (MALDI) or using electron spray mass spectrometry (ESI).
- 23.A method according to Claims 1 through 5 comprising the following steps;
 - a) obtaining a biological sample containing genomic DNA
 - b) extracting the genomic DNA
 - c) digesting the genomic DNA comprising at least one or more CpGs of the genes ABL1, ABL1, APAF1, APC, AR, ARHI, BAK1, BAX, BCL2, CASP10, CASP8, CASP9, CCND2, CDC2,

CDC25A, CDH1, CDH3, CDK 4, CDKN1A, CDKN1B (p27 Kip1), CDKN1C, CDKN2a, CDKN2B, CSNK2B, DAPK1, EGR4, ELK1, ESR1, FOS, GPIb beta, GPR37, GSK3ß, GSTP1, HIC-1, HOXA5, IGF2, MDR1, MGMT, MLH1, MOS, Humos, MPL, MYC, MYCL1, MYOD1, N33, PITX2, PML, PMS2, PRAME, PTEN, RB1, RBL2, SDC4, SFN, TCL1A, TGFBR2, TP73, WT1, N-MYC, L-MYC, C-ABL, ELK1, Tubulin, CSF1, CD1R3, CSNK2B, Me491/TD63, AR, CDK 4, Humos, CDC25A, CMYCex3 with one or more methylation sensitive restriction enzymes

- d) detection of the DNA fragments generated in the digest of step c)
- 24.A method according to Claim 23,

wherein the DNA digest is amplified prior to Step d.

- 25. The method as recited in Claim 24, characterized in that the amplification is carried out by means of the polymerase chain reaction (PCR).
- 26. The method as recited in one of the Claims 24 and/or 25, characterized in that the amplification of more than one DNA fragments is carried out in one reaction vessel.
- 27. The method as recited in one of the Claims 24 through 26, characterized in that the polymerase is a heat-resistant DNA polymerase.
- 28.An isolated nucleic acid of a pretreated genomic DNA according to one of the sequences taken from the group comprising SEQ ID NO: 95 to SEQ ID NO: 386 and sequences complementary thereto.
- 29.An oligomer, in particular an oligonucleotide or peptide nucleic acid (PNA)-oligomer, said oligomer comprising at least one base sequence of at least 10 nucleotides which hybridizes to or

- is identical to a pretreated genomic DNA according to one of the SEQ ID NO:95 to SEQ ID NO: 386 according to Claim 28.
- 30. The oligonucleotide as recited in Claim 29; wherein the base sequence includes at least one CpG or TpG dinucleotide sequence.
- 31. The oligonucleotide as recited in Claim 30; characterized in that the cytosine of the at least one CpG or TpG dinucleotide is/are located approximately in the middle third of the oligomer.
- 32. An oligomer, in particular an oligonucleotide or peptide nucleic acid (PNA)-oligomer, according to one of the sequences taken from the group comprising SEQ ID NO: 535 to SEQ ID NO: 1258.
- 33.A set of oligonucleotides, comprising at least two oligonucleotides according to any of Claims 29 through 32.
- 34.A set of oligonucleotides, comprising at least two oligonucleotides according to SEQ ID NO: 1245, 1245, 1246, 1246, 1247, 1247, 1248, 1248, 1249, 1249, 1250, 1250, 1251, 1251, 1252, 1252, 1253, 1253, 1254, 1254, 1255, 1255, 1256, 1256, 1257, 1257, 1258, 1258.
- 35. One or more isolated nucleic acid(s) taken from the group according to SEQ ID NO: 88 94.
- 36.A set of oligonucleotides, comprising at least two oligonucleotides according to SEQ ID NO: 1211 1231, 1231, 1232, 1232, 1233, 1233, 1234, 1234, 1235, 1235, 1236, 1236, 1237, 1237, 1238, 1238, 1239, 1239, 1240, 1240, 1241, 1241, 1242, 1242, 1243, 1243, 1244, 1244.
- 37. One or more isolated nucleic acid(s) taken from the group according to SEQ ID NO: 74 87.
- 38. A set of oligomers, peptide nucleic acid (PNA)-oligomers and/or isolated nucleic acids as recited in Claims 33 through 35 and 37, comprising oligomers for detecting the methylation state

of all CpG dinucleotides within one or more of the sequences according to SEQ ID NO: 1 to SEQ ID NO: 73 and sequences complementary thereto.

- 39.Use of a set of oligomers or peptide nucleic acid (PNA)-oligomers according to any of claims 29 through 34 and 36 as probes for determining the cytosine methylation state and/or single nucleotide polymorphisms (SNPs) of sequences according to SEQ ID NO: 1 to SEQ ID NO: 73 and sequences complementary thereto.
- 40.Use of a set of oligonucleotides according to Claim 34 or nucleic acid(s) according to Claim 35 for the differentiation between healthy hematopoietic cells and proliferative disorder hematopoietic cells.
- 41.Use of a set of oligonucleotides according to Claim 36 or nucleic acid(s) according to Claim 37 for the differentiation between acute lymphocytic leukemia and acute myelogenous leukemia.
- 42.A set of at least two oligonucleotides or peptide nucleic acid (PNA)-oligomers as recited in Claim 29, as primer oligonucleotides for the amplification of DNA sequences of one of SEQ ID NO: 95 to SEQ ID NO: 386 according to Claim 28 and/or sequences complementary thereto and segments thereof.
- 43. Use of a pretreated genomic DNA according to claim 28 for the determination of the methylation status of a corresponding genomic DNA and/or detection of single nucleotide polymorphisms (SNPs).
- 44. A set of oligonucleotides or peptide nucleic acid (PNA)-oligomers as recited in Claims 33, 34 or 36 characterized in that at least one oligonucleotide is bound to a solid phase.
- 45. A set of oligonucleotides or peptide nucleic acid (PNA)-oligomers as recited in Claims 33, 34 or 36, characterized in that all members of the set are bound to a solid phase.

- 46.A method for manufacturing an arrangement of different oligomers or peptide nucleic acid (PNA)-oligomers (array) for analyzing diseases associated with the corresponding genomic methylation status of the CpG dinucleotides within one of the SEQ ID NO: 1 to SEQ ID NO: 73 and sequences complementary thereto, wherein at least one oligomer according to any of the Claims 33, 34 or 36 is coupled to a solid phase.
- 47.An arrangement of different oligomers or peptide nucleic acid (PNA)-oligomers (array) obtainable according to claims 44 and 45.
- 48. An array of different oligonucleotide- and/or PNA-oligomer sequences as recited in Claim 47, characterized in that these are arranged on a plane solid phase in the form of a rectangular or hexagonal lattice.
- 49. A nucleic acid or peptide nucleic acid array for the analysis of hematopoietic cell proliferative disorders associated with the methylation state of genes comprising at least one nucleic acid according to one of the preceding claims.
- 50. The array as recited in any of the Claims 47 through 49, characterized in that the solid phase surface is composed of silicon, glass, polystyrene, aluminium, steel, iron, copper, nickel, silver, or gold.
- 51.A kit comprising a bisulfite (= disulfite, hydrogen sulfite) reagent as well as oligonucleotides and/or PNA-oligomers according to one of the Claims 29 through 37.
- 52. The use of oligonucleotides or peptide nucleic acid (PNA)-oligomers according to SEQ ID NO: 74 to SEQ ID NO: 94 and SEQ ID NO: 535 to SEQ ID NO: 1258 for the detection of a predisposition to, differentiation between subclasses, diagnosis, prognosis, treatment and/or

monitoring of hematopoietic cell proliferative disorders.

53.A DNA sequence according to one of the sequences taken from the group comprising SEQ ID NO: 95 to SEQ ID NO: 386 and sequences complementary thereto for use in the analysis of cytosine methylation within said nucleic acid for the detection of a predisposition to, differentiation between subclasses, diagnosis, prognosis, treatment and/or monitoring of hematopoietic cell proliferative disorders.

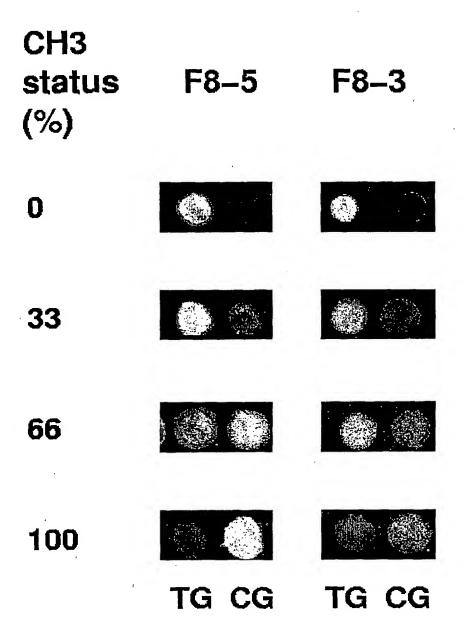


Figure 1A

Figure 1B



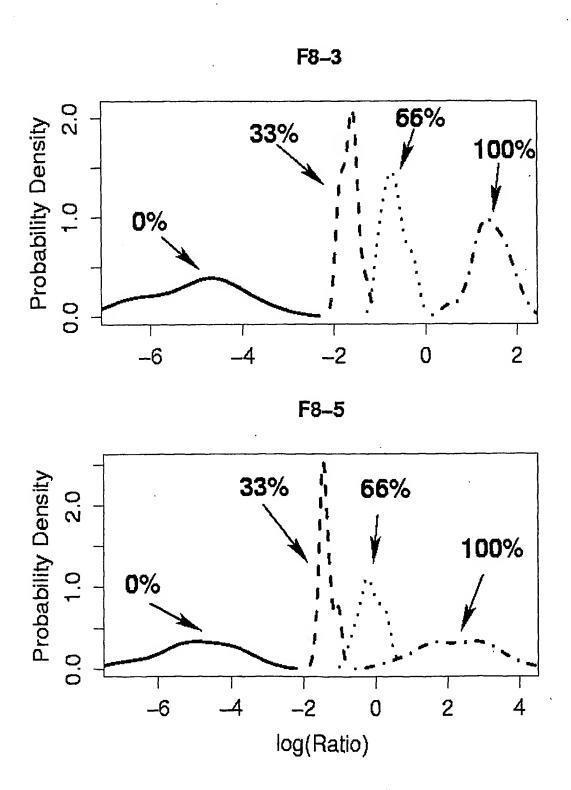


Figure 2A

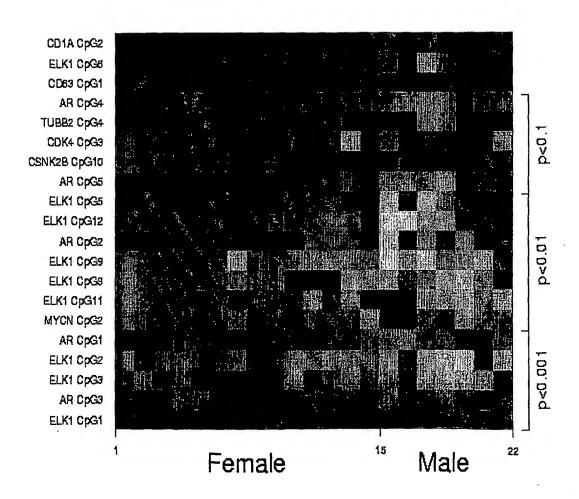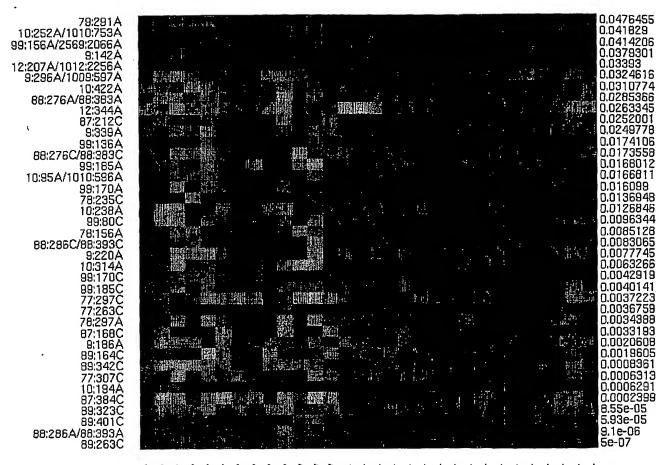

4/18



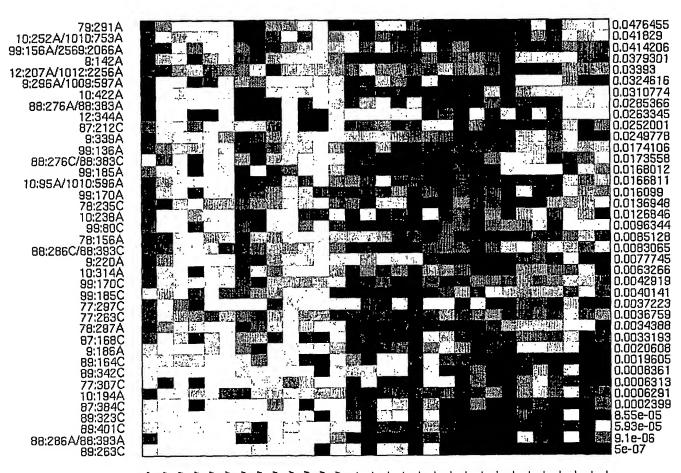


Figure 3A

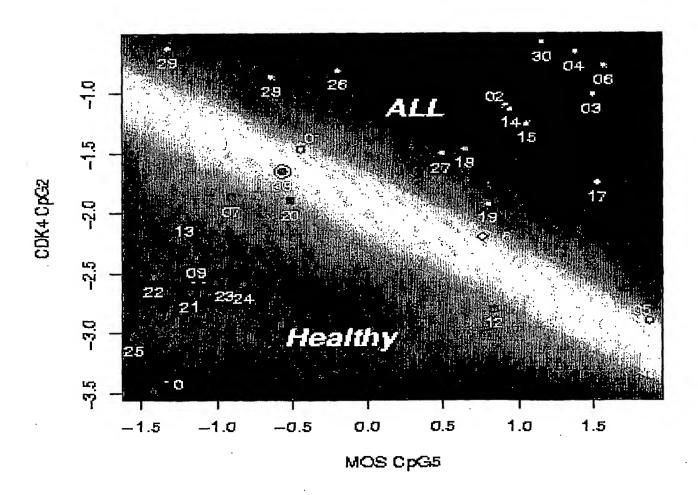
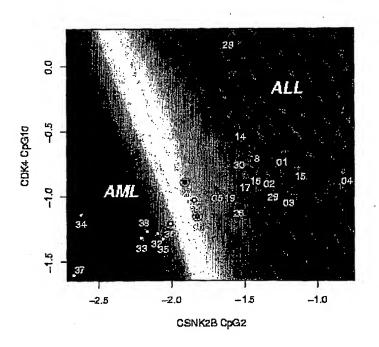


Figure 3B



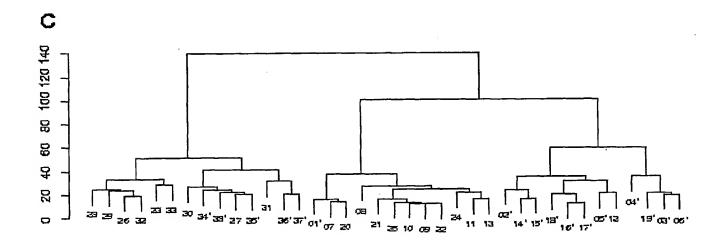


Figure 3C

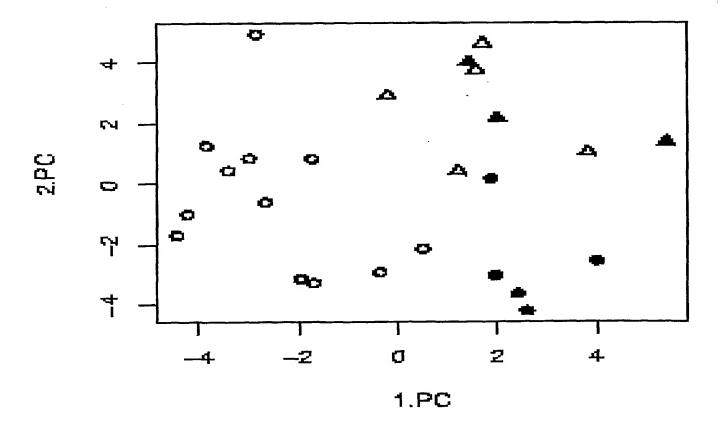


Figure 4A

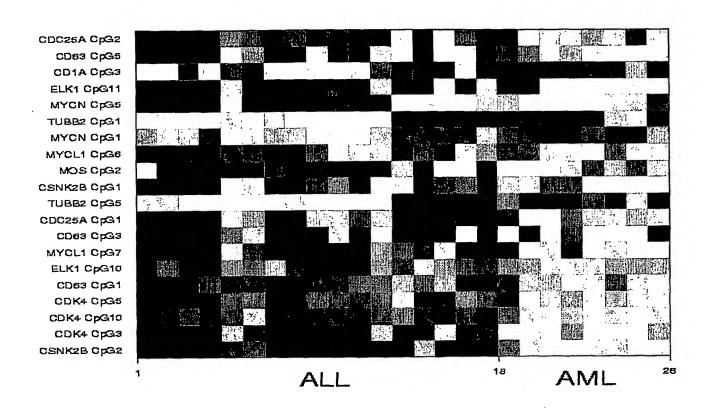


Figure 4B

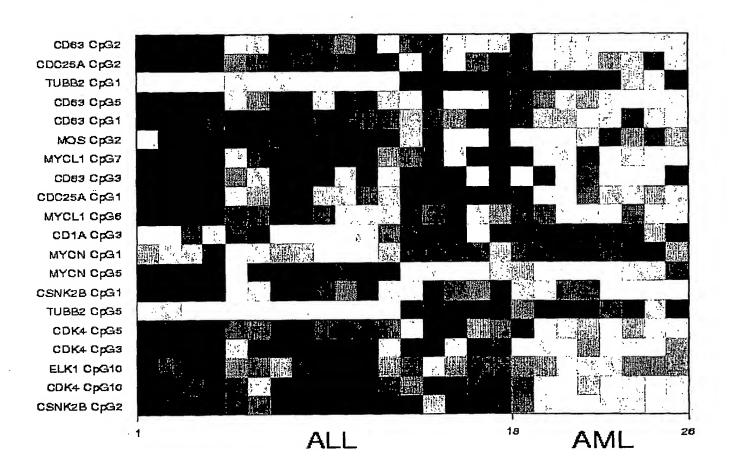


Figure 4C

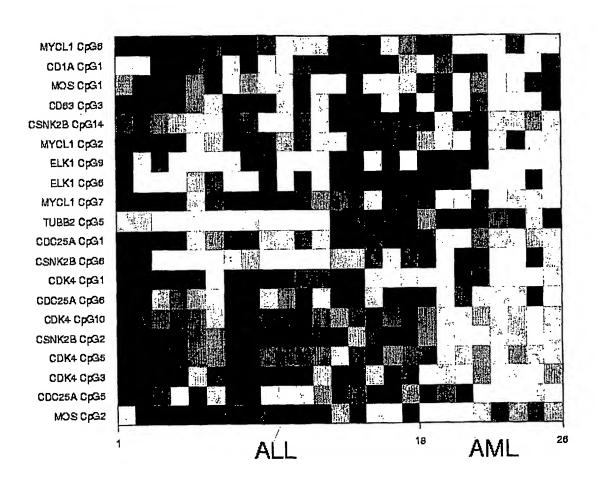


Figure 4D

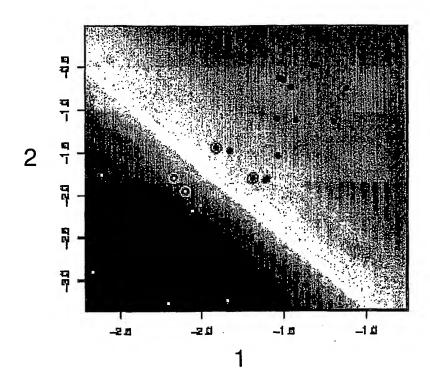


Fig. 5

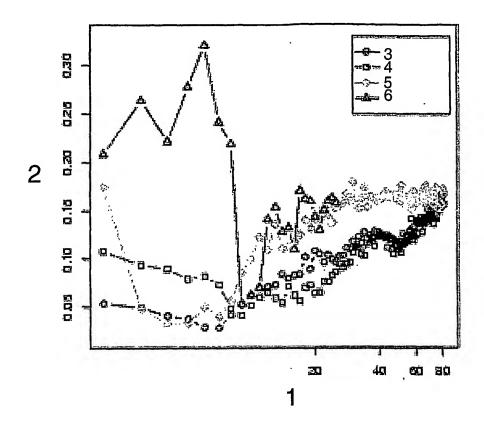


Fig. 6


15/18

77:263C 78:297A 99:80C 99:136C 10:304A 9:186A 12:80A 10:194A 67:361C 79:77A 12:344A 78:67C 79:302A 77:191C 12:107A/1012:2156A 78:201C 99:156A/2569:2066A 10:381A 99:136A 10:95A/1010:596A 9:220A 89:164C 79:390A/1079:479A 79:175A 79:291A 88:348A/88:456A 9:142A 78:67A 88:433A/88:540A 88:276A/88:383A 67:251C 12:207A/1012:2256A 67:155C 78:156A 88:433C/88:540C



0.0490278 0.0412449 0.0386675 0.0354139 0.0354059 0.0333188 0.0310073 0.0303154 0.0287562 0.0263738 0.0248328 0.0230314 0.0158705 0.0146667 0.014391 0.0120971 0.0110565 0.0038815 0.0032655 0.0031618 0.002885 0.0027582 0.0023383 0.0022992 0.0017879 0.0013198 0.0011476 0.0007615 0.0006906 0.0006357 0.0002499 0.0002047 0.00014 4.39e-05 2.8e-05

77:263C 78:297A 89:80C 99:136C 10:304A 9:186A 12:80A 10:194A 67:361C 79:77A 12:344A 78:67C 79:302A 77:191C 12:107A/1012:2156A 78:201C 99:156A/2569:2066A 10:381A 88:136A 10:95A/1010:596A 9:220A 89:164C 78:390A/1079:479A 79:175A 79:291A 88:349A/88:456A 9:142A 78:67A 88:433A/88:540A 88:276A/88:383A 67:251C 12:207A/1012:2256A 67:155C 78:156A 88:433C/88:540C

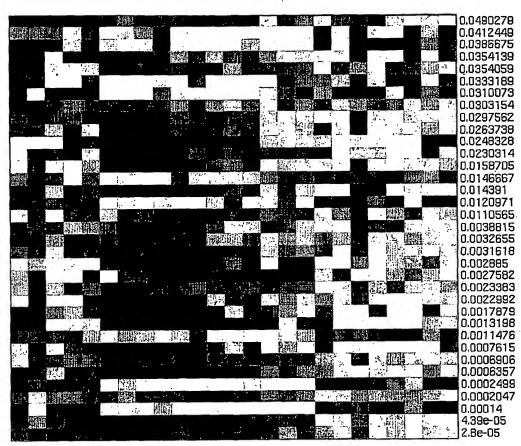

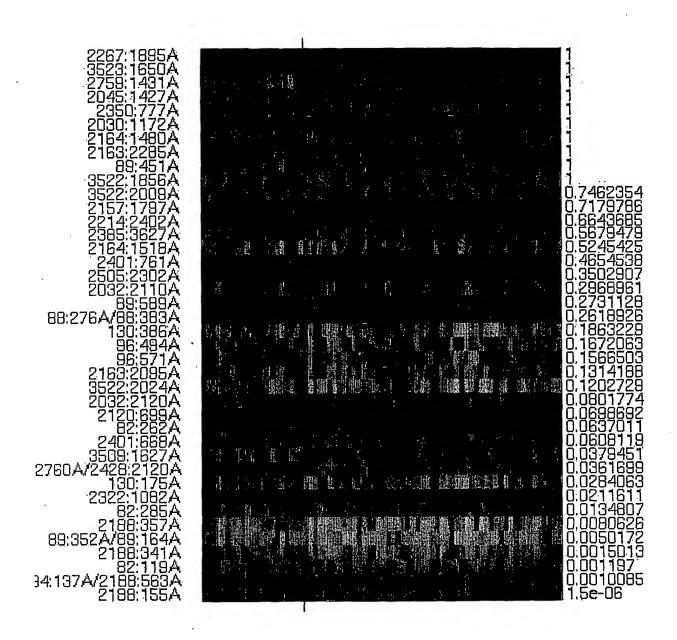


Figure 8A

